



TrioDocs

Version: 0.4.0

Date: July 03, 2026

Download the latest version on:
<https://triodocs.org>


Insulin Pump

Insulin Pump

This section allows you to choose your insulin pump. Most options are self-explanatory. For more information on compatible insulin pumps, please see the following [link](#)

Step 1: Add Pump

The first step in setting up your insulin pump on Trio is to tap the "Add Pump" button in the Devices Menu

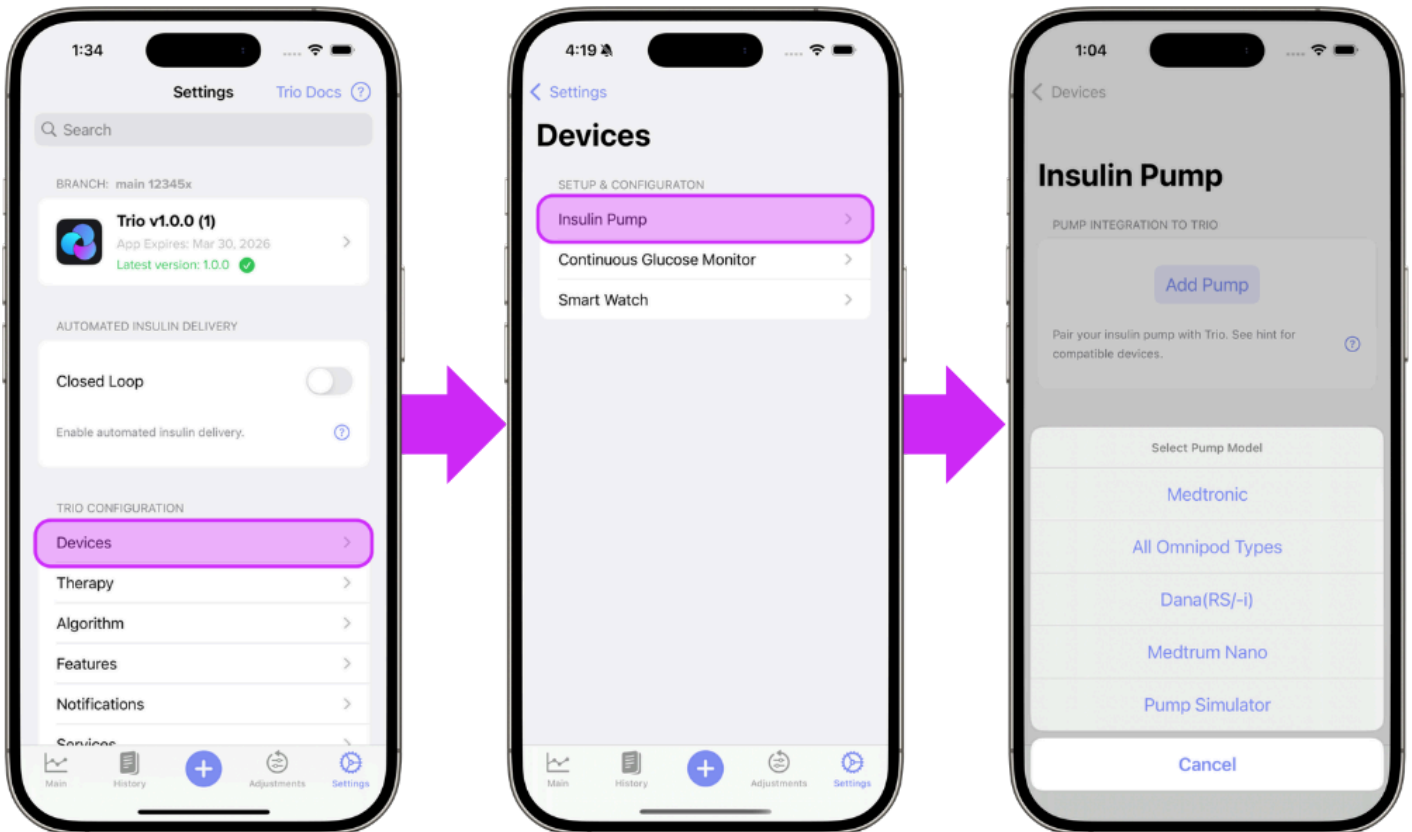

Step 2: Select Your Pump

Select your pump from the in-app menu and from the options below for step-by-step instructions. The links below will guide you through the connection instructions for your specific pump.

- [All Omnipod Types](#)
 - [Omnipod Eros](#)
 - [Omnipod Dash](#)
 - [Omnipod 5](#)
- [Dana\(RS/-i\)](#)
- [Medtrum Nano](#)
- [Medtronic](#)
- [Pump Simulator](#)

All Omnipod Types

Step 3 Select which Omnipod type you will use with Trio

The screenshot shows a mobile app interface for selecting a pod type. At the top, the status bar shows the time 1:09 and battery level. The app header has a back arrow, the text 'Insulin Type', and a 'Cancel' button. The main title is 'Pod Type'. Below it is an instructional text box: 'Select the particular Omnipod Pod type to use with Trio. Be sure to select the correct Pod type or Trio will not be able to communicate with it.' There are three options, each with a radio button: 'Omnipod Classic' (described as also known as Eros Pods, requiring a RileyLink), 'Omnipod DASH' (described as having a blue needle tab), and 'Omnipod 5' (described as having a clear needle tab and SmartAdjust technology). At the bottom is a large blue button labeled 'Confirm Pod Type'.

1:09

< Insulin Type Cancel

Pod Type

Select the particular Omnipod Pod type to use with Trio. Be sure to select the correct Pod type or Trio will not be able to communicate with it.

Omnipod Classic
Omnipod Classic Pods (also known as Eros Pods) have a clear needle tab with a 6-character LOT number starting with 'L'. These Pods require the use of a RileyLink to communicate with the iPhone.


Omnipod DASH
Omnipod DASH Pods have a blue needle tab with a 12-character LOT number typically starting with 'PD1'.

Omnipod 5
Omnipod 5 Pods have a clear needle tab with a 12-character LOT number typically starting with 'PH1'. The Pod's "SmartAdjust" technology will not be used for closed loop control.


Confirm Pod Type

Omnipod Eros

Step 4 Begin the process of configuring your Eros pod by tapping the blue "Continue" button at the bottom

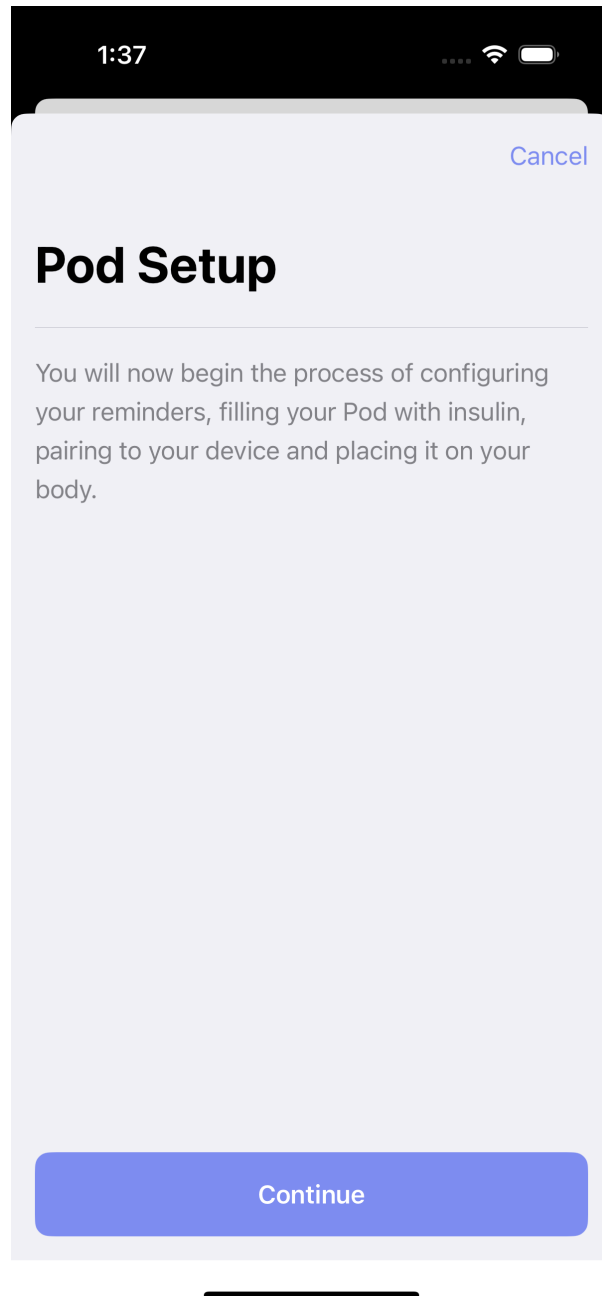

Step 5 Choose the expiration reminder that you prefer and tap "Next"

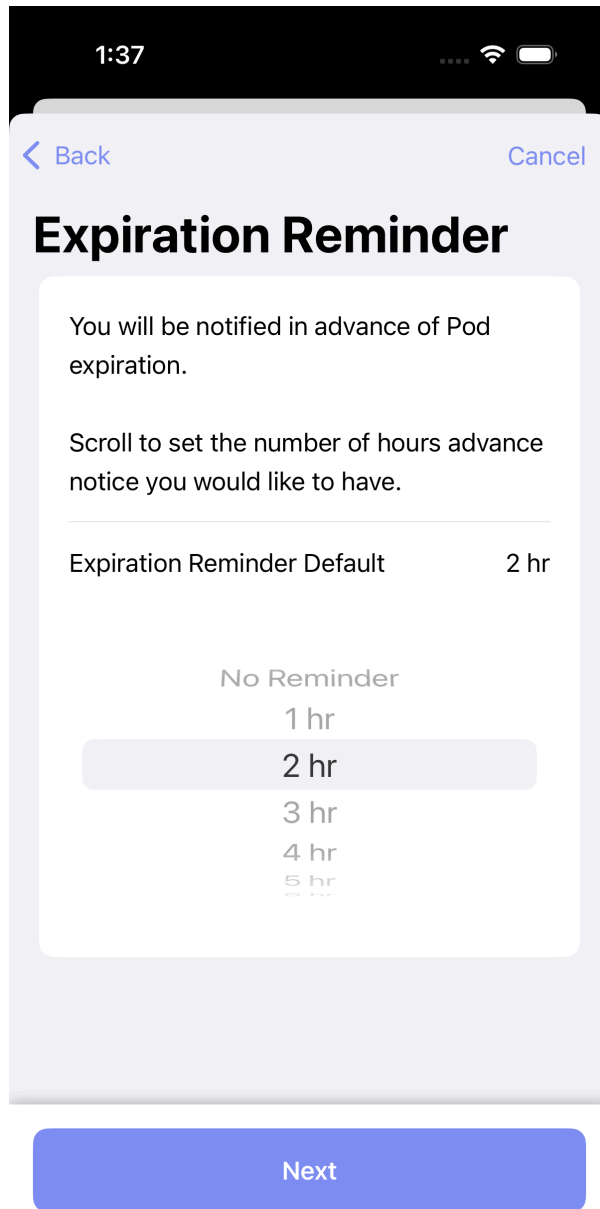

Step 6 Choose the low reservoir reminder that you prefer and tap "Next"

1:38

< Back

Cancel

Low Reservoir

You will be notified when the amount of insulin in the Pod reaches a selected level.

Scroll to set the number of units (1 to 50) at which you would like to be reminded.

Low Reservoir

10 U

6 U

7 U

8 U

9 U

10 U

11 U

12 U

13 U

14 U

Next

Step 7 Choose your insulin type and tap "Continue"

1:38

< Back

Cancel

Insulin Type

Select the type of insulin that you will be using in this pod.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒




Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐

Continue


Step 8 Connect your RileyLink, OrangeLink, or EmaLink and tap "Continue"

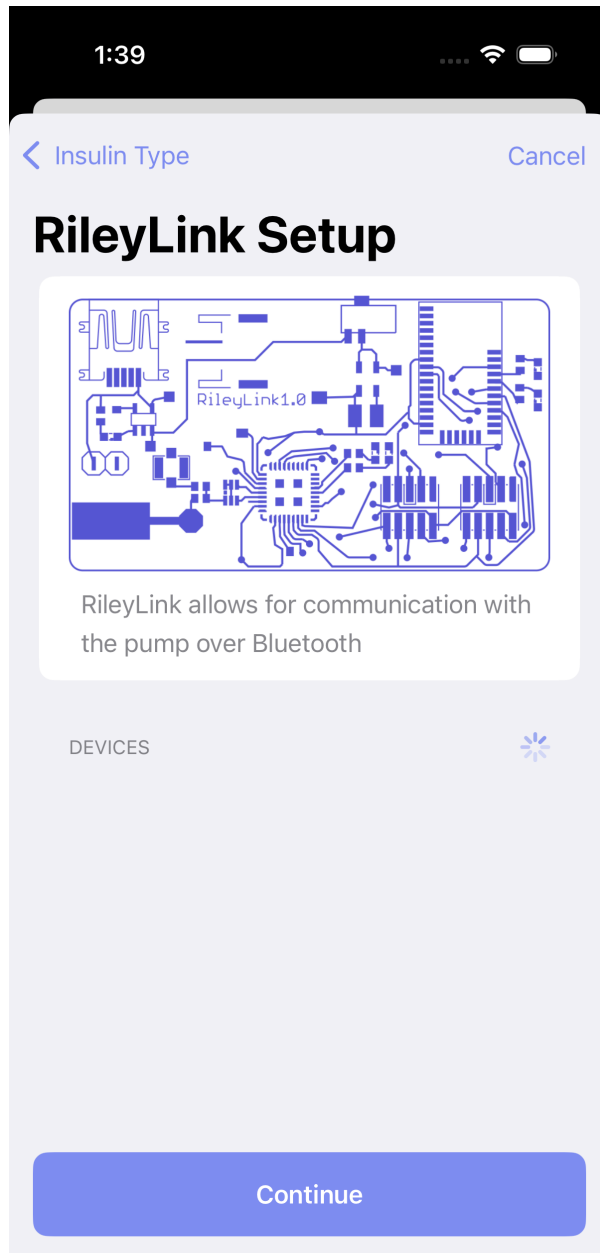


Step 8 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Omnipod Dash

Step 4 Begin the process of configuring your Dash pod by tapping the blue "Continue" button at the bottom

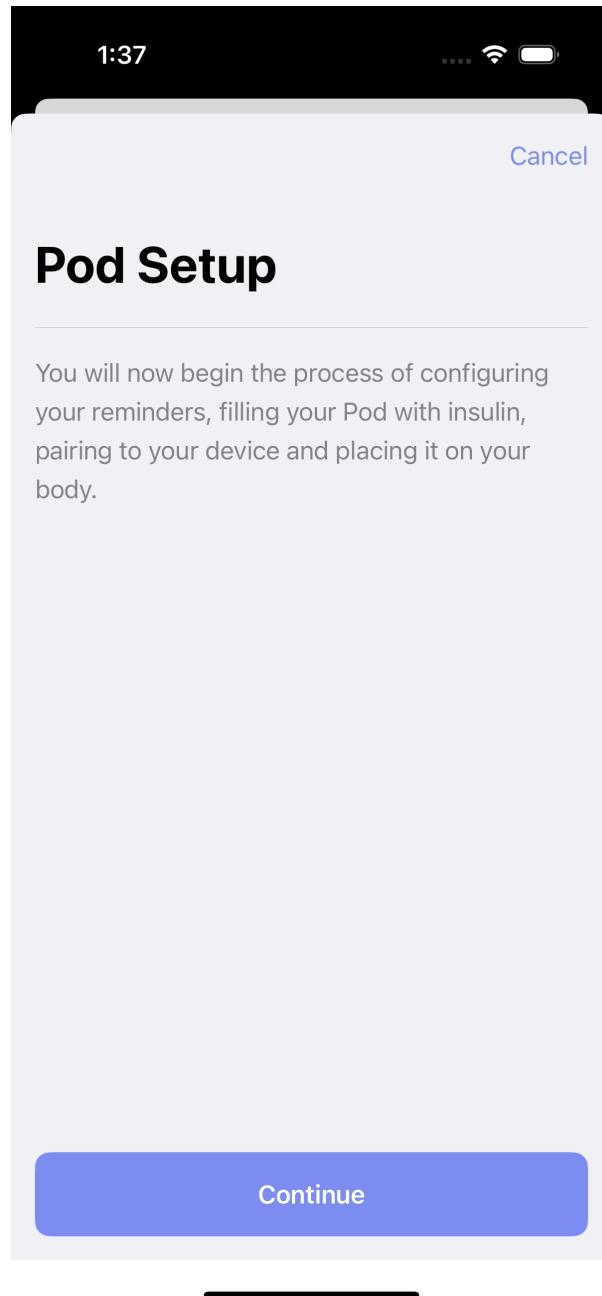

Step 5 Choose the expiration reminder that you prefer and tap "Next"

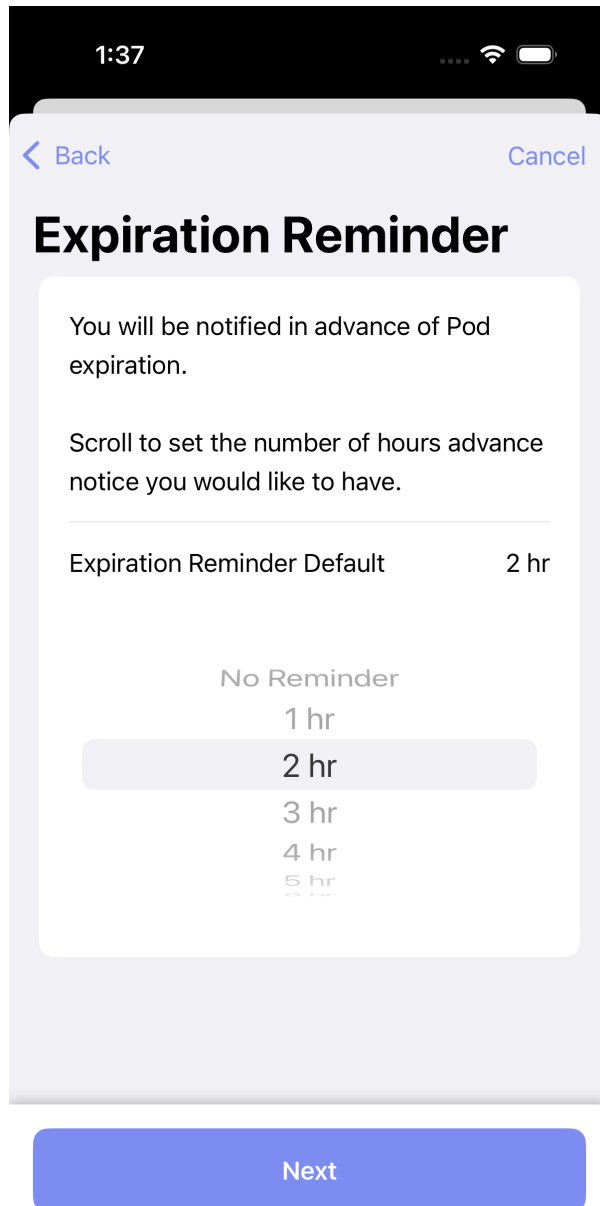

Step 6 Choose the low reservoir reminder that you prefer and tap "Next"

1:38

< Back

Cancel

Low Reservoir

You will be notified when the amount of insulin in the Pod reaches a selected level.

Scroll to set the number of units (1 to 50) at which you would like to be reminded.

Low Reservoir

10 U

6 U

7 U

8 U

9 U

10 U

11 U

12 U

13 U

14 U

Next

Step 7 Choose your insulin type and tap "Continue"

1:38

< Back Cancel

Insulin Type

Select the type of insulin that you will be using in this pod.



Apidra (insulin glulisine)
Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)
Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)
NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒



Fiasp
Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk


☐

Continue


Step 8 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Omnipod 5

Step 4 Begin the process of configuring your Omnipod 5 pod by selecting Omnipod 5 and tap the purple "Continue" button

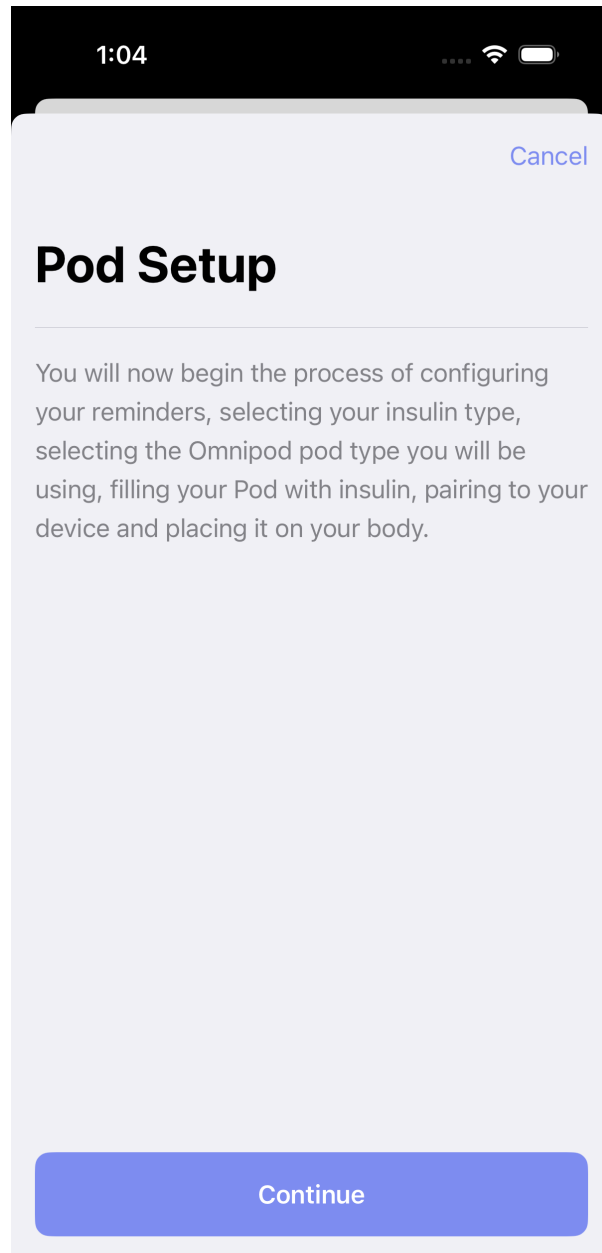

Step 5 Choose the expiration reminder that you prefer and tap "Next"

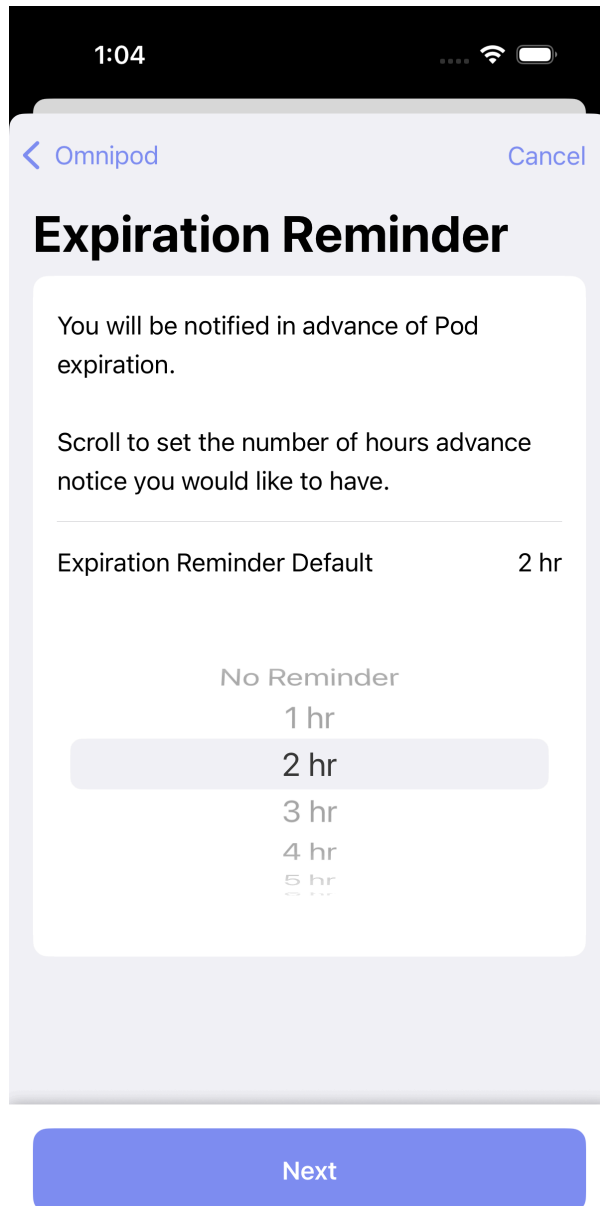

Step 6 Choose the low reservoir reminder that you prefer and tap "Next"

1:04

< Back

Cancel

Low Reservoir Default

You can be optionally notified when the amount of insulin remaining in the reservoir reaches a selected level (1 - 50 units).

Set the default low reservoir alert level to configure when pairing a Pod.

Low Reservoir Default

10

7

8

9

10

11

12

13

Next

Step 7 Choose your insulin type and tap "Continue"

1:04

< Back Cancel

Insulin Type

Select the type of insulin that you will be using in this pod.



Apidra (insulin glulisine)
Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)
Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)
NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒




Fiasp
Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk


☐

Continue


Step 8 Omnipod 5 requires a certificate to be pulled from a database online. This step requires internet access, but only for the very first Omnipod 5 pod. All future pods will not need this step or internet connection.

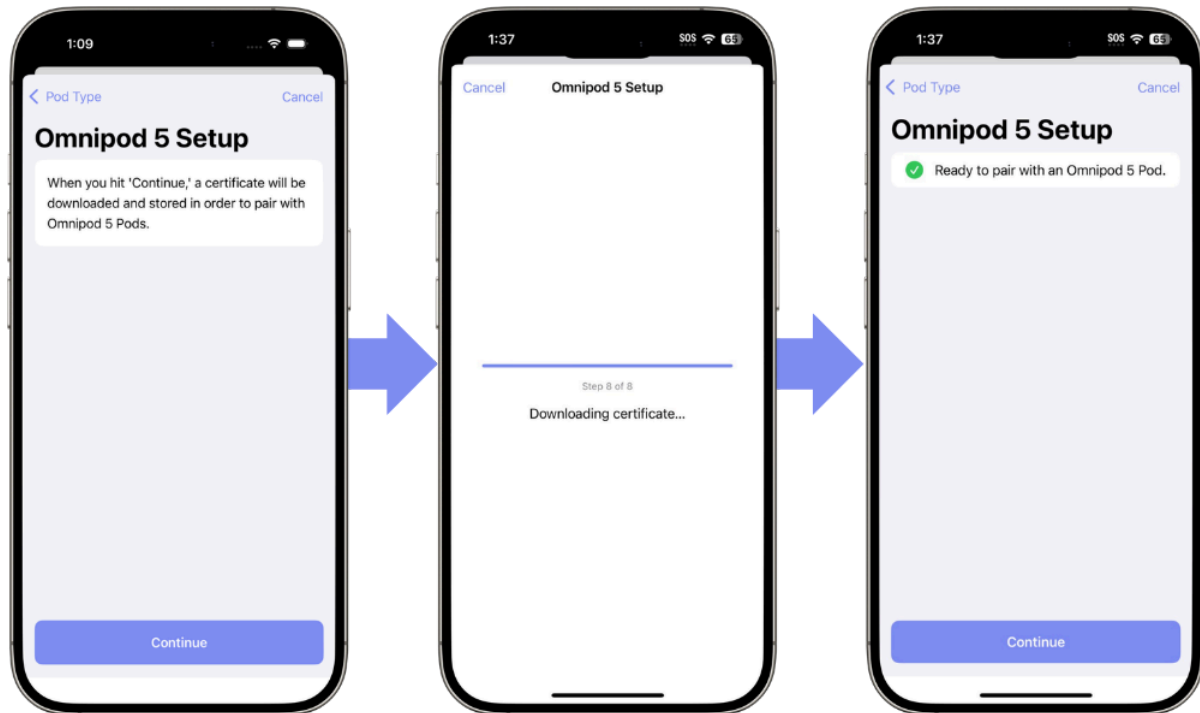


Step 9 Follow the on-screen instructions to pair your Omnipod 5 pod

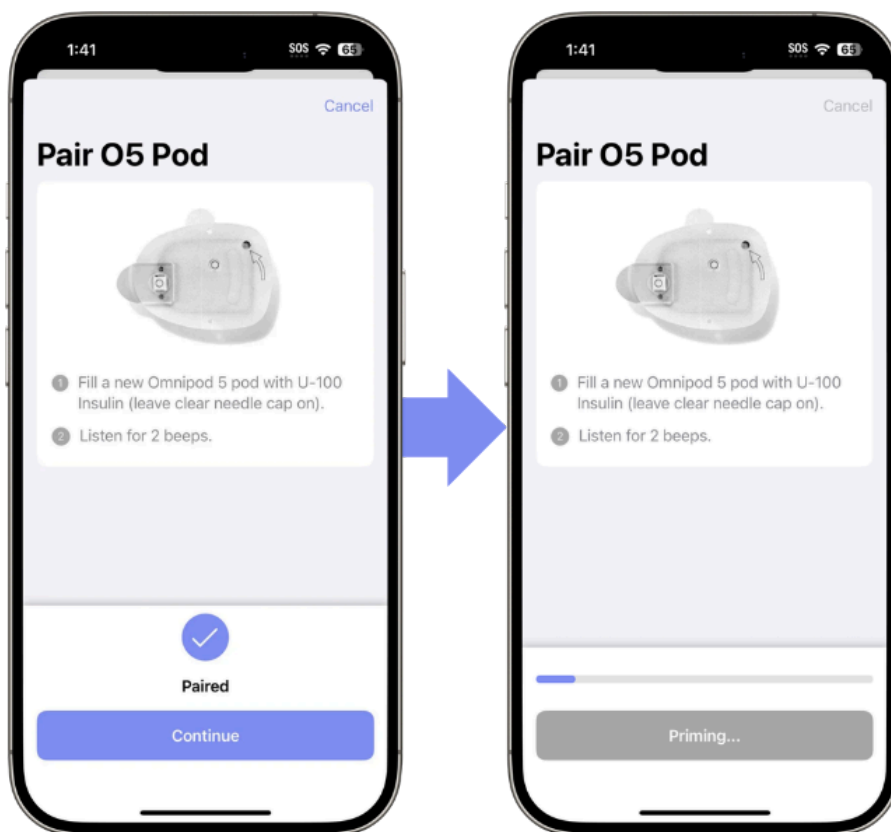

Step 10 Proceed through the steps to insert the cannula

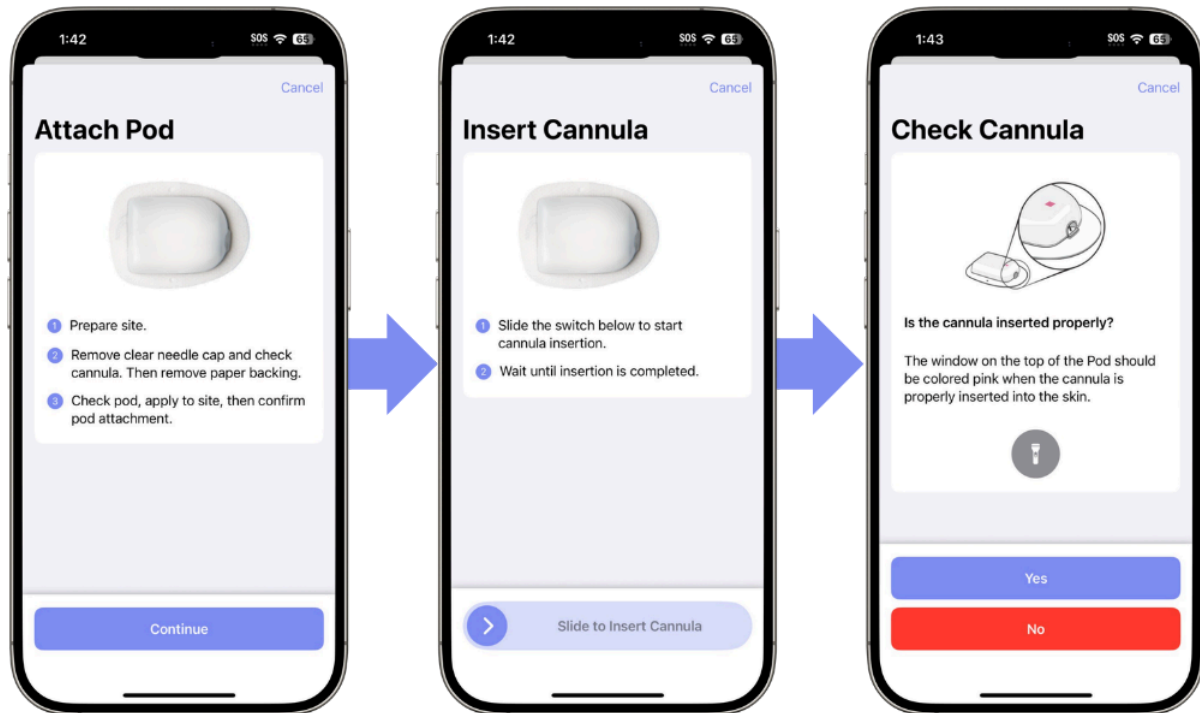

Step 11 Setup complete!

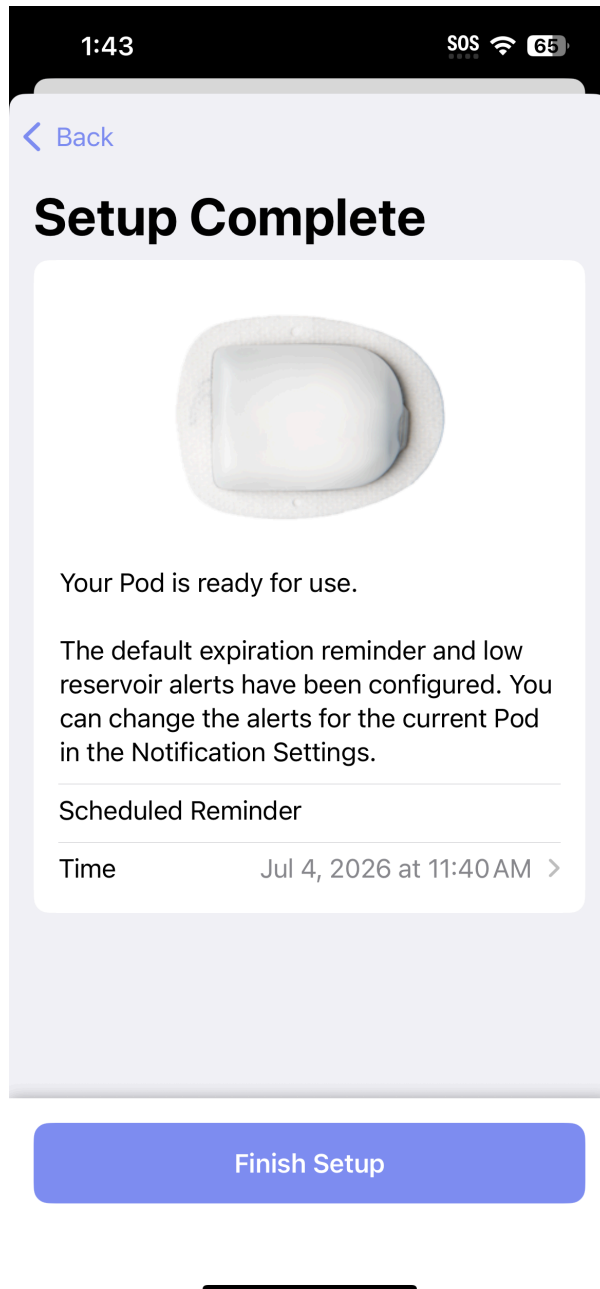


Step 12 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Dana(RS/-i)

Step 3 Select your Dana pump model from the list and tap "Continue"



Step 4 Read the in-app information and tap "Continue" when ready to proceed

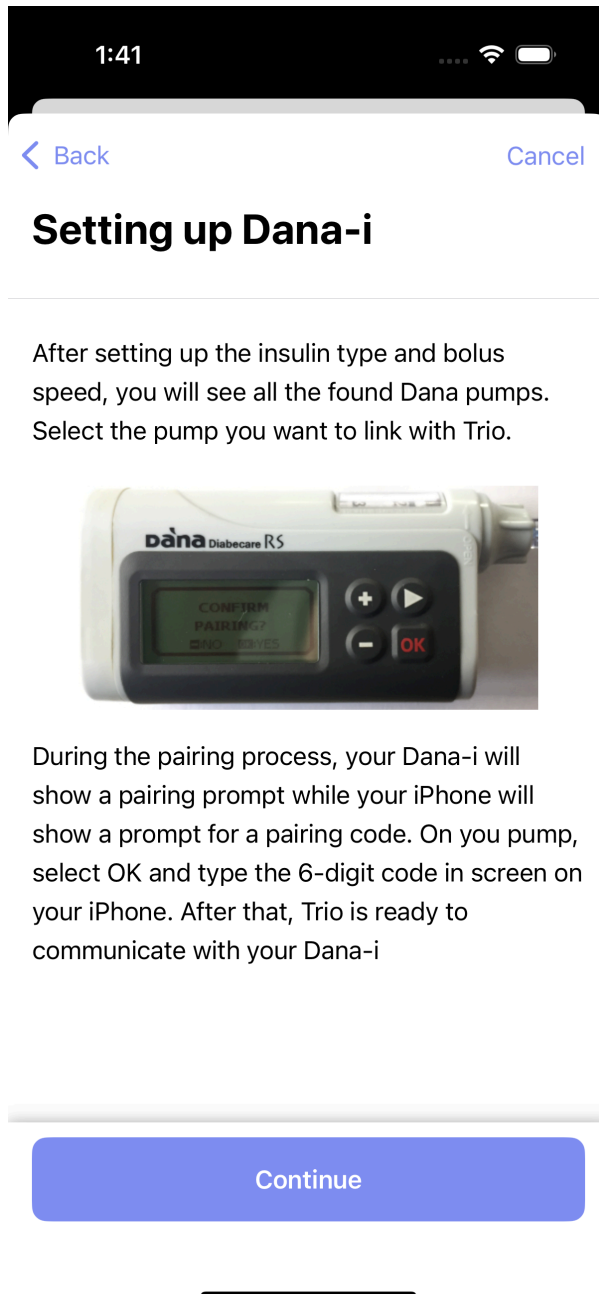

Step 5 Choose your insulin type and tap "Continue"

1:41

< Back

Cancel

Select insulin type

Select the type of insulin that you will be using in this pump.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☒



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☐



Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐




Lyumjev

Lyumjev is a mealtime insulin lispro formulation with the addition of citrate and


☐

Continue

Step 6 Choose your insulin delivery speed and tap "Continue"

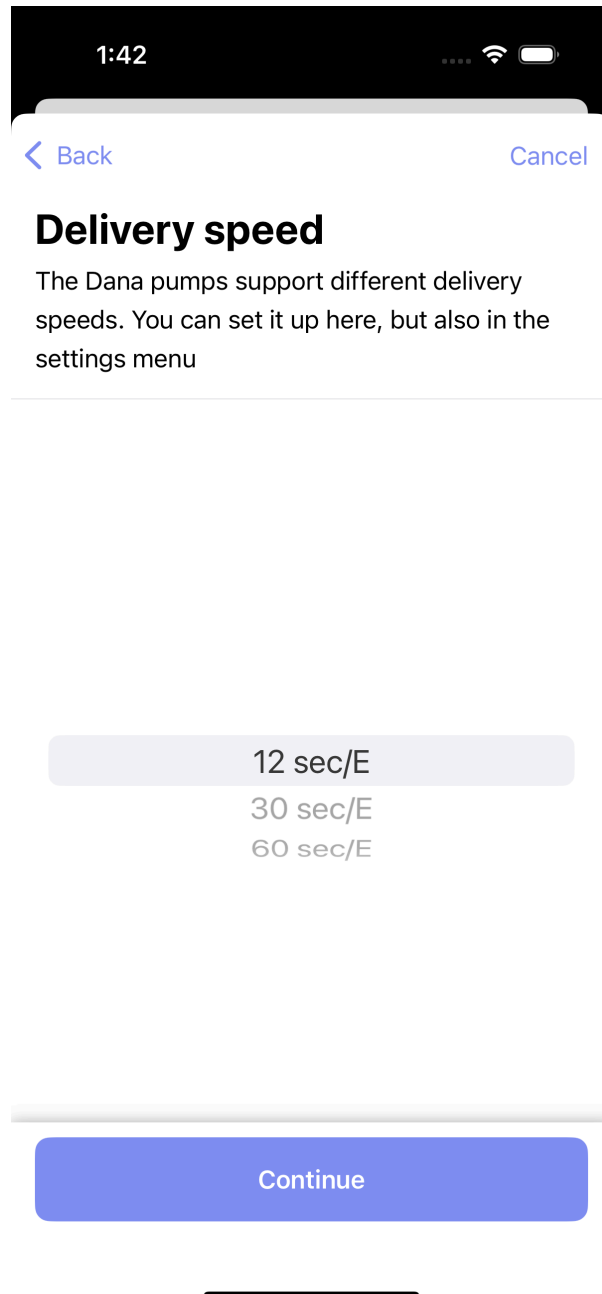

Step 7 Select your pump from the list of found pumps

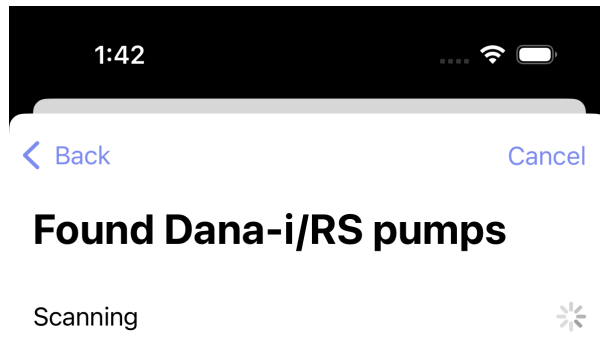

Step 8 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Medtrum Nano

 **Setup Instructions Coming Soon**

Medtronic

Step 3 Choose your insulin type and tap "Continue"

1:43

Insulin Type

Cancel

Select the type of insulin that you will be using in this pump.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒




Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐

Continue

Step 4 Connect your RileyLink, OrangeLink, or EmaLink and tap "Continue"

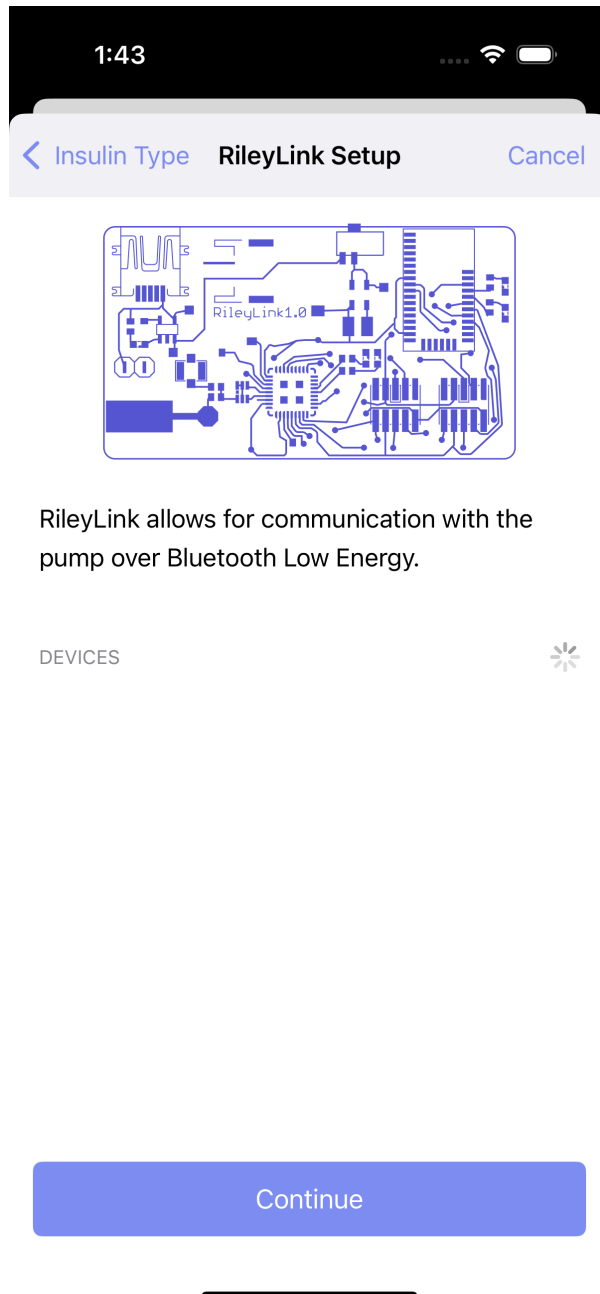

Step 5a Enter your pump region and color

Step 5b Enter your 6-digit pump ID


Then, tap the "Connect" button when it turns blue at the bottom.

1:44

< Back


Pump Setup

Cancel

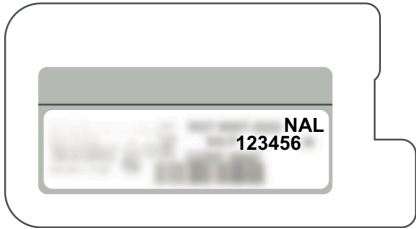

REGION AND COLOR

Enter the pump region

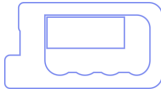

The pump region and color are denoted as the last 3 letters of the the model number (labeled as REF).

PUMP ID

Enter the 6-digit pump ID

The pump ID is the 6-digit numerical portion of the serial number (labeled as SN or S/N).

Connect

Step 6 Continue to [Connect CGM](#) OR return to [New User Setup](#)


Pump Simulator

Warning

If using a pump simulator, it is important to understand:

- You will only experience the user interface of Trio.
- Using a pump simulator does not indicate how the app will perform, nor will it give accurate guidance or suggestions for insulin dosing.
- Trio may not run consistently in the background when using a pump simulator.
- If you use a pump simulator and switch to a live pump, you must delete all insulin dosing history by uninstalling the Trio app and reinstalling it.
- **Only use a pump simulator if you understand the conditions above.**

Step 3 Choose from the variety of simulator pump options and tap "Done" in the upper right corner to activate the pump simulator

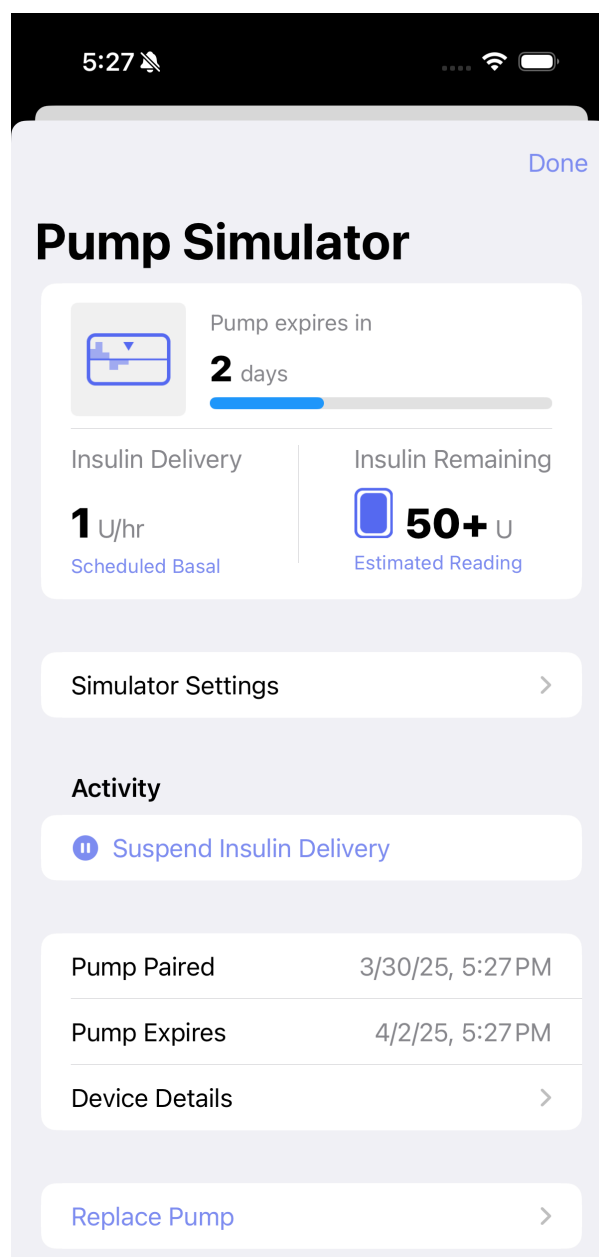

Step 4 Continue to [Connect CGM](#) *OR* return to [New User Setup](#)
