



TrioDocs

Version: 0.4.0

Date: December 29, 2025

Download the latest version on:
<https://triodocs.org>


Insulin Pump

Insulin Pump

This section allows you to choose your insulin pump. Most options are self-explanatory. For more information on compatible insulin pumps, please see the following [link](#)

Step 1: Add Pump

The first step in setting up your insulin pump on Trio is to tap the "Add Pump" button in the Devices Menu

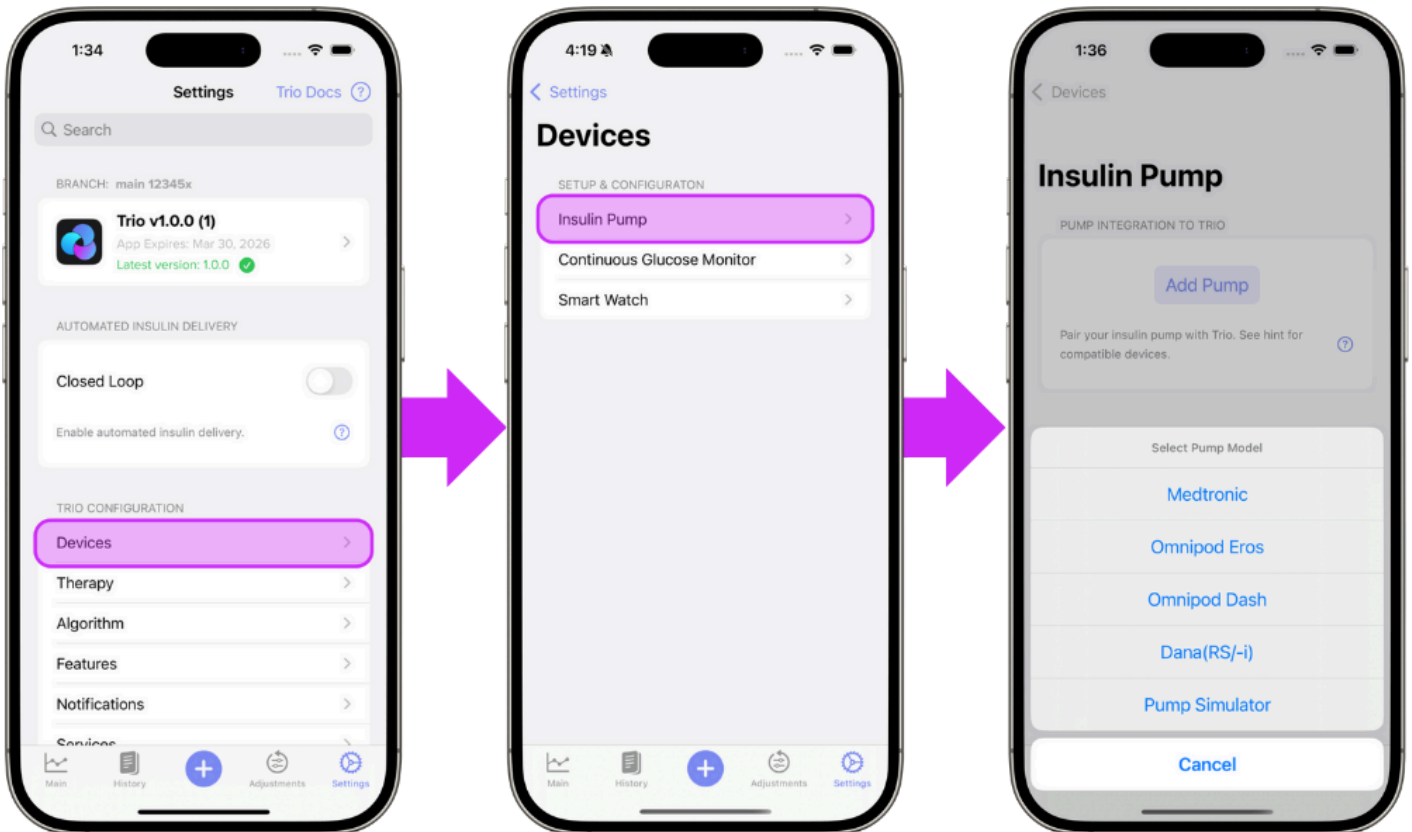

Step 2: Select Your Pump

Select your pump from the in-app menu and from the options below for step-by-step instructions. The links below will guide you through the connection instructions for your specific pump.

- [Medtronic](#)
- [Omnipod Eros](#)
- [Omnipod Dash](#)
- [Dana\(RS/-i\)](#)
- [Pump Simulator](#)

Medtronic

Step 3 Choose your insulin type and tap "Continue"

1:43

Insulin Type

Cancel

Select the type of insulin that you will be using in this pump.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒




Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐

Continue


Step 4 Connect your RileyLink, OrangeLink, or EmaLink and tap "Continue"

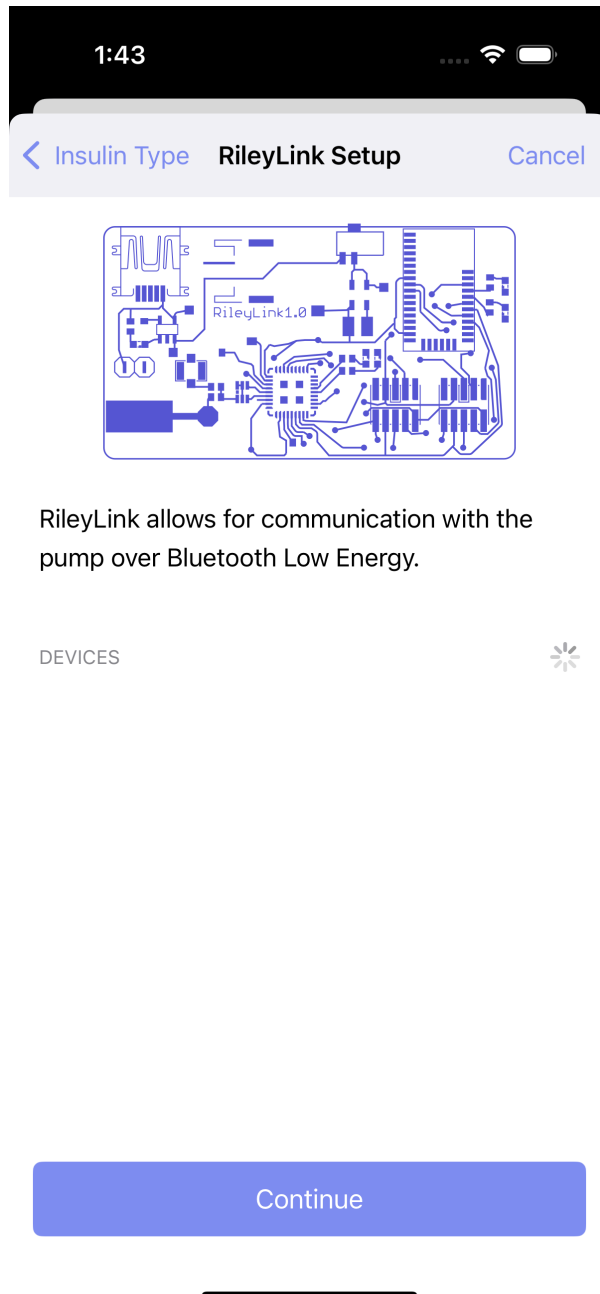

Step 5a Enter your pump region and color

Step 5b Enter your 6-digit pump ID


Then, tap the "Connect" button when it turns blue at the bottom.

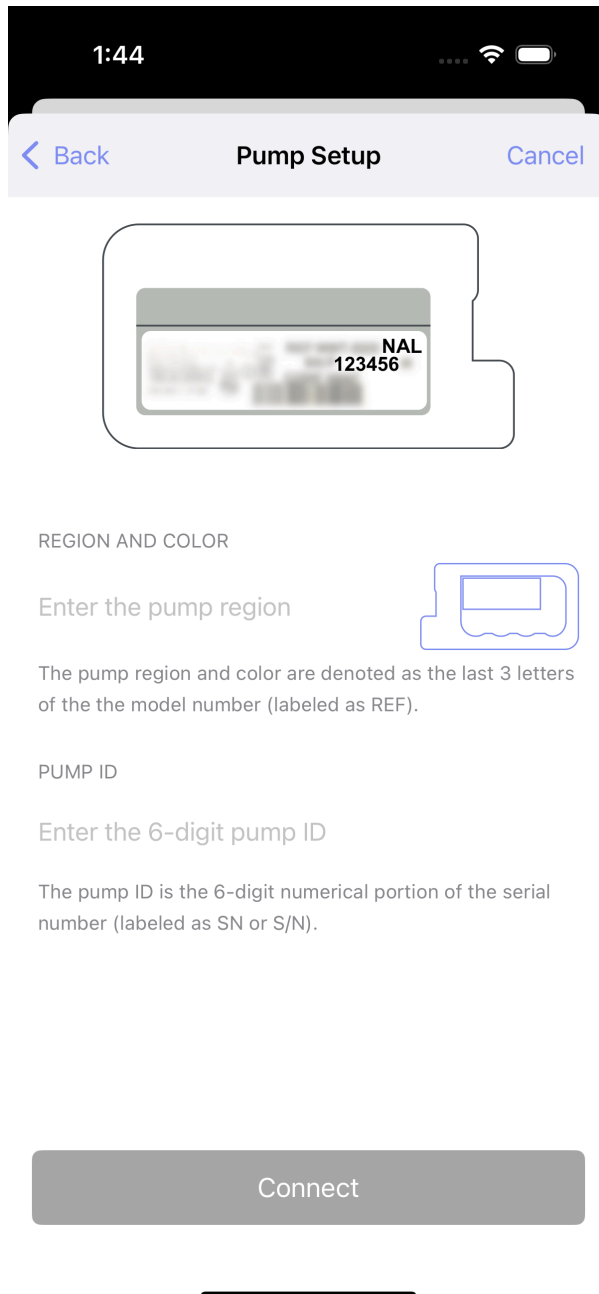


Step 6 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Omnipod Eros

Step 3 Begin the process of configuring your Eros pod by tapping the blue "Continue" button at the bottom

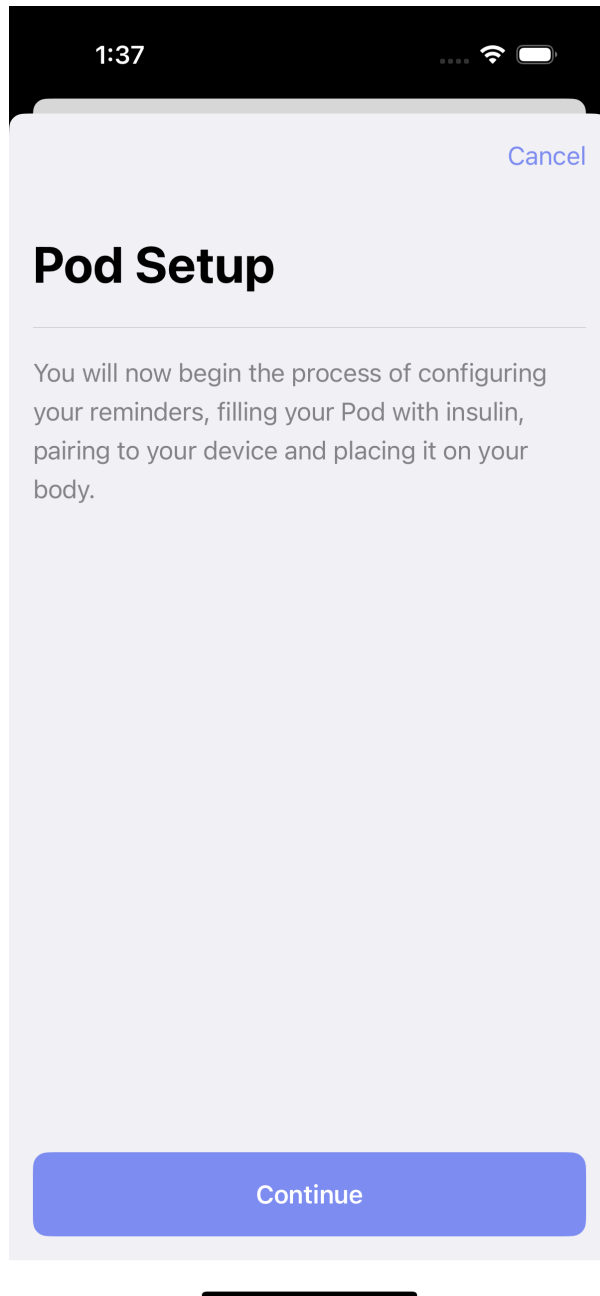

Step 4 Choose the expiration reminder that you prefer and tap "Next"

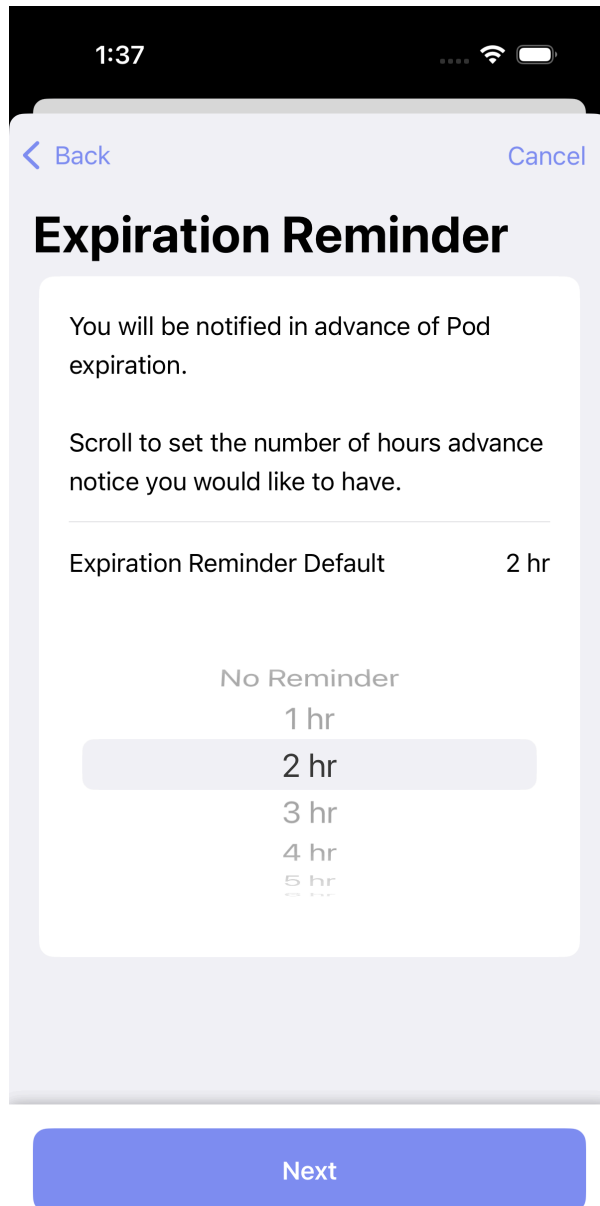

Step 5 Choose the low reservoir reminder that you prefer and tap "Next"

1:38

< Back

Cancel

Low Reservoir

You will be notified when the amount of insulin in the Pod reaches a selected level.

Scroll to set the number of units (1 to 50) at which you would like to be reminded.

Low Reservoir

10 U

6 U

7 U

8 U

9 U

10 U

11 U

12 U

13 U

14 U

Next

Step 6 Choose your insulin type and tap "Continue"

1:38

< Back

Cancel

Insulin Type

Select the type of insulin that you will be using in this pod.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒




Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐

Continue


Step 7 Connect your RileyLink, OrangeLink, or EmaLink and tap "Continue"

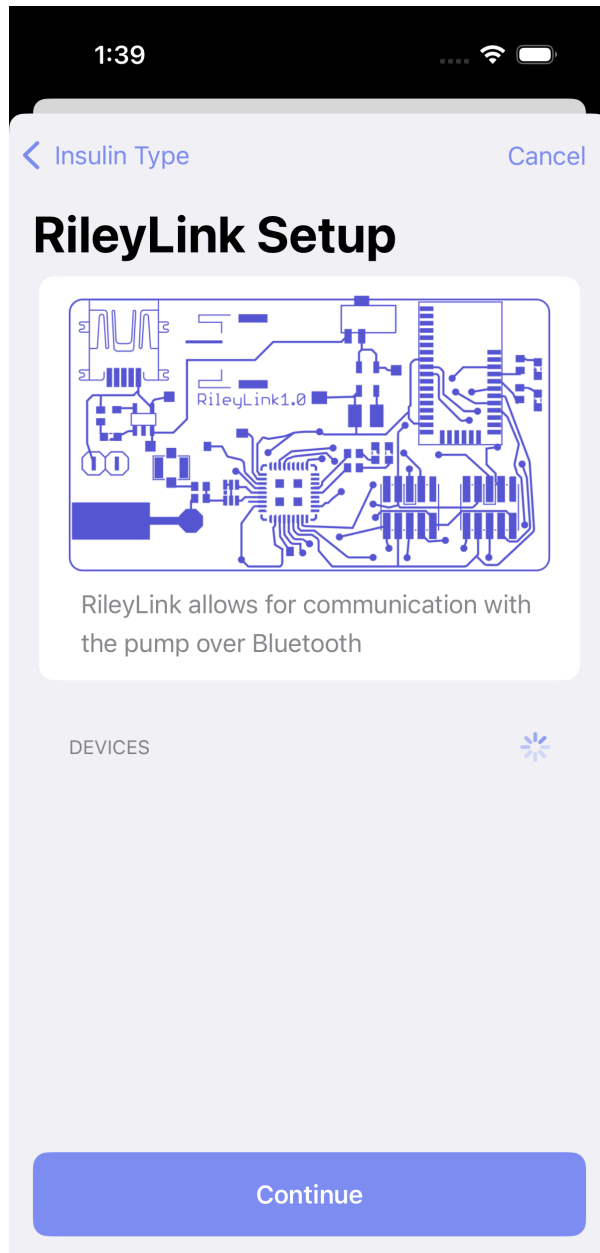


Step 8 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Omnipod Dash

Step 3 Begin the process of configuring your Dash pod by tapping the blue "Continue" button at the bottom

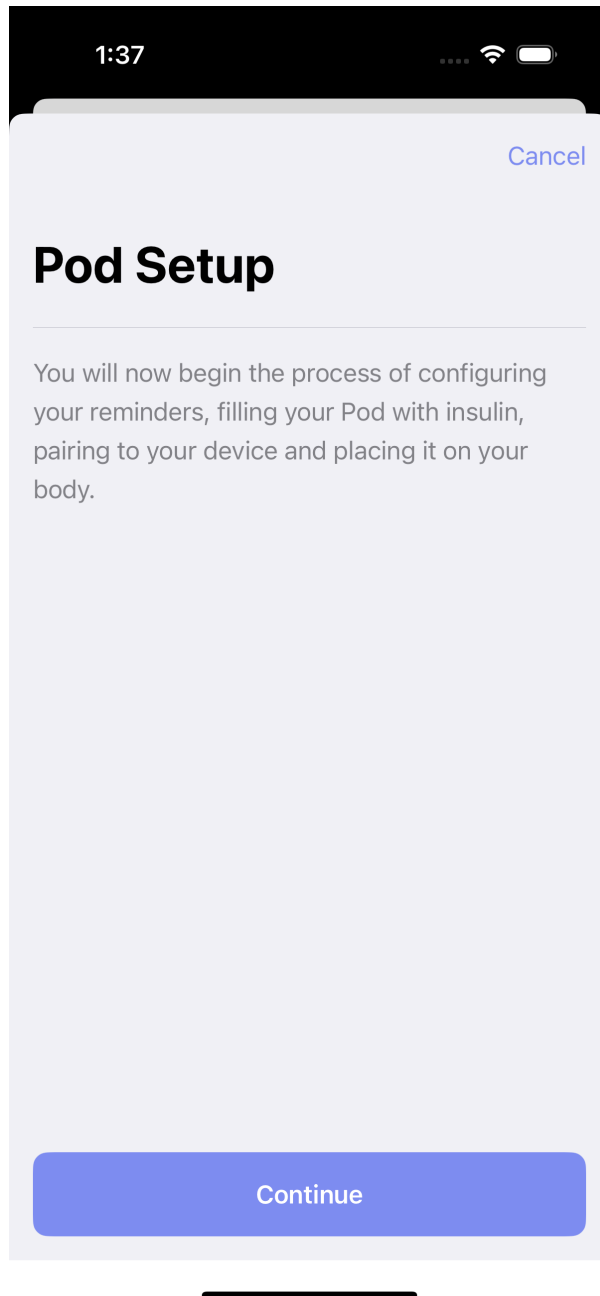

Step 4 Choose the expiration reminder that you prefer and tap "Next"

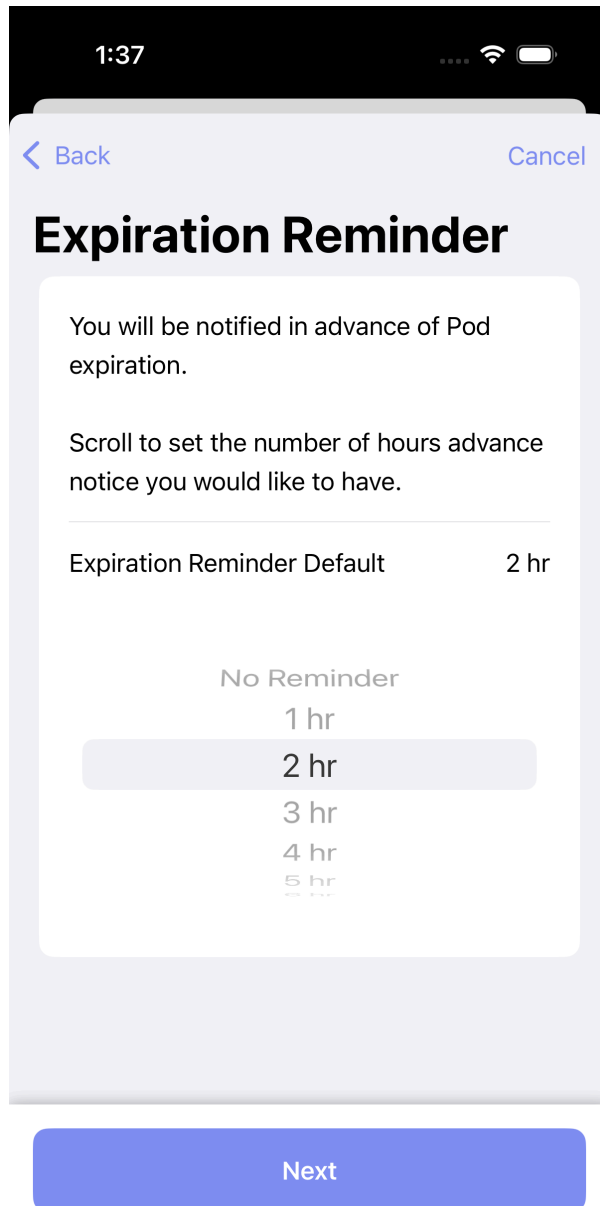

Step 5 Choose the low reservoir reminder that you prefer and tap "Next"

1:38

< Back

Cancel

Low Reservoir

You will be notified when the amount of insulin in the Pod reaches a selected level.

Scroll to set the number of units (1 to 50) at which you would like to be reminded.

Low Reservoir

10 U

6 U

7 U

8 U

9 U

10 U

11 U

12 U

13 U

14 U

Next

Step 6 Choose your insulin type and tap "Continue"

1:38

< Back Cancel

Insulin Type

Select the type of insulin that you will be using in this pod.



Apidra (insulin glulisine)
Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☐



Humalog (insulin lispro)
Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)
NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☒



Fiasp
Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐

Continue


Step 7 Continue to [Connect CGM](#) OR return to [New User Setup](#)

Dana(RS/-i)

Step 3 Select your Dana pump model from the list and tap "Continue"



Step 4 Read the in-app information and tap "Continue" when ready to proceed

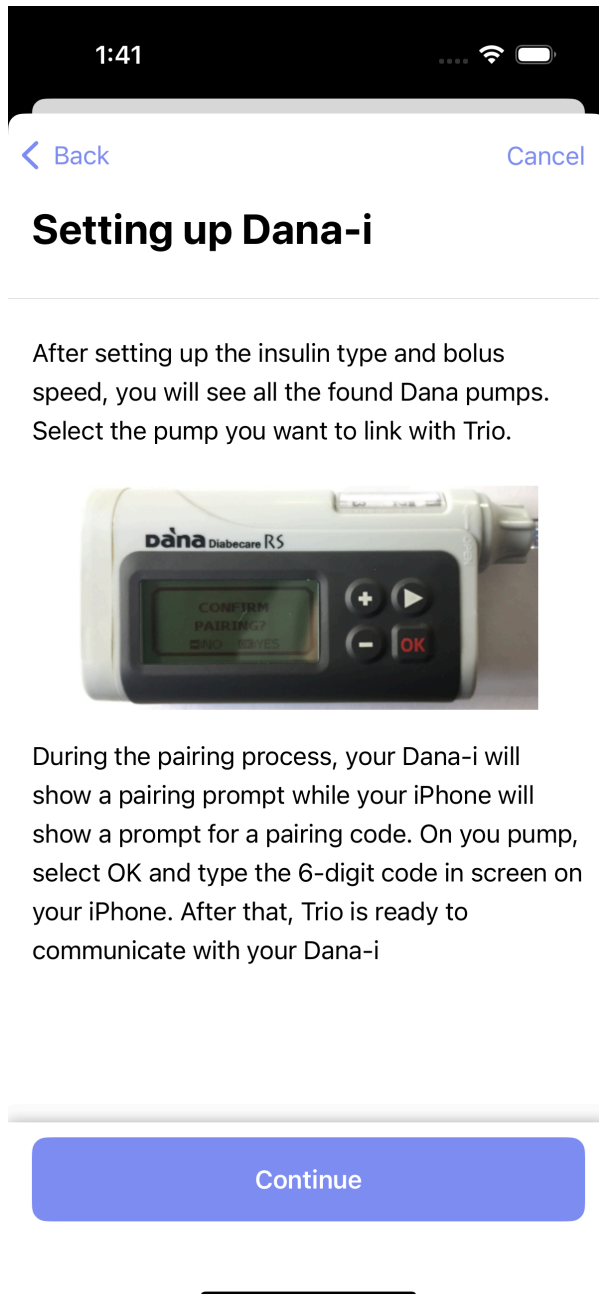

Step 5 Choose your insulin type and tap "Continue"

1:41

< Back

Cancel

Select insulin type

Select the type of insulin that you will be using in this pump.



Apidra (insulin glulisine)

Apidra (insulin glulisine) is a rapid-acting insulin made by Sanofi-aventis

☒



Humalog (insulin lispro)

Humalog (insulin lispro) is a rapid-acting insulin made by Eli Lilly

☐



Novolog (insulin aspart)

NovoLog (insulin aspart) is a rapid-acting insulin made by Novo Nordisk

☐



Fiasp

Fiasp is a mealtime insulin aspart formulation with the addition of nicotinamide (vitamin B3) made by Novo Nordisk

☐




Lyumjev

Lyumjev is a mealtime insulin lispro formulation with the addition of citrate and


☐

Continue

Step 6 Choose your insulin delivery speed and tap "Continue"

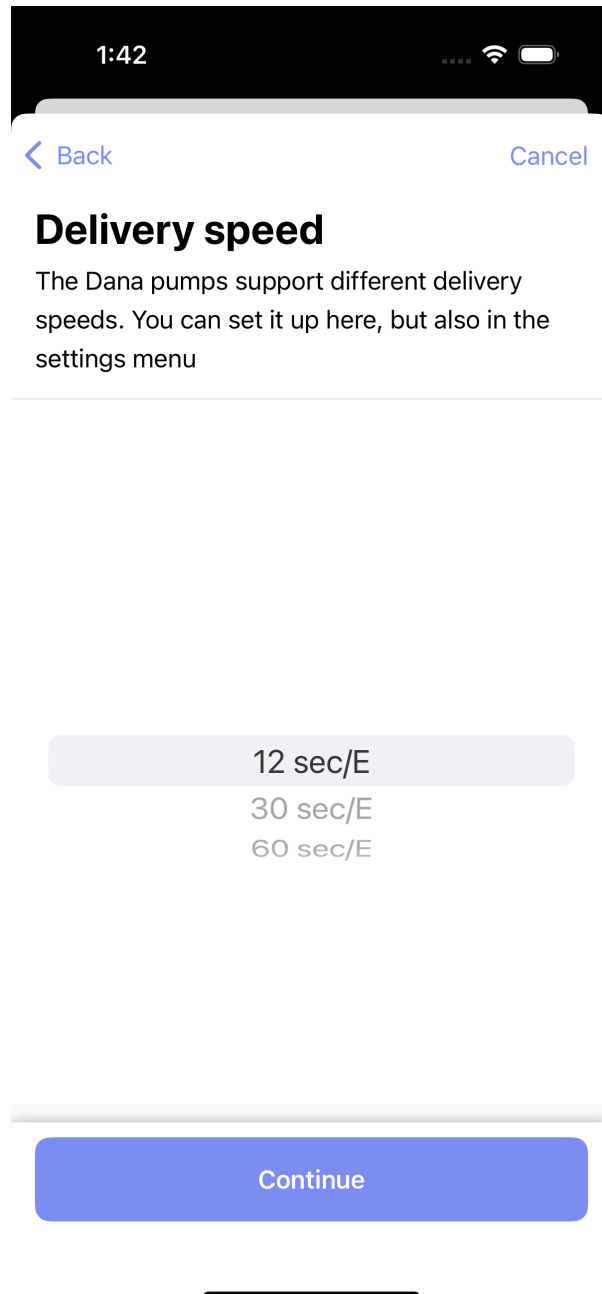

Step 7 Select your pump from the list of found pumps

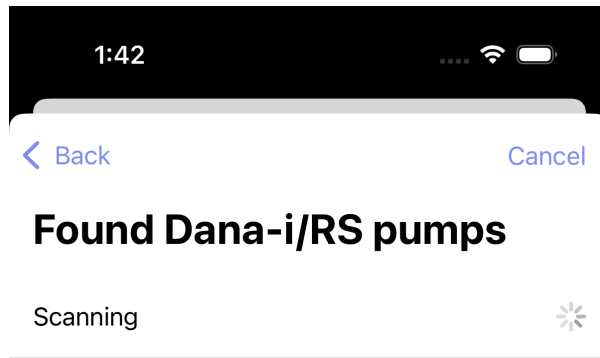

Step 8 Continue to [Connect CGM](#) OR return to [New User Setup](#)


Pump Simulator

 Warning

If using a pump simulator, it is important to understand:

- You will only experience the user interface of Trio.
- Using a pump simulator does not indicate how the app will perform, nor will it give accurate guidance or suggestions for insulin dosing.
- Trio may not run consistently in the background when using a pump simulator.
- If you use a pump simulator and switch to a live pump, you must delete all insulin dosing history by uninstalling the Trio app and reinstalling it.
- **Only use a pump simulator if you understand the conditions above.**

Step 3 Choose from the variety of simulator pump options and tap "Done" in the upper right corner to activate the pump simulator

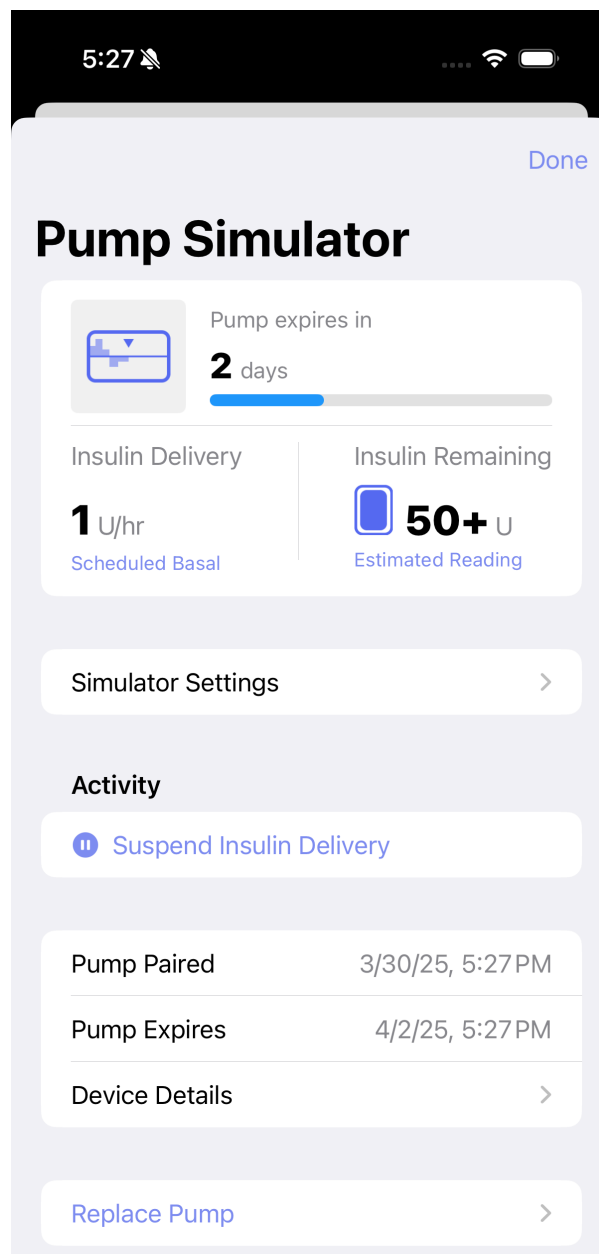

Step 4 Continue to [Connect CGM](#) OR return to [New User Setup](#)
