



TrioDocs

Version: 0.4.0

Date: May 17, 2025

Download the latest version on:  
<https://triodocs.org>

# Glucose Targets

# Glucose Targets

Trio will target this value when calculating insulin needs. It should be set to the blood glucose you would like to reach when corrections are made. A recommended value is between 90-110mg/dL or 5-6mmol/L

## Tip

You can set different targets at different times.


Some users only have one target.

Others have a lower target during the day and a higher target during the night to avoid lows.

## How to Enter Your Glucose Targets Into Trio


### Step 1

Enter the Glucose Targets screen

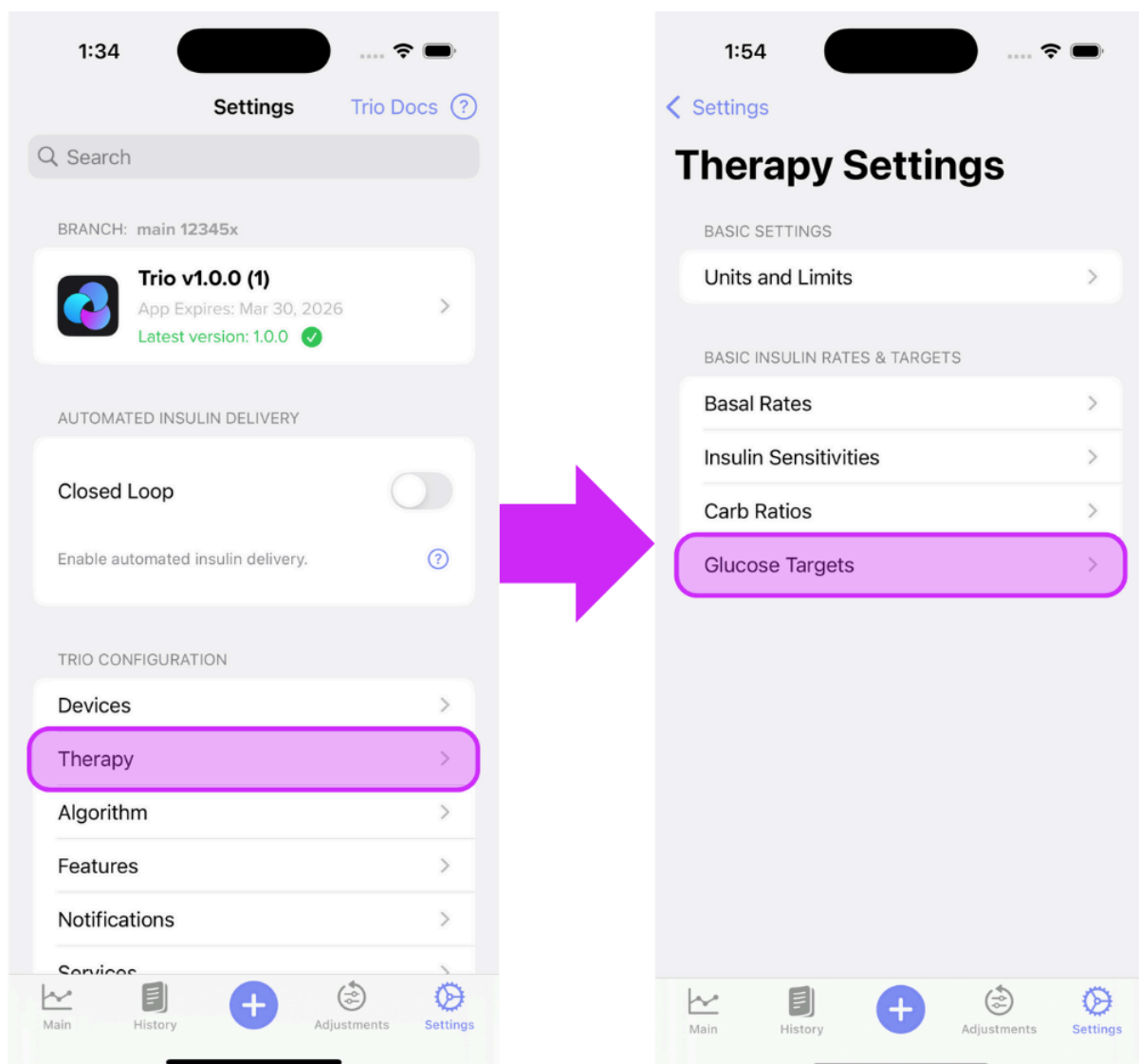



## Step 2

Tap the "Add Rate +" on the top right of the screen until you have the number of targets you require. Then, edit each target by tapping the arrow to the right of the Glucose Target.

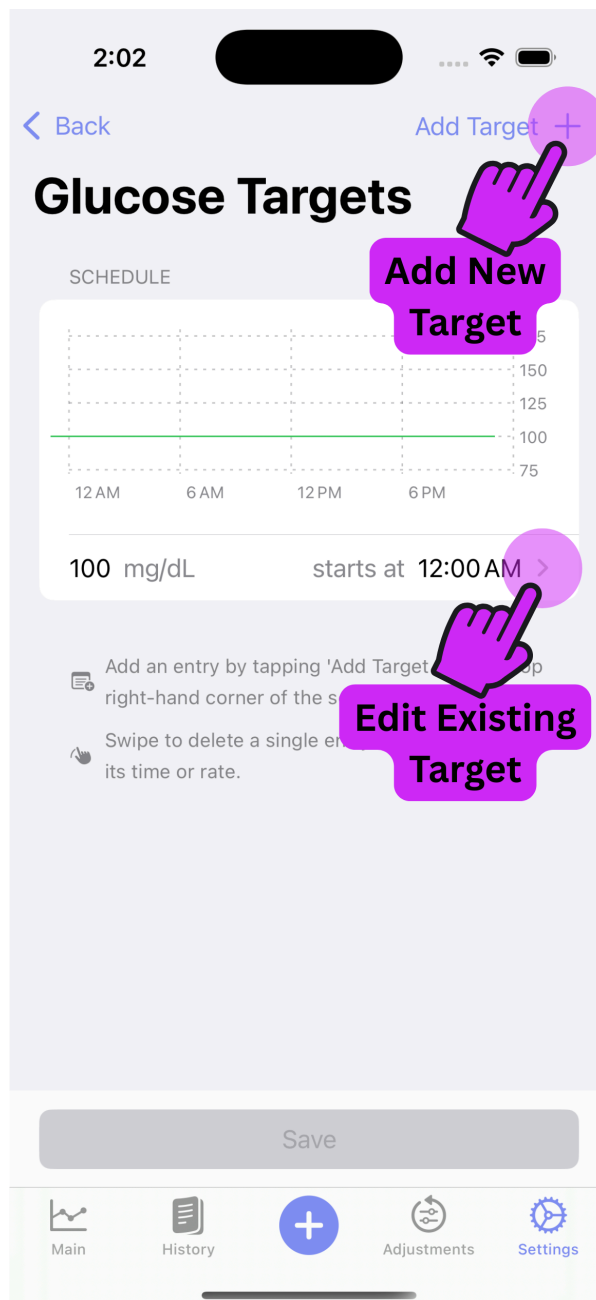


## Step 3

Adjust the target glucose value

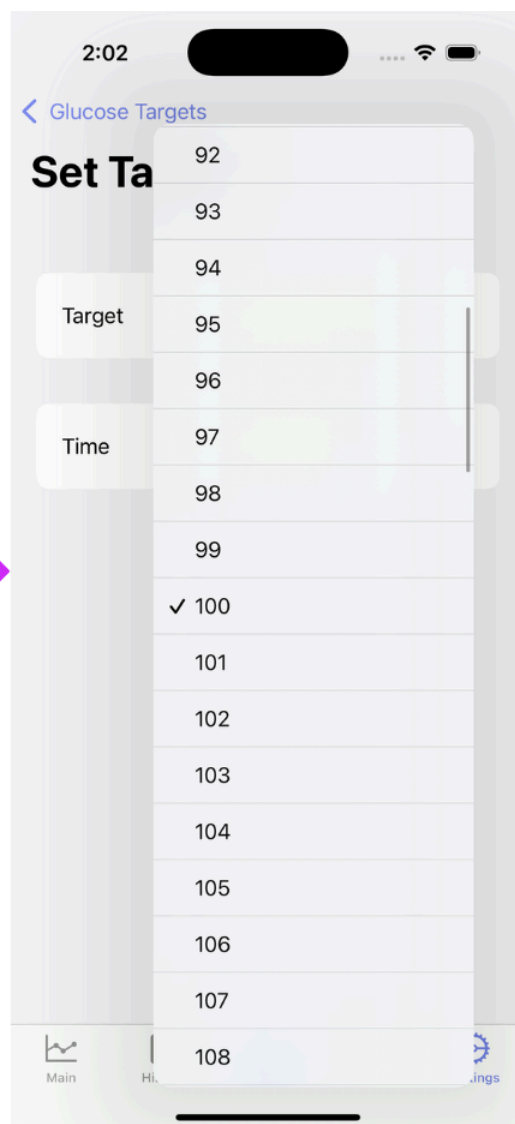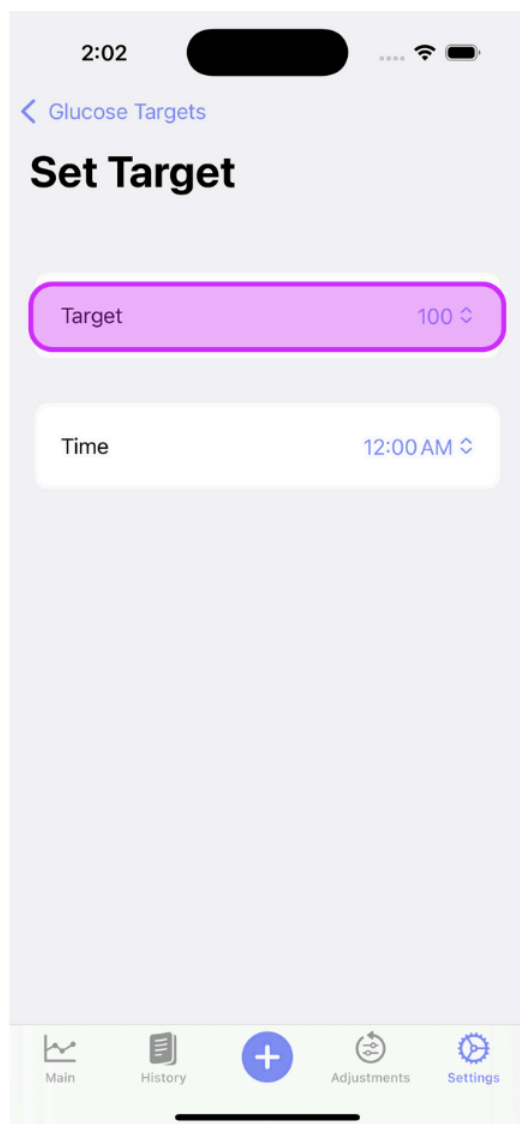

## Step 4

Adjust the time

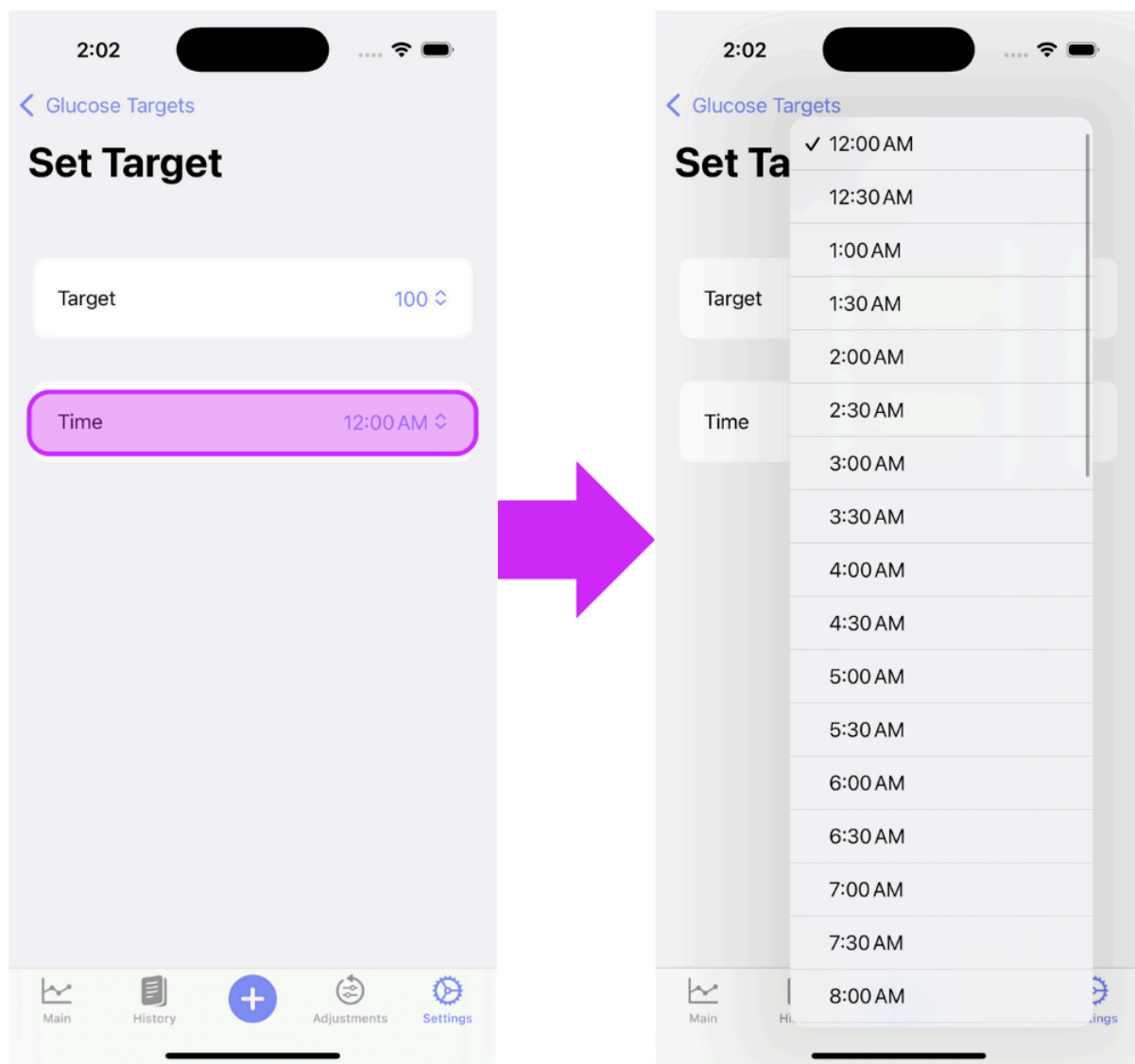


## Step 5

Repeat Steps 2, 3, and 4 until all CRs are set


## Delete a Glucose Target Entry

Should you need to delete an CR entry, just swipe left on the rate you need to remove.

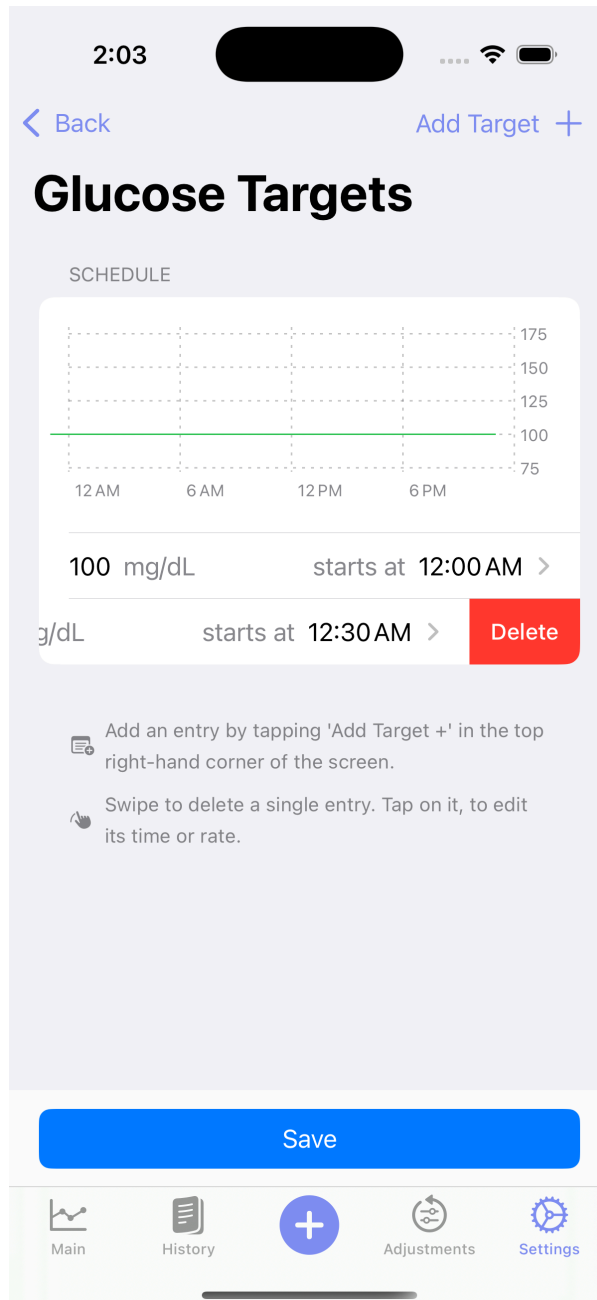

## Step 6 **IMPORTANT**

Save your changes!

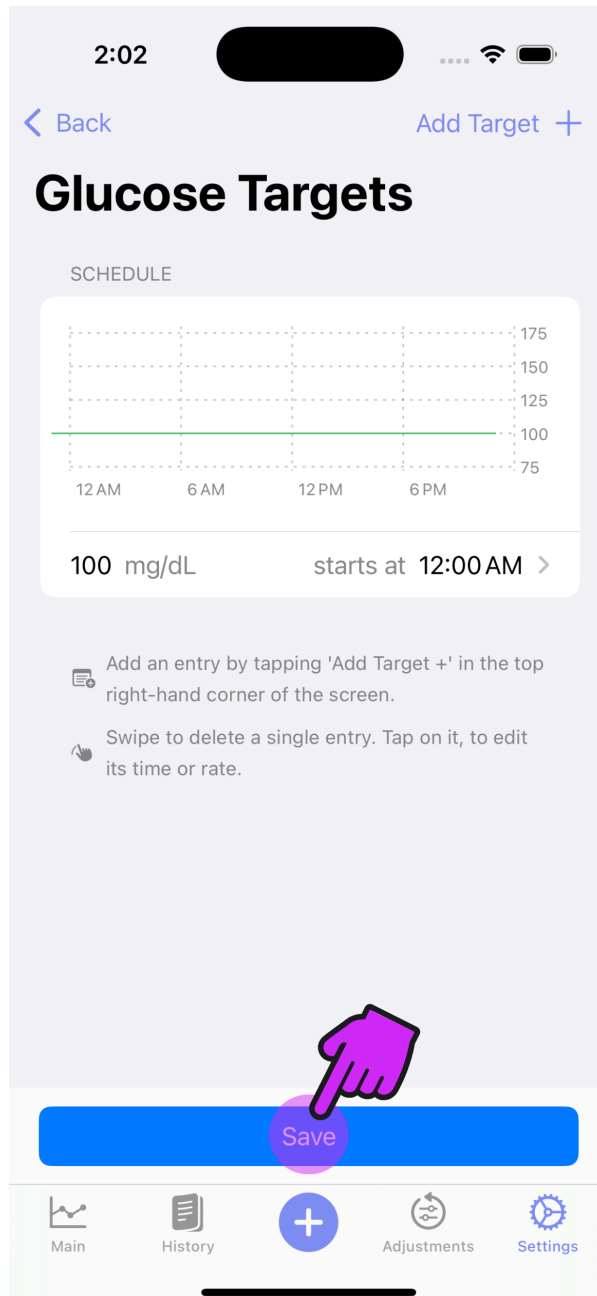

## Step 7

Proceed to [Basal Rates](#) or return to [New User Setup](#)