



TrioDocs

Version: 0.4.0

Date: April 22, 2026

Download the latest version on:  
<https://triodocs.org>

## Install Trio from TestFlight

## General Installation Information

**The *Trio* app must be built at least every 90 days when using a browser to build.**

The automatic build actions should occur every Sunday.

It is recommended you manually install the new build using *TestFlight* at your convenience.

After you [Build the Trio App](#) with a browser and you get the email that it is available in *TestFlight*, you are ready to install on as many phones as you and your family members need.


- If you later need to add an adult family member to your list, refer to [TestFlight Overview](#).
- Children (under 13 in US, varies by country) cannot use *TestFlight* with their ID. When you use [TestFlight for a Child](#), you will need to use your ID on their phone (not the whole phone - just the Media & Purchase portion), so send the *TestFlight* invitation to the email associated with your ID.

## Install *TestFlight*

If you already have the *TestFlight* app installed on your phone, skip ahead to [Install App with TestFlight](#).

To install *TestFlight*, refer to the GIF below:

- On the phone, open the App Store and Search for *TestFlight*
- Install or Download to that phone *TestFlight*
  - Hint: On child's phone, do this while logged in as yourself for Media & Purchase
  - Logging in as an adult is explained in [TestFlight for a Child](#)

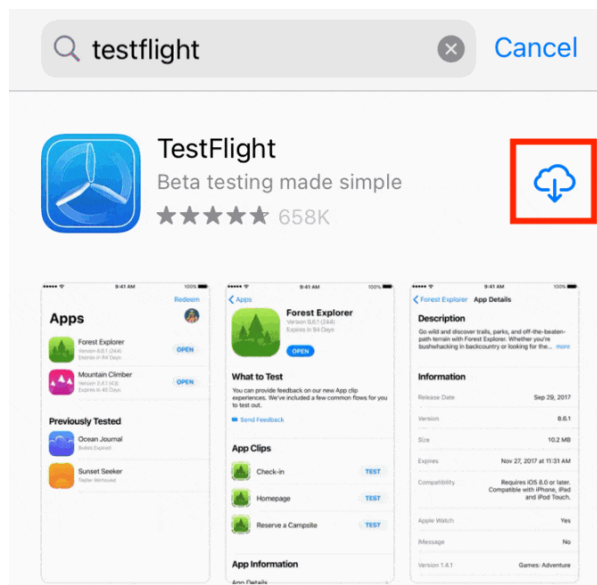

## Install App with *TestFlight*

Once you get an email that your app is available to test on iOS and watchOS, you can install the app on your phone from *TestFlight*.

- Note this can be half-hour to an hour after the build displays the green check mark on your browser
- Once the app is on your phone, you can choose to install the watch version using your phone *Watch* app
- If you did not already add your email to your *TestFlight* group for your app, go and do it now using these [instructions](#) ↗

The first time you use *TestFlight* on any phone associated with a given email, you must *Redeem* the code sent to that email inviting you to test the app. The GIF below is for someone who has never used *TestFlight*; the example is for *Loop* but *Trio* acts the same.

- Initial screen indicates there are no Apps available to test, tap on Redeem
- Enter your code and tap redeem to enter it
- Click on OK to acknowledge
- Click on Install

1:01



Redeem

## No Apps Available to Test

You currently have no apps to test.  
Tap on the public link or email  
invitation again to join the beta, or  
contact the developer for an  
invitation code.

---


If you already have the *Trio* app on the phone, you'll see the warning about possible loss of data. Don't worry, all your settings remain. Go ahead with the installation.

### Subsequent Times on Phone

- Open the *TestFlight* app and find the name you used for your *Trio* app in the [Configure Trio App](#) step
- If you have previously used *TestFlight* on this phone and don't see the latest build, sometimes quitting and restarting *TestFlight* will bring that new build into the list of available apps
- Most people will just tap install to get the most recent build
  - If you have more than one version number available in *TestFlight* you can choose which version to install by tapping on the [Previous Builds](#) row and then selecting the desired version
- Tap on Install
  - If you already have the *Trio* app installed on this phone, you will be warned that the app already exists on your phone and that you might lose data

- Click Install again (your pump connection and all your data will be fine)
- Choose Open
- Make sure the *Trio* app is operating as expected

This GIF shows an example of installing *Loop* from *TestFlight*.

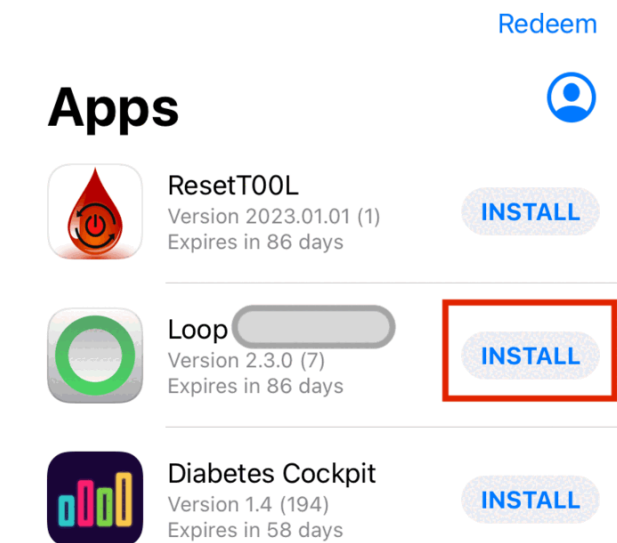

## Automatic Update, Build, Install

If you followed all the instructions for setting up Build with Browser, your GitHub account will automatically take the following actions:

- Update the version of your fork within a week of a new release
  - Automatically create a new build and upload it to *TestFlight*
  - This is only for the `default` branch, typically `main`
- Build the app at least once a month and upload it to *TestFlight*

You need to take this manual action for your *Trio* app in *TestFlight*, at every 90 days, but better to do it monthly.

- Install each new build from *TestFlight* on the phone as soon as it is detected
- That's fine for some apps, like *LoopFollow*, but you do not want an app that controls your insulin delivery to install when you are not paying attention

## Recommendation

Recommended settings:

- Allow automatic update of your fork
- Allow automatic build and upload to *TestFlight*
- **Disable automatic installation** on phone from *TestFlight*


If you are running the development code, you may prefer to turn off the automatic update, but keep the automatic build. To read more about modifying automatic update and build options, please read [LoopDocs: Automatic Update & Build](#).

## Disable Automatic Install from *TestFlight*

Once the app is available in *TestFlight*, you can adjust whether it is automatically installed when updated versions are available. We recommend you disable automatic installation so you can choose when to switch to a newer build, which in some cases, may be a newer version.

Go back to the *TestFlight* app on your phone and tap on your app name in the list to see an expanded screen similar to the graphic below. The row to enable or disable automatic updates is highlighted in the graphic, which shows the feature disabled. This is recommended for all users.

- If you leave automatic update enabled (default), then whenever a new build is created and uploaded to *TestFlight*, it will be installed immediately (see [Unexpected \*TestFlight\* Beta Expiration](#))
- WARNING: If you switch between Building with Browser and Mac, you must disable automatic update or Xcode will not be able to install to your phone

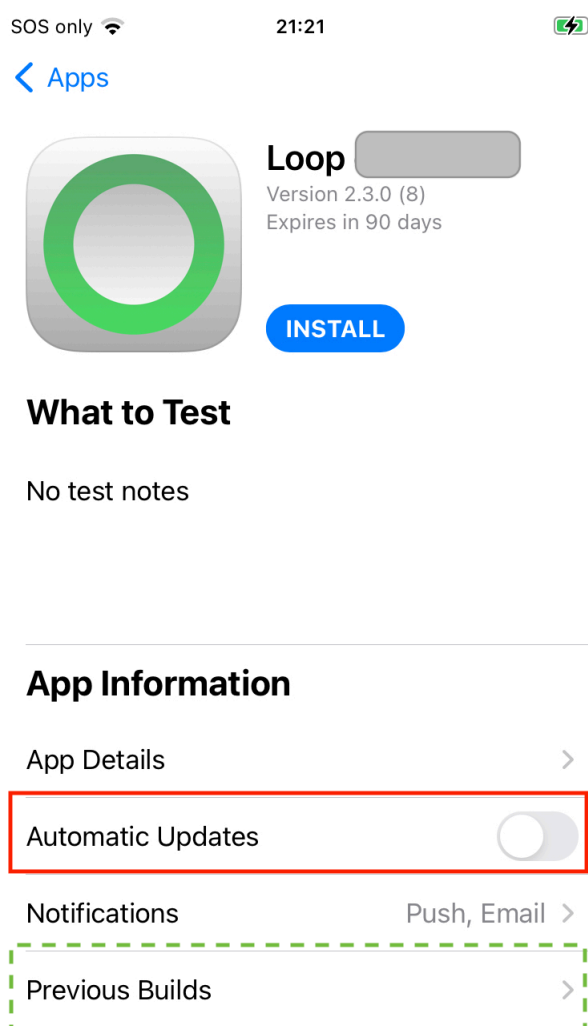

When you are ready to install, just open the *TestFlight* app and click Install to get the most recent build and then click Open when it completes the installation. All your settings and connections to CGM and Pump are maintained.

## Previous Builds

If you tap on the row that says `Previous Builds`, highlighted by the dashed-green rectangle in the graphic above, you can view and choose an older (or lower version number) build (as long as it has not expired).

- In some cases, you need to do this to see the newest build if the most recent build has a lower version number than an older build

## Unexpected `TestFlight` Beta Expiration

### ! **Apple `TestFlight` Bug**

It doesn't happen to everyone but it has happened a few times over the years and only for people who have Automatic Update enabled for their `TestFlight` app.

Symptom:

- App has plenty of time before expiration
- A new build becomes available, for example from the automatic monthly rebuild
- User is told *Trio* Beta is not available

Solution:

- Open `TestFlight` and turn off automatic installation for your app
- Install the previous version of the app on the phone manually (with monthly build, it should still be ok)
- Manually start a new build (if there were updates you want to get)
- Manually install the new update

More information:

- [Link 1: Apple forum](#)
- [Link 2: Stackoverflow thread](#)

## `TestFlight` for a Child

Minor children are not allowed to install or use the `TestFlight` app.

An adult, who is a member of the [Internal `TestFlight` Group](#) ↗ can log into Media & Purchase (see steps below) without affecting the Apple ID associated with a phone (and thus the health records used by the *Trio* app for the minor child). After the adult installs or updates the app using `TestFlight`, they probably should reverse those steps to remove their credentials from Media & Purchase.

Media & Purchase affects access to the App Store, Books, Music and Podcasts.

On the Child phone:

- Tap on Settings

- At the very top of Settings, tap on the Name of the phone, for example, `my kids phone`
- Apple ID Screen appears
  - Tap on Media & Purchases
  - Tap on Sign Out, and confirm
  - Sometimes the phone requires a reboot before you can sign in with a different ID
- Sign in with the adult Apple ID email and password
- Install or Update the app from TestFlight on child phone
- Repeat the process to sign out the adult and (if needed) sign back in the child

## Change the App Store Connect Name

Suppose you really don't like the name you picked initially for the *Trio* app that shows up in the TestFlight app.

You can change it.

Open this link: [App Store Connect Apps](#) and log in as needed.

- Click on your app name.
- Click on App Information on the left side (make browser wider if you don't see this).
- Modify the Name under Localizable Information and click on the Save button (upper left)
- If you chose a name that is in use, you'll see the warning screen - try again