



TrioDocs

Version: 0.4.0

Date: May 28, 2026

Download the latest version on:
<https://triodocs.org>

Prepare App

New Trio Builders

Create the Trio App Group

If you already have a Trio App Group

- You can skip this step - your existing App Groups are found at this link: [App Group List](#)
- If your Trio App Group was created from a Mac with Xcode, you may choose to edit the Description to make the **NAME** match

If you do not have a Trio App Group :

- Go to [Register an App Group](#) on the Apple developer site and use the table below to help you create one.
- Replace TEAMID with your Apple Developer ID.

NAME	Xcode version (NAME)	IDENTIFIER
Trio App Group	group org nightscout TEAMID trio trio-app-group	group.org.nightscout.TEAMID.trio.trio-app-group


Table of Identifiers

These are the Identifiers created by running the [GitHub](#) action " Add Identifiers ". Only 3 of them need to have the Trio App Group added.

NAME	Xcode version (NAME)	IDENTIFIER
Trio	XC org nightscout TEAMID trio	org.nightscout.TEAMID.trio
Trio LiveActivity	-	org.nightscout.TEAMID.trio.LiveActivity
Trio Watch App	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp
Trio Watch Complication	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp.TrioWatchComplication

Add Trio App Group to Identifiers

- Open the [App IDs Identifier page for your Apple Developer Account](#).
 - Click on the Trio Identifier and assign the Trio App Group to the Identifier - see graphic below.

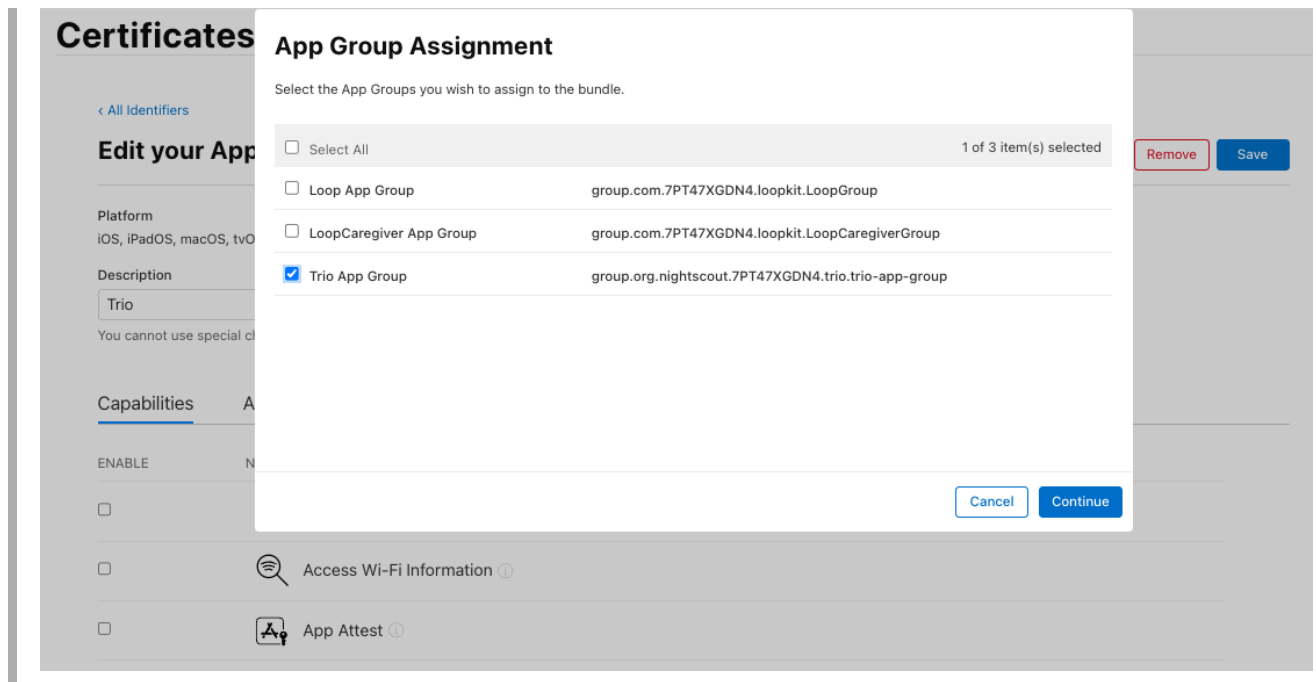

- Repeat this for the other 2 Identifiers that need to have an App Group assigned: Trio Watch App and Trio Watch Complication

Configure Trio App

Follow the directions in [LoopDocs](#), but use the Trio Bundle ID

- In App Store Connect, the Bundle ID for Trio will be: org.nightscout.TEAMID.trio

If you want to take it a step at a time, continue to [Create Certificates](#).

If you are feeling confident, you can go straight to [Build App](#) which includes the create certificates step.

Navigation: [← Back: Add Identifiers](#) | [Next: Prepare TestFlight Users →](#)