



TrioDocs

Version: 0.4.0

Date: September 20, 2025

Download the latest version on:
<https://triodocs.org>

Build Beta with Browser

Use Browser Build for Public Beta Version

This page is provided to help people currently running the public beta version of Trio.

! Public Beta

The public beta is built from the `dev` branch of Trio.

This code is undergoing testing and receives rapid updates. Be sure to follow along in Trio Discord and to update frequently (at least weekly).

- [Trio Discord Invitation](#)

To ensure weekly updates, make the `dev` branch the default branch in your fork on GitHub. Remember, this will just automatically trigger a build of the latest `dev` to TestFlight every Wednesday at 08:43 UTC. You are still responsible for manually installing the update via the TestFlight app on your phone.

Once this beta is released to `main` as 1.0, then you should return to building from the `main` branch.

Update to Trio `dev` branch from Trio 0.2.x main

! A Note on Compatibility

Upgrading to Trio 0.5.x from 0.2.x is smooth and straightforward. The new Onboarding Wizard will guide you step by step. Your pump, CGM, therapy settings, and the last 24 hours of treatment and glucose history will be brought over automatically.

Be Aware

- **Once you have upgraded to 0.5.x, going back to 0.2.x is not supported**
 - If you choose to downgrade, you will need to set everything up again from scratch
 - If you accidentally build 0.2.x over 0.5.x, and then return to 0.5.x:
 - You will see "Oops? Some data didn't make it over" every time you restart the app
 - But all your "stuff" is there
 - See [Remove Annoying Message](#)
- Some guardrails have changed compared to 0.2.x
- **Your saved Temp Targets and Overrides will not be maintained**
 - The storage method in 0.5.x is different from 0.2.x, so you will need to recreate them
 - Before you update from 0.2.x, capture a screenshot of each named Temp Target or Override that you want to add to 0.5.x

The Browser Build documentation is under construction for Trio version 0.5.x and newer.

This temporary set of instructions is intended for the subset of users who **previously built Trio 0.2.x using Browser Build**. Because this is aimed at experienced builders, some of the steps are abbreviated and will be expanded as the documents are updated.

If you are new to Trio — please wait until the documents have been expanded and preferably until the release of Trio 1.0.

Mac-Xcode to Browser Build

! Experienced Mac-Xcode Builder

If you are an experienced Trio 0.2.x user who wants to join the open beta testing, but also wants to switch from Mac-Xcode build to Browser Build — welcome. For the time being, please do this:

- Use the [0.2.x: Build Trio with GitHub](#) instructions to build Trio 0.2.x and wait for it to show up in your TestFlight. Do not install it — you just want to make sure you can succeed with a Browser Build.
- Then return to this page and follow the directions; you will need to run `Add Identifiers` and `Create Certificates` again when updating to 0.5.x.

Transition from Other Apps to Trio Browser Build

! Coming from Other Apps

If you are an experienced Browser Build person who has not built Trio before, you will need to add the Trio App Group to 3 `Identifiers` and create a Trio App in App Store Connect.

The steps are summarized in [New Trio Builders](#).

If you are not experienced with Browser Build, we suggest you follow the 0.2.x instructions: [0.2.x: Build Trio with GitHub](#) first. Then return to this page to update to 0.5.x; you will need to run `Add Identifiers` and `Create Certificates` again when updating to 0.5.x.

Summary of Tasks to Build Trio 0.5.x

These steps assume:

- You previously built Trio 0.2.x using [GitHub Actions](#) (Browser Build)
- You confirmed your *Apple Developer* license agreements are up to date

These are the new steps for you to follow:

1. [Configure Browser Build Certificate Automation](#)
2. [Configure Fork with dev branch](#)
 - [Configure the dev branch as default](#) (while public beta is ongoing)
3. [Update Identifiers](#)
4. [Update Certificates](#)

5. Build Trio 0.5.x

Browser Build Certificate Automation

Browser Build Certificate Automation was added to Trio 0.2.3 - if you have not added the `ENABLE_NUKE_CERTS` variable, you should add it now. Otherwise, skip ahead to [Configure Fork with dev branch](#).


In order to utilize the new automatic certificate renewal feature, you'll need to add a new Variable. Variables are located in *GitHub*, in the same location as your Secrets. The exact location will depend upon whether you build using a *GitHub organization* or a personal account.

If you use a personal account, click on your Trio repository. If you have other repositories, just follow these same steps for each of them.

If you build using a *GitHub organization* and have already added this variable to your organization, there's nothing for you to do. All repositories in your organization are covered. Otherwise, click on your organization name.

The numbered steps correspond to numbers in the graphics below:

1. Choose settings
2. Scroll down to select Secret and variables
3. Choose Actions
4. Choose Variable
5. Tap on "Create new organization variable" or "New repository variable"
6. For the Name, type `ENABLE_NUKE_CERTS`
7. For the Value, type `true`
8. Tap on Add Variable.

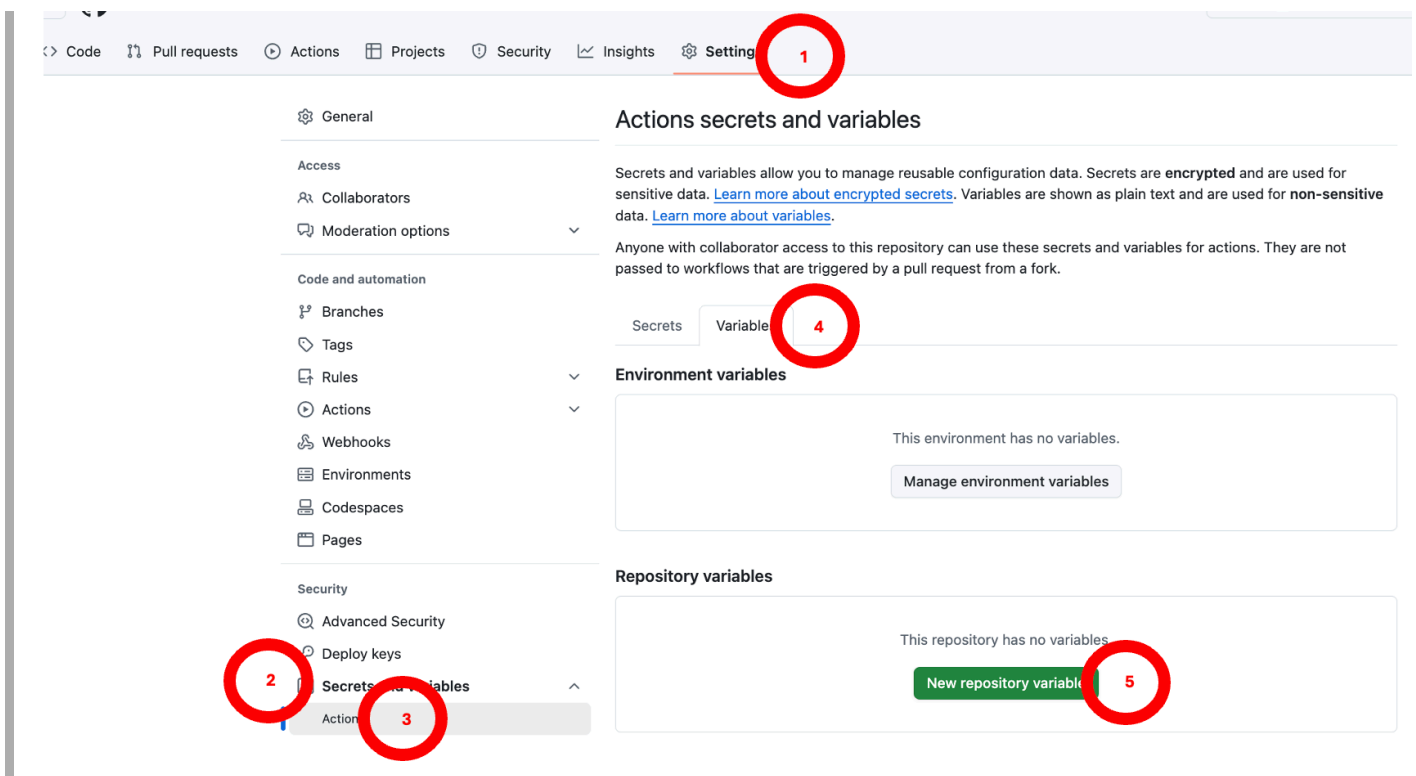

Actions variables / New variable

Note: Variable values are exposed as plain text. If you need to encrypt and mask sensitive information, [create a secret](#) instead.

Name *

ENABLE_NUKE_CERTS

6

- Variable names may only contain alphanumeric characters ([a-z], [A-Z], [0-9]) or underscores (_).
- Variable names cannot start with a number.
- Variable names cannot start with GITHUB_ prefix.

Value *

true

7

Add variable

8


Configure Fork with dev branch

The open beta testing for Trio uses the dev branch.


- If you do not have a dev branch you must first configure one following these directions.
- If you already have a dev branch, skip ahead to [Update Branch](#) and be sure to select the dev branch when you update your fork .

Add Branch

First, you need to get to your Trio repository and tap on the branch icon. (Do not worry about how many branches you have.):

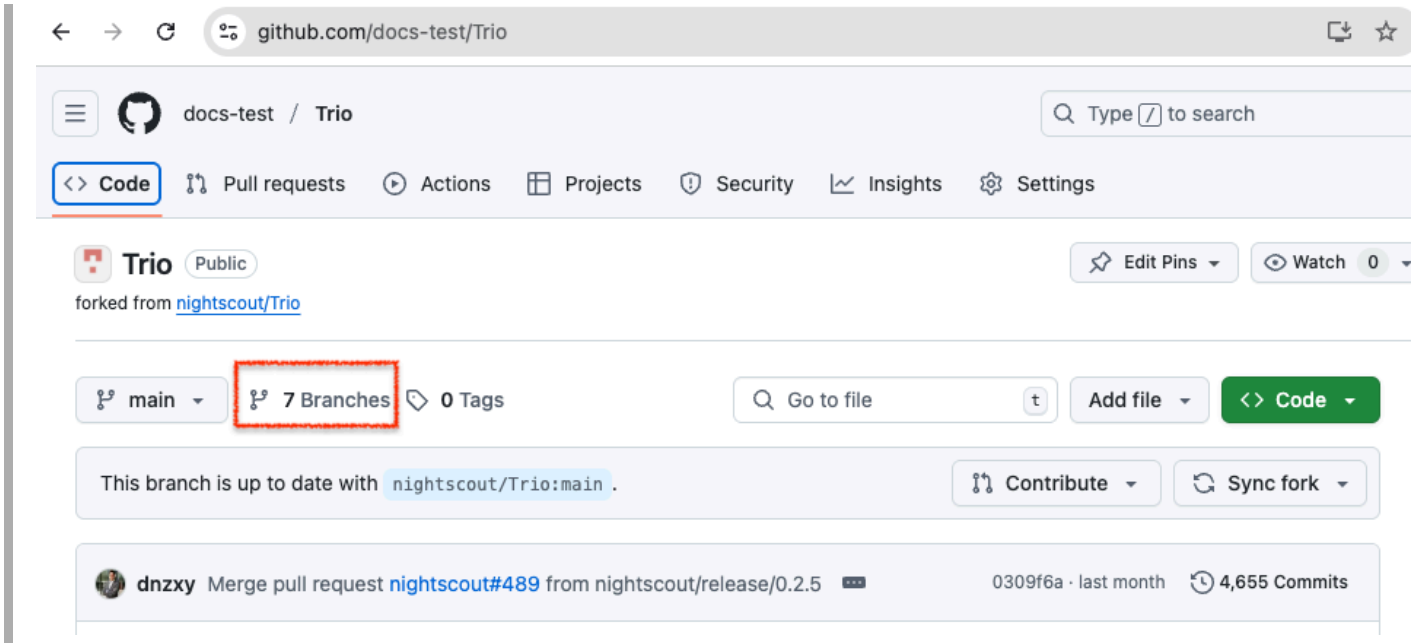


The next screen displays the branches you currently have. It shows a **New branch** button in the upper right. If you don't see that, you are not logged in. Tap on the **New branch** button.

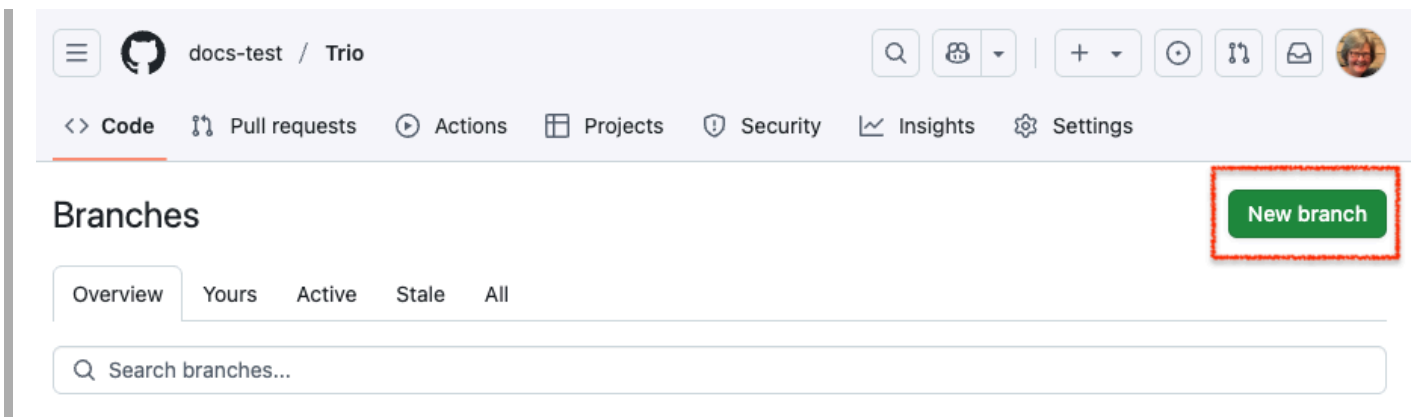

Each step in the list below matches with the number in the graphic. On the left half of the graphic is the default selections for your **fork** when **main** is the default branch. The right side shows the display after making the indicated selections.:

1. Click on the drop-down menu labeled 1 in the graphic and choose nightscout/Trio
2. Click on the drop-down menu labeled 2 in the graphic and choose **dev**
3. Click on the **Branch** name box labeled 3 in the graphic and type **dev**
 - The **branch** name in your **fork** should always match the **branch** name you are adding; check that you type it correctly
4. Review the dialog items to make sure everything is correct, and then tap on **Create branch**

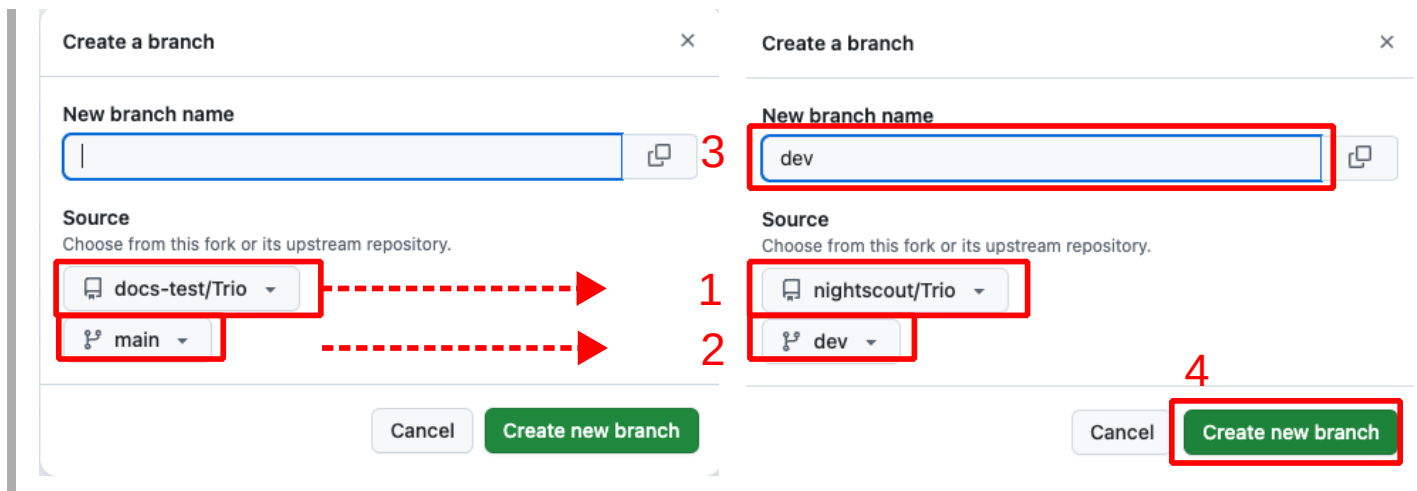

Configure the `dev` branch as default

This recommendation is only while the public beta is ongoing. Typically only developers configure the `dev` branch as their default branch.


1. Once you commit to the 0.5.x version of Trio, you want to stay with it
2. By making the `dev` branch your default, you will get automatic build updates weekly

These are the steps to modify the default branch.

For this example, we show how to change from a default branch of `main` to a default branch of `dev`. Note — only the owner of the repository can take this action and they must be logged in. Otherwise the Settings tab does not appear.

For the numbered steps below, refer to the graphic found under each group of steps.


1. Click on the Settings Icon near the top right of your Trio repository
 - You may need to scroll down to see the `Default Branch` as shown in the graphic
 - Do not tap on the Branches tab to the left under Code and Automation, that is not the correct menu

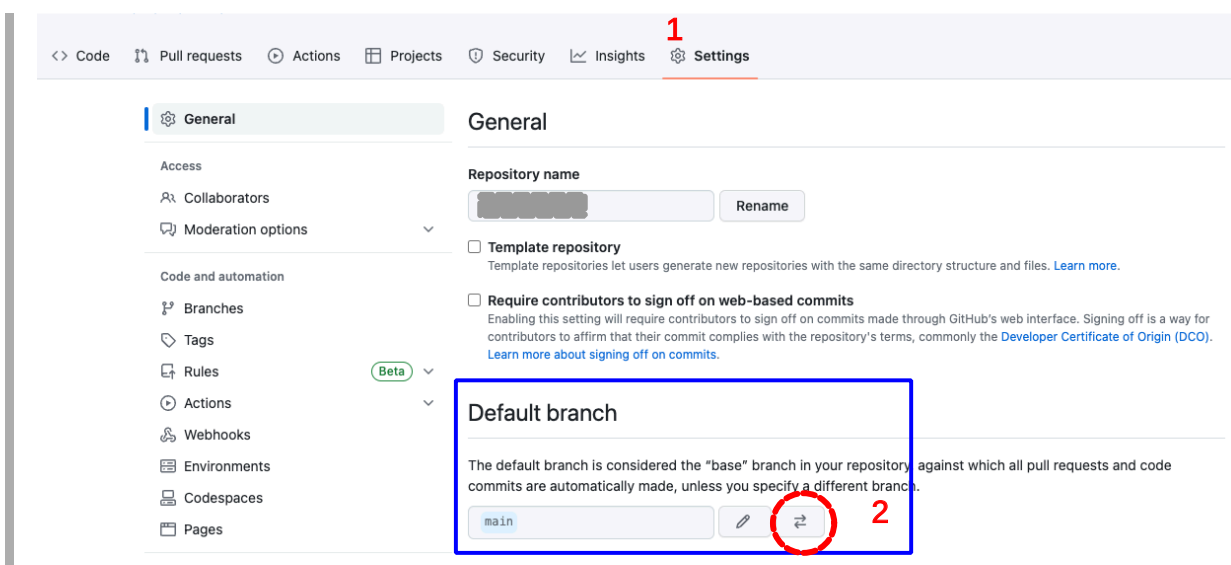

2. To the right of the default branch name there is a pencil and a left-right arrow icon (↔)

- Tap on the left-right arrow icon (↔) to bring up the `Switch default branch to another branch` dialog

3. Click on the dropdown next to the current default branch, in this example, `main`

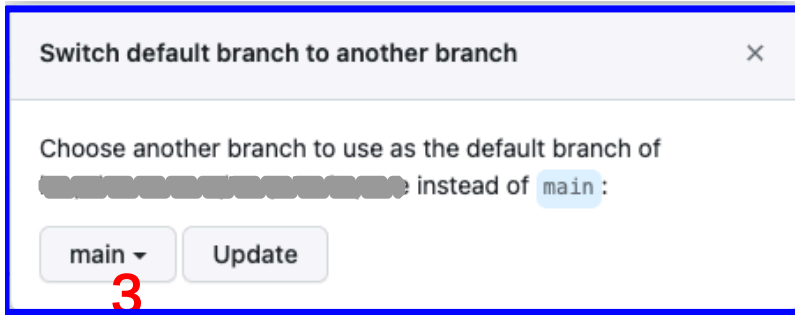

1. Select the desired default branch, for the public beta, choose `dev`

2. Click on the `Update` button

3. You will be presented with an are-you-sure question.

- Click on the red `I understand, update the default branch.` button


Your default branch has been changed.

Update Branch

The `dev` branch of the `nightscout/Trio` repository will be updated frequently. This is how you know if your fork needs to be updated as well. As long as the automatic runs are happening weekly for your repository – you will not need these manual instructions, unless you want to manually build an update sooner.

Tap the `Code` button (upper left) and ensure this branch in your fork is up to date.

- Select the desired branch in the dropdown menu (this graphic shows `main` branch, to get 0.5.x, you must choose `dev` branch)
- If the message indicates this branch is "behind", tap on the sync fork button and then the `Update branch` button

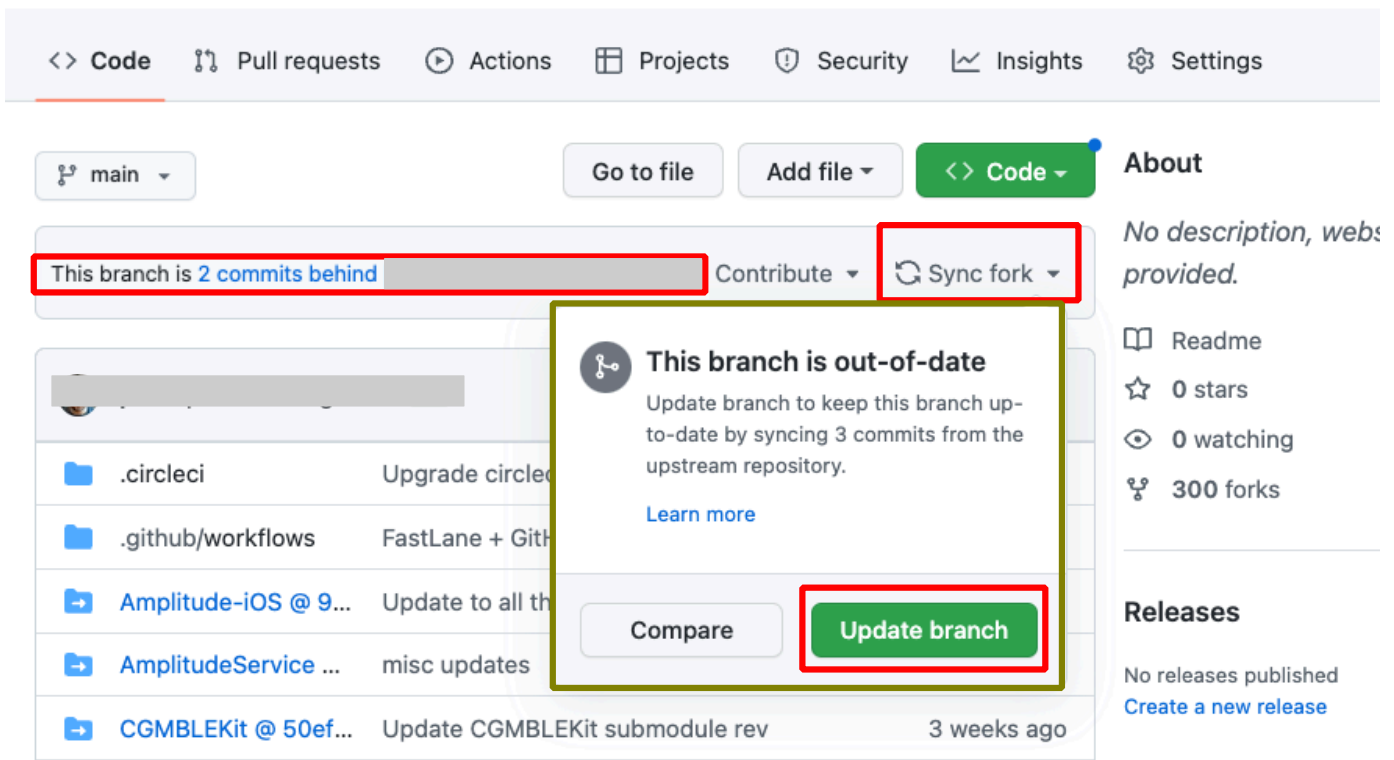


Update Identifiers

For Trio 0.5.x and newer, you must have Apple Push Notification enabled to build the app. This capability is added to the existing Trio Identifier by running the Action: `Add Identifiers` after you update your `fork`.

Refer to the graphic below for the numbered steps:

1. Click on the `Actions` tab of your Trio repository
2. On the left side, click on `2. Add Identifiers`
3. On the right side, click `Run Workflow` to show a dropdown menu
 - You will see your default `branch` (typically this is `main`)
 - **To update the `Identifiers` so you can build the Trio 0.5.x, you must select `dev`**
4. Tap the green button that says `Run workflow`.

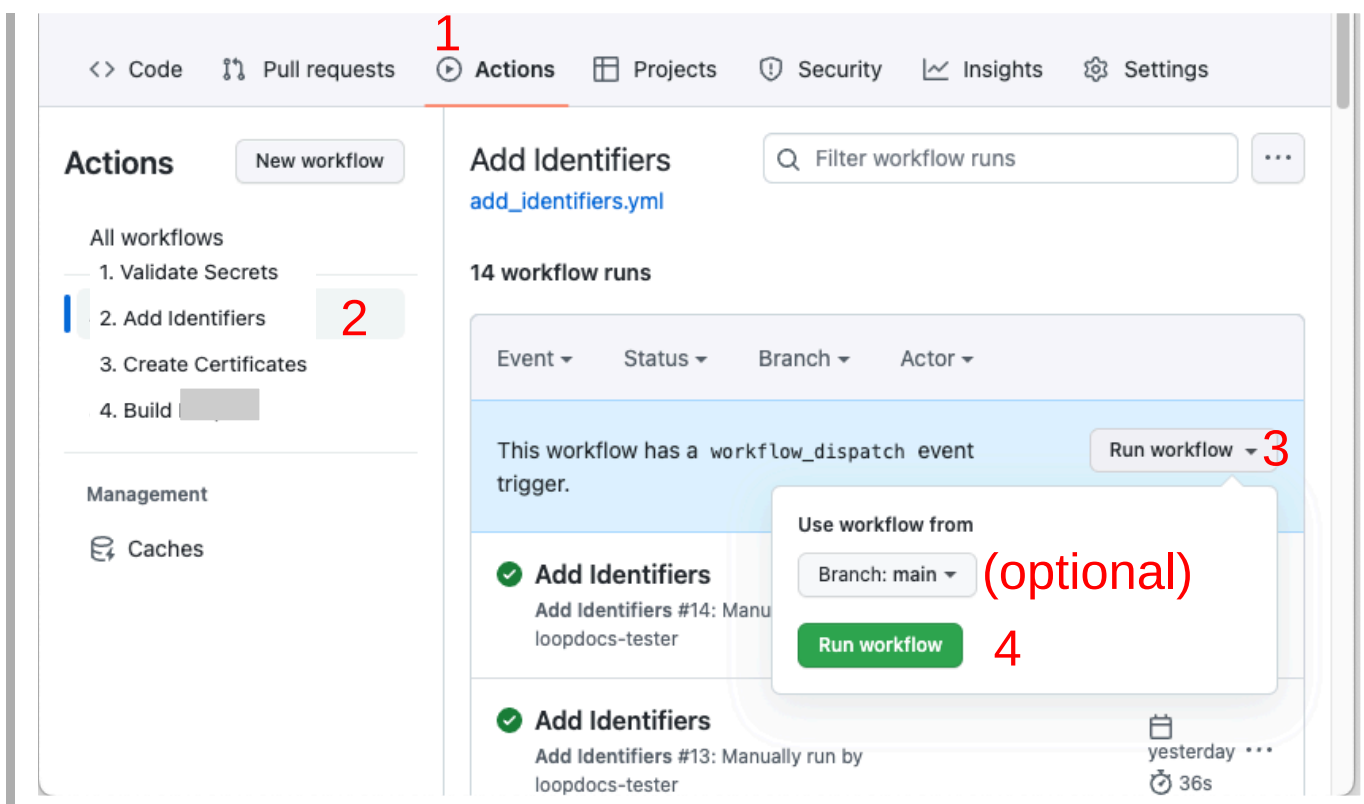

Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes


The `Add Identifiers` Action should succeed or fail in a few minutes. Do not continue to the next step until this one succeeds.

Update Certificates

For Trio 0.5.x and newer, you must have Apple Push Notification enabled to build the app. After you add this capability to the Trio Identifier, you must update the Certificate by running the Action: `Create Certificates`.

Refer to the graphic below for the numbered steps:

1. Click on the `Actions` tab of your Trio repository
2. On the left side, click on `3. Create Certificates`
3. On the right side, click `Run workflow` to show a dropdown menu
 - You will see your default branch (typically `main`)
 - **Because you plan to build Trio 0.5.x, select `dev`**
4. Tap the green button that says `Run workflow`.

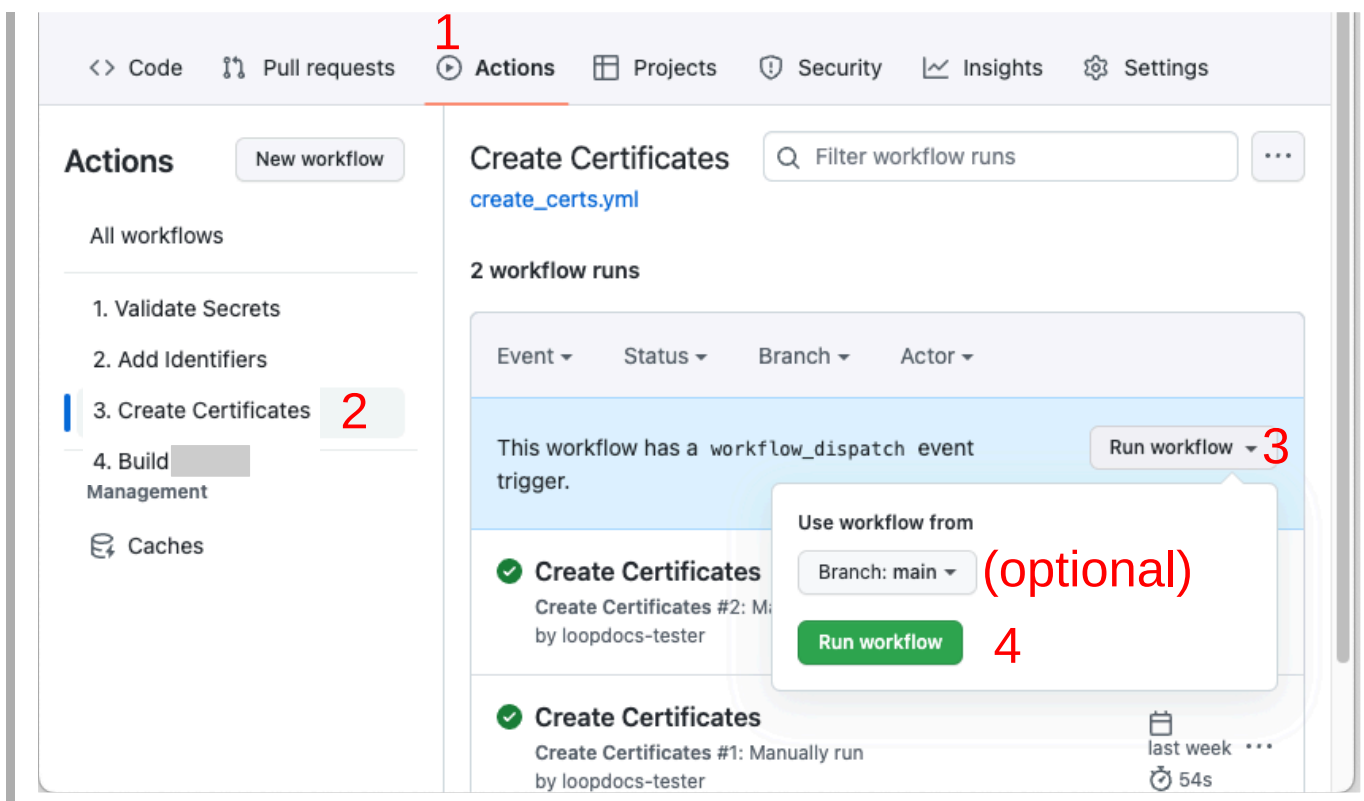

5. Wait a minute or two for the action to finish

Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Once you see the green check mark by `Create Certificates`, the next step is to Build the app.

Build the App

If you completed all the steps on this page successfully (got a green checkmark ) , you are ready to run Action: Build Trio.

We recommend that while participating in the public beta, you should **Configure the dev branch as default**. This will ensure weekly updates to **TestFlight**. You must still manually install updates to your phone from **TestFlight**.

If you choose to build a different branch than your default branch, there is an extra step when you **Build Trio**. In addition to the normal steps 1, 2, and 3 in the graphic below, you must also do the (optional) step. Select the **dev** branch in the **branch** dropdown menu before continuing to step 4 and tapping on the green **Run workflow** button:

The screenshot shows the GitHub Actions interface. The top navigation bar has tabs for Code, Pull requests, Actions (highlighted with a red '1'), Projects, Security, Insights, and Settings. On the left, under 'Actions', there is a list of workflows: '1. Validate Secrets', '2. Add Identifiers', '3. Create Certificates', '4. Build Trio' (highlighted with a red '2'), '8. DONT RUN Add ... Disabled', and '8. DONT RUN close... Disabled'. Below this is a 'Management' section with 'Caches'. The main area shows the '4. Build Trio' workflow details, including a search bar 'Filter workflow runs' and a '133 workflow runs' section. A modal is open showing 'This workflow has a workflow_dispatch event trigger.' and a 'Run workflow' button (highlighted with a red '3'). Below this, there is a 'Use workflow from' section with a 'Branch: dev' dropdown (highlighted with a red '(optional)') and a 'Run workflow' button (highlighted with a red '4').

Refresh, Do Not Repeat

Hit Refresh


After you tap the green **Run workflow** button, **GitHub** can be slow to update.

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Disable the DO NOT RUN Actions

You may notice some **GitHub** actions that have **DONT RUN** or **DO NOT RUN** in their names. Those are special actions that are used by the developers. They are skipped in your **Fork**, but you may see logs for them.

If they already have the notation **disabled** beside them, so won't see them run. Otherwise, click on each one, click on the three dots at the upper right, and then select **Disable workflow**. Then you will avoid the annoyance of seeing logs that say the action was automatically run and then skipped.

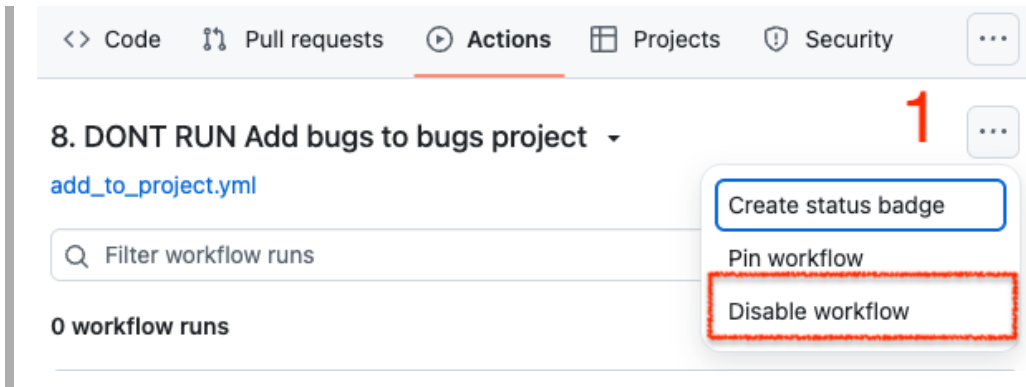

Update Build Errors

The most likely build error is that you did not [Update Identifiers and Certificates](#). In that case, you will see this **error** in the **build action annotations**.

Error



If you decide to look at the action **log**, instead of reading the annotations, you may see an **error** similar to the two graphics below.



```

518 [23:27:00]: ▶ ::error file=/Users/runner/work/Trio-dev/Trio-dev/Trio.xcodeproj::Provisioning profile "match AppStore org.nightscout.***.trio" doesn't
support the Push Notifications capability. (in target 'Trio' from project 'Trio')
519 [23:27:00]: ▶ /Users/runner/work/Trio-dev/Trio-dev/Trio.xcodeproj: error: Provisioning
profile "match AppStore org.nightscout.***.trio" doesn't include the aps-environment
entitlement. Profile qualification is using entitlement definitions that may be out of
date. Connect to network to update. (in target 'Trio' from project 'Trio')
520 [23:27:00]: ▶ note: Disabling previews because SWIFT_VERSION is set and
SWIFT_OPTIMIZATION_LEVEL=-O, expected -Onone (in target 'Trio' from project 'Trio')
521 [23:27:00]: ▶ note: Run script build phase 'Swiftformat' will be run during every build
because the option to run the script phase "Based on dependency analysis" is unchecked.

```

Solution

1. Update [Identifiers](#) - make sure you select `dev` branch
2. Update [Certificates](#) - make sure you select `dev` branch

New Trio Builders

For experienced browser builders, this minimal set of instructions might be sufficient.

Warning - there is no hand holding in these directions.

Every app you build will use the same 6 Secrets.

- Fork from: <https://github.com/nightscout/Trio>
- Use the `Trio App Group` for Trio, see [Create the Trio App Group](#)
- [Identifiers](#) for Trio: see [Table of Identifiers](#)
- Add the `Trio App Group` to these [Identifiers](#):
 - `Trio`
 - `Trio Watch`
 - `Trio WatchKit Extension`
- In `App Store Connect`, the `Bundle ID` for Trio will be: `org.nightscout.TEAMID.trio`

Create the `Trio App Group`

If you already have a `Trio App Group`

- You can skip this step - your existing App Groups are found at this link: [App Group List](#)
- If your `Trio App Group` was created from a Mac with `Xcode`, you may choose to edit the Description to make the **NAME** match

If you do not have a `Trio App Group`:

- Go to [Register an App Group](#) on the Apple developer site and use the table below to help you create one.
- Replace `TEAMID` with your *Apple Developer* ID.

NAME	Xcode version (NAME)	IDENTIFIER
Trio App Group	group org nightscout TEAMID trio trio-app-group	group.org.nightscout.TEAMID.trio.trio-app-group


Table of Identifiers

These are the Identifiers created by running the GitHub action "Add Identifiers". Only 3 of them need to have the Trio App Group added. Be sure that you selected the dev branch when you start the "Add Identifiers" action. For more details, see Update Identifiers , but then return here when that action succeeds.

NAME	Xcode version (NAME)	IDENTIFIER
Trio	XC org nightscout TEAMID trio	org.nightscout.TEAMID.trio
Trio LiveActivity	-	org.nightscout.TEAMID.trio.LiveActivity
Trio Watch App	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp
Trio Watch Complication	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp.TrioWatchComplication

Add Trio App Group to Identifiers

- Open the App IDs Identifier page for your Apple Developer Account.
- Click on the Trio Identifier and assign the Trio App Group to the Identifier - see graphic below.

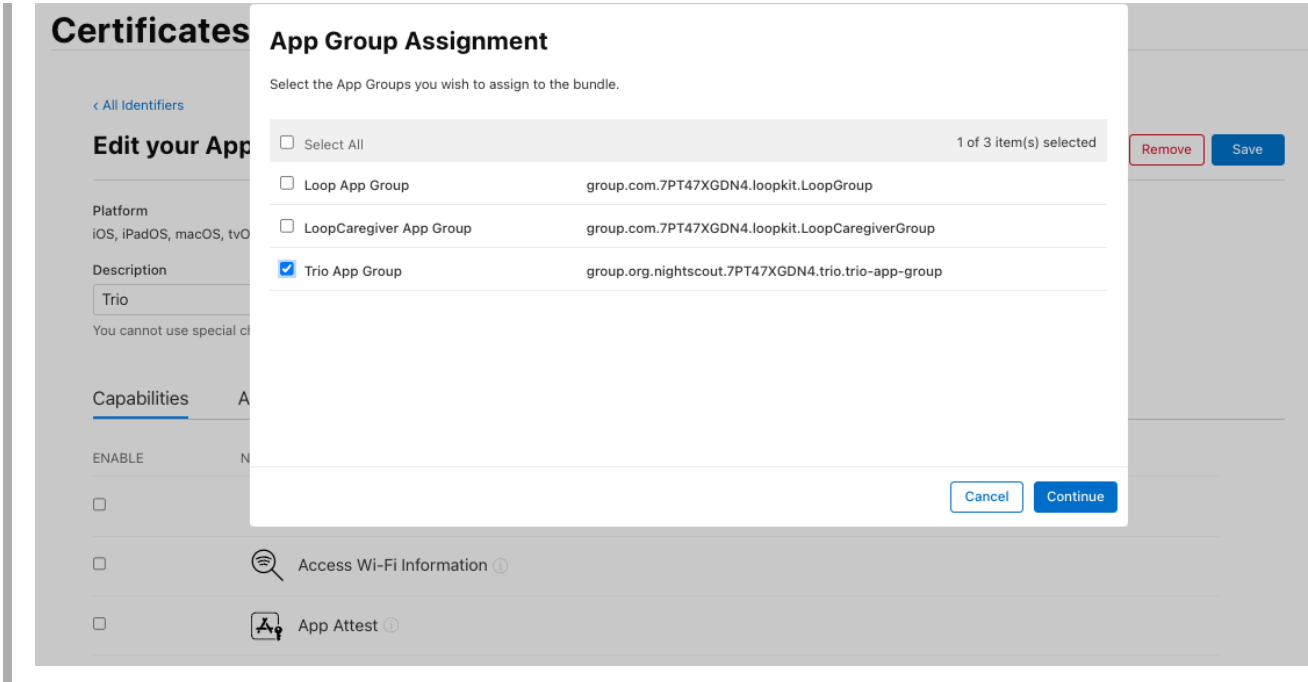

- Repeat this for the other 2 `Identifiers` that need to have an `App Group` assigned: `Trio Watch` and `Trio WatchKit Extension`

Configure Trio App

Follow the directions in [LoopDocs](#), but use the Trio `Bundle ID`

- In `App Store Connect`, the `Bundle ID` for Trio will be: `org.nightscout.TEAMID.trio`

Return to the main set of instructions on this page, [Update Certificates](#), and keep going until you have a successful build.

Extra Tips

Remove Annoying Message

If you accidentally build 0.2.x over 0.5.x, and then return to 0.5.x:

- You will see `"Oops? Some data didn't make it over"` every time you restart the app
- But all your "stuff" is there

You see that because some special files are labeled with `.migrated` by 0.5.x when you complete the onboarding action when you upgrade from 0.2.x to 0.5.x.

- When you accidentally build 0.2.x over 0.5.x, the old names for the files are recreated
- When you restore to 0.5.x, it sees the duplicate files, one set has `.migrated` in the file name and one set does not
- You are warned that migration can't happen - but for this case, you don't need it

You can get rid of the annoying message by taking these steps:

First quit (stop running, swipe up) the Trio app

- Search for Files on your phone
- Open the `Trio` folder on your phone
- Open the `enact` folder
 - you may see 4 files, instead of 2
 - delete the ones that do **not** have `.migrated` in the name
- Open the `monitor` folder
 - you will probably see many files, but 3 of them have `.migrated` in the name
 - find the files with the same name except for `.migrated` and delete them

Open the Trio app and it should just say "Getting everything ready for you" and just keep going. Annoying message is gone.