

Update with Browser 



TrioDocs

Version: 0.4.0

Date: May 19, 2025

Download the latest version on:  
<https://triodocs.org>

Update with  
Browser 

# Update to Trio dev branch from Trio 0.2.x main

## ! A Note on Compatibility

Upgrading to Trio `v0.5.0` from `0.2.x` is smooth and straightforward. The new Onboarding Wizard will guide you step by step. Your pump, CGM, therapy settings and the last 24 hours of treatment and glucose history will be brought over automatically.

### Be Aware

- Once you have upgraded to `v0.5.0`, going back to `0.2.x` is not supported.
- If you choose to downgrade, you will need to set everything up again from scratch.

We have designed this release to be a stable and long-term foundation—so when you are ready, you can upgrade with confidence.

The Browser Build documentation is under construction for versions Trio 0.5.x and newer.

This temporary set of instructions is intended for the subset of users who **previously built Trio 0.2.x using Browser Build**. Because this is aimed at experienced builders, some of the steps are abbreviated and will be expanded as the documents are updated.

If you are new to Trio - please wait until the documents have been expanded and preferably until the release of Trio 1.0.

If you are an experienced Trio 0.2.x user who wants to join the open beta testing, but also wants to switch from Mac-Xcode build to Browser Build - welcome. For the time being, please do this:

- Use the [0.2.x: Build Trio with GitHub](#) instructions to build Trio 0.2.5 and wait for it to show up in your TestFlight. Do not install it - you just want to make sure you can succeed with a Browser Build.
- Then return to this page and follow the directions

## Summary of Tasks to Build Trio 0.5.x

These steps assume:

- You previously built Trio 0.2.x using [GitHub Actions](#) (Browser Build)
- You confirmed your *Apple* Developer license agreements are up to date

These are the new steps for you to follow:

1. [Configure Browser Build Certificate Automation](#)
2. [Configure Fork with dev branch](#)
3. [Update Identifiers](#)
4. [Update Certificates](#)
5. [Build Trio 0.5.x](#)

# Browser Build Certificate Automation

Browser Build Certificate Automation was added to Trio 0.2.3 - if you have not added the `ENABLE_NUKE_CERTS` variable, you should add it now. Otherwise, skip ahead to [Configure Fork with dev branch](#).


In order to utilize the new automatic certificate renewal feature, you'll need to add a new Variable. Variables are located in [GitHub](#), in the same location as your Secrets. The exact location will depend upon whether you build using a [GitHub organization](#) or a personal account.

If you use a personal account, click on your Trio repository. If you have other repositories, just follow these same steps for each of them.

If you build using a [GitHub organization](#) and have already added this variable to your organization - there's nothing for you to do. All repositories in your organization are covered. Otherwise, click on your organization name.

The numbered steps correspond to numbers in graphics below:

1. Choose settings
2. Scroll down to select Secret and variables
3. Choose Actions
4. Choose Variable
5. Tap on "Create new organization variable" or "Create new repository variable"
6. Add the variable `ENABLE_NUKE_CERTS`
7. Type in the value `true`
8. Tap on Add Variable.

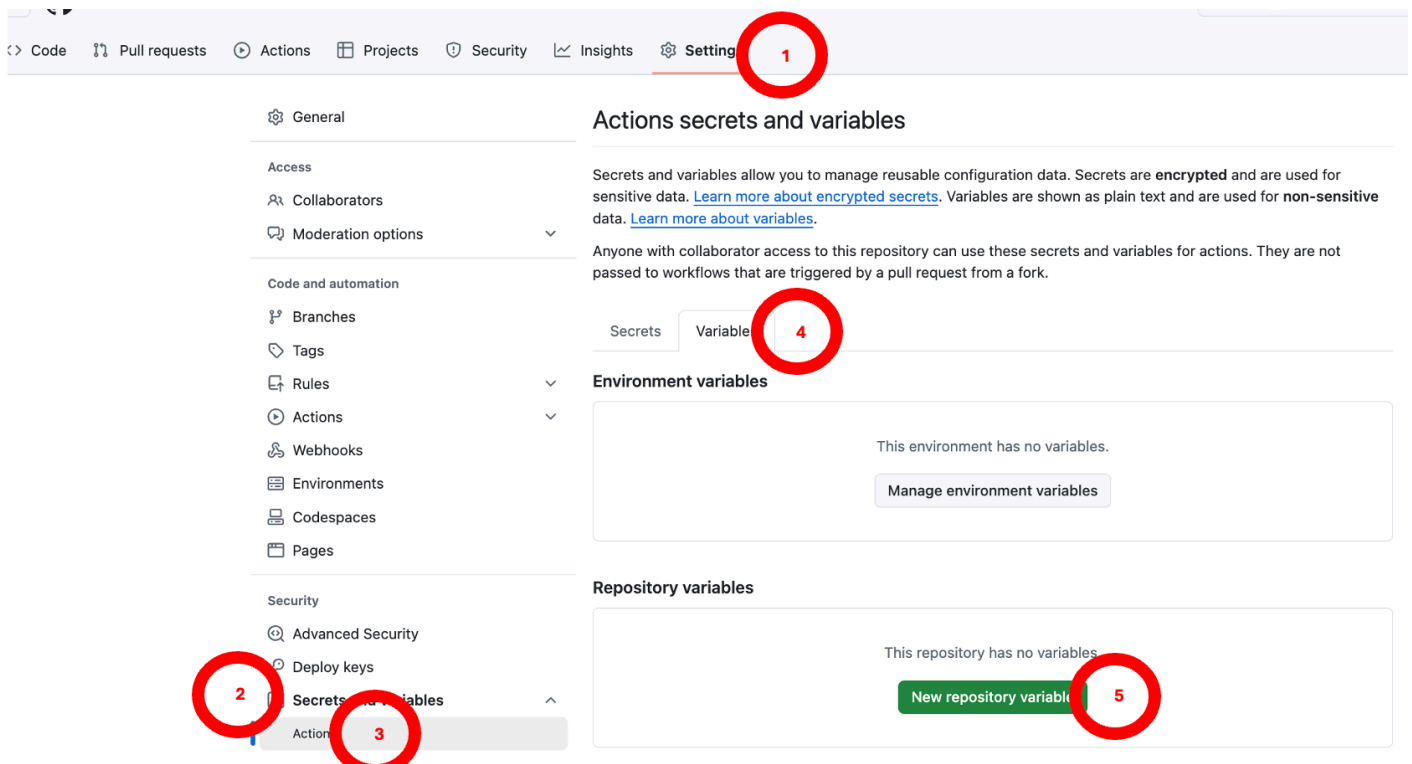

## Actions variables / New variable

Note: Variable values are exposed as plain text. If you need to encrypt and mask sensitive information, [create a secret](#) instead.

Name \*

ENABLE\_NUKE\_CERTS

6

- Variable names may only contain alphanumeric characters ([a-z], [A-Z], [0-9]) or underscores (\_).
- Variable names cannot start with a number.
- Variable names cannot start with GITHUB\_ prefix.

Value \*

true

7

Add variable

8

## Configure Fork with dev branch

The open beta testing for Trio uses the dev branch.

- If you do not have a dev branch you must first configure one following these directions.
- If you already have a dev branch, skip ahead to [Update Branch](#) and be sure to select the dev branch when you update your fork.

### Add Branch

Each step in the list below matches with the number in the graphic. In the top half of the graphic, the left side shows the initial display and the right side shows the display after making the indicated selections.:

1. Click on the drop down menu labeled 1 in the graphic and choose nightscout/Trio similar to the top right graphic
2. Click on the drop down menu labeled 2 in the graphic and choose dev
3. Click on the Branch name box labeled 3 in the graphic and type dev
  - The branch name in your fork should always match the branch name you are adding; check that you type it correctly
4. Review the dialog items to make sure everything is correct and then tap on Create branch

Create a branch

New branch name

Source  
Choose from this fork or its upstream repository.

docs-test/Trio 1

main 2

Cancel Create new branch

Create a branch

New branch name

Source  
Choose from this fork or its upstream repository.

nightscout/Trio 1

main 2

Cancel Create new branch

Create a branch

New branch name

dev 3

Source  
Choose from this fork or its upstream repository.

nightscout/Trio

main


Cancel Create new branch 4

## Update Branch

Be aware that building from the `dev` branch will upgrade to Trio 0.5.x. Once you have upgraded to 0.5.0, going back to 0.2.x is not supported.

Tap the `Code` button (upper left) and ensure this branch in your fork is up to date.

- Select the desired branch in the dropdown menu (this graphic shows `main` branch, to get 0.5.0, you must choose `dev` branch)
- If the message indicates this branch is "behind", tap on the sync fork button and then the Update branch button

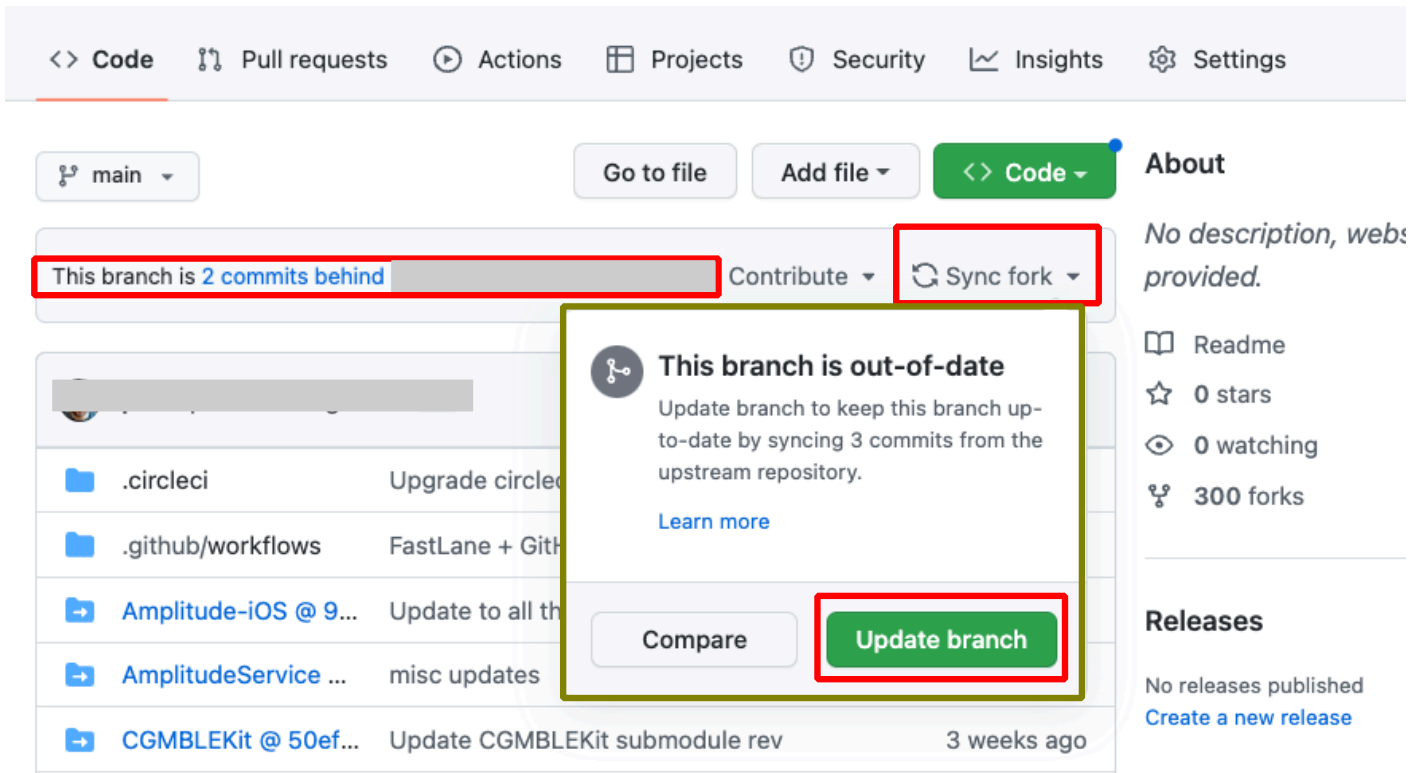


## Update Identifiers

For Trio 0.5.x and newer, you must have *Apple* Push Notification enabled to build the app. This capability is added to the existing Trio Identifier by running the Action: Add Identifiers after you update your fork.

Refer to the graphic below for the numbered steps:

1. Click on the `Actions` tab of your `Trio` repository
2. On the left side, click on 2. `Add Identifiers`
3. On the right side, click `Run Workflow` to show a dropdown menu
  - You will see your default branch (typically this is `main`)
  - To update the Identifiers so you can build the Trio 0.5.x, you must select `dev`
4. Tap the green button that says `Run workflow`.

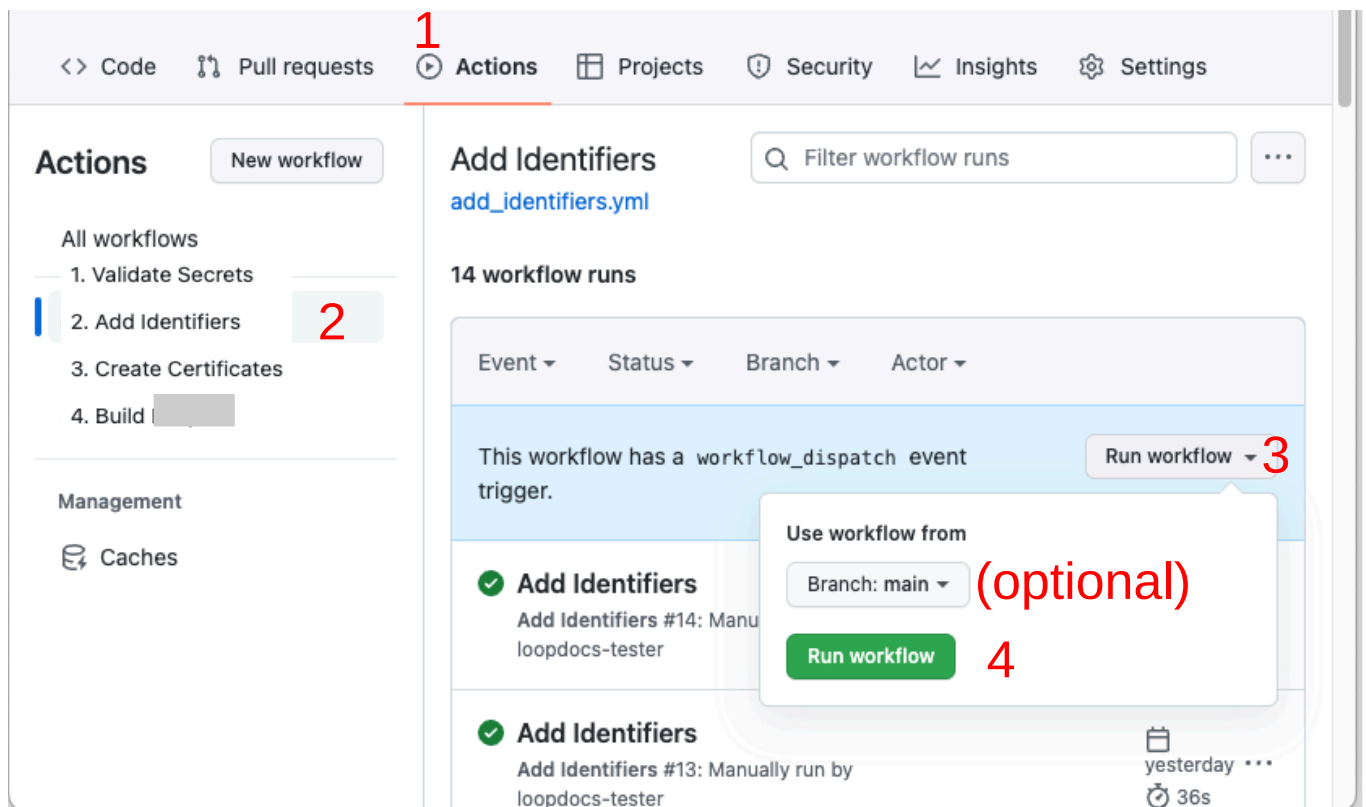

### Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes


The `Add Identifiers` Action should succeed or fail in a few minutes. Do not continue to the next step until this one succeeds.

## Update Certificates

For Trio 0.5.x and newer, you must have *Apple* Push Notification enabled to build the app. After you add this capability to the Trio Identifier, you must update the Certificate by running the Action: Create Certificates.

Refer to the graphic below for the numbered steps:

1. Click on the "Actions" tab of your `Trio` repository
2. On the left side, click on "Create Certificates"
3. On the right side, click "Run Workflow" to show a dropdown menu
  - You will see your default branch (typically `main`)
  - Because you plan to build Trio 0.5.x, select `dev`
4. Tap the green button that says "Run workflow".

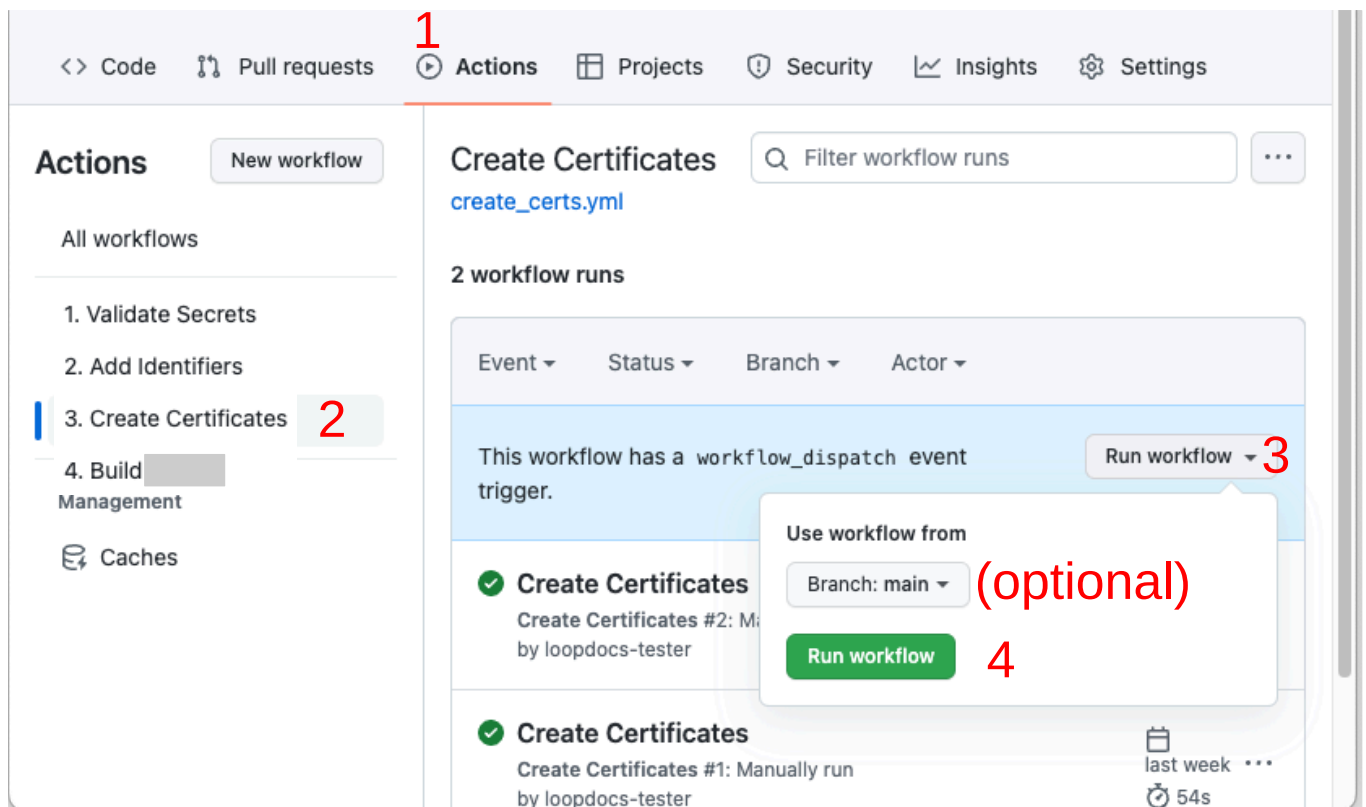

5. Wait a minute or two for the action to finish

#### Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Once you see the green check mark by `Create Certificates`, the next step is to Build the app.

## Build the App

If you completed all the steps on this page successfully (got a green checkmark ) , you are ready to run Action: Build Trio.

**We recommend most users leave their default branch as `main` . That means you must remember to select `dev` branch each time you build Trio 0.5.x**

If you have one branch as default, for example `main` , and choose to build a different branch, there is an extra step when you `Build Trio` . In addition to the normal steps 1, 2 and 3 in the graphic below, you must also do the (optional) step. Select the `dev branch` in the `branch dropdown` menu before continuing to step 4 and tapping on the green Run workflow button:



<> Code

Pull requests

**1** Actions

Projects

Security

Insights

Settings

Actions

New workflow

All workflows

1. Validate Secrets

2. Add Identifiers

3. Create Certificates

**4. Build Trio** **2**

8. DONT RUN Add ... Disabled

8. DONT RUN close... Disabled

Management

Caches

4. Build Trio

build\_trio.yml

133 workflow runs

Event Status Branch Actor

This workflow has a workflow\_dispatch event trigger.

Run workflow **3**

Use workflow from

Branch: dev (optional)

Run workflow **4**

Build Trio (bast\_fe

4. Build Trio #133: Sched

Refresh, Do Not Repeat

Hit Refresh

After you tap the green Run workflow button, GitHub can be slow to update.

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Update Build Errors

The most likely build error is that you did not Update Identifiers and Certificates. In that case, you will see this error in the build action annotations.

Error

Annotations

2 errors and 1 warning

Build

Process completed with exit code 1.

Build: Trio.xcodeproj#L0

Provisioning profile "match AppStore org.nightscout.\*\*\*.trio" doesn't support the Push Notifications capability. (in target 'Trio' from project 'Trio')

Show less

If you decide to look at the log, instead of reading the annotations, you may see an error similar to the two graphics below.

```
691 [23:27:00]: Called from Fastfile at line 138
692 [23:27:00]: ```
693 [23:27:00]: 136:      )
694 [23:27:00]: 137:
695 [23:27:00]: => 138:      gym(
696 [23:27:00]: 139:          export_method: "app-store",
697 [23:27:00]: 140:          scheme: "Trio",
698 [23:27:00]: ```
699 [23:27:00]: Error building the application - see the log above
700
701 +-----+
702 | fastlane summary |
```

```
from project 'Firebase' )
518 [23:27:00]: ▶ ::error file=/Users/runner/work/Trio-dev/Trio-
dev/Trio.xcodeproj::Provisioning profile "match AppStore org.nightscout.***.trio" doesn't
support the Push Notifications capability. (in target 'Trio' from project 'Trio')
519 [23:27:00]: ▶ /Users/runner/work/Trio-dev/Trio-dev/Trio.xcodeproj: error: Provisioning
profile "match AppStore org.nightscout.***.trio" doesn't include the aps-environment
entitlement. Profile qualification is using entitlement definitions that may be out of
date. Connect to network to update. (in target 'Trio' from project 'Trio')
520 [23:27:00]: ▶ note: Disabling previews because SWIFT_VERSION is set and
SWIFT_OPTIMIZATION_LEVEL=-O, expected -Onone (in target 'Trio' from project 'Trio')
521 [23:27:00]: ▶ note: Run script build phase 'Swiftformat' will be run during every build
because the option to run the script phase "Based on dependency analysis" is unchecked.
```

## Solution

[Update Identifiers and Certificates.](#)