

Update with Browser 



TrioDocs

Version: 0.4.0

Date: June 29, 2025

Download the latest version on:
<https://triodocs.org>

Update with
Browser 

Update to Trio dev branch from Trio 0.2.x main

! A Note on Compatibility

Upgrading to Trio 0.5.x from 0.2.x is smooth and straightforward. The new Onboarding Wizard will guide you step by step. Your pump, CGM, therapy settings, and the last 24 hours of treatment and glucose history will be brought over automatically.

Be Aware

- **Once you have upgraded to 0.5.x, going back to 0.2.x is not supported**
 - If you choose to downgrade, you will need to set everything up again from scratch
- Some guardrails have changed compared to 0.2.x
- **Your saved Temp Targets and Overrides will not be maintained**
 - The storage method in 0.5.x is different from 0.2.x, so you will need to recreate them
 - Capture a screenshot of each named Temp Target or Override that you want to add to 0.5.x

We have designed this release to be a stable and long-term foundation—so when you are ready, you can upgrade with confidence.

The Browser Build documentation is under construction for versions Trio 0.5.x and newer.

This temporary set of instructions is intended for the subset of users who **previously built Trio 0.2.x using Browser Build**. Because this is aimed at experienced builders, some of the steps are abbreviated and will be expanded as the documents are updated.

If you are new to Trio - please wait until the documents have been expanded and preferably until the release of Trio 1.0.

Mac-Xcode to Browser Build

! Experienced Mac-Xcode Builder

If you are an experienced Trio 0.2.x user who wants to join the open beta testing, but also wants to switch from Mac-Xcode build to Browser Build - welcome. For the time being, please do this:

- Use the [0.2.x: Build Trio with GitHub](#) instructions to build Trio 0.2.x and wait for it to show up in your TestFlight. Do not install it - you just want to make sure you can succeed with a Browser Build.
- Then return to this page and follow the directions; you will need to run **Add Identifiers** and **Create Certificates** again when updating to 0.5.x.

Transition from Other Apps to Trio Browser Build

! Coming from Other Apps

If you are an experienced Browser Build person who has not built Trio before, you will need to add the Trio App Group to 3 Identifiers and create a Trio App in App Store Connect.

The steps are summarized in New Trio Builders.

If you are not experienced with Browser Build, we suggest you follow the 0.2.x instructions: 0.2.x: Build Trio with GitHub first. Then return to this page to update to 0.5.x; you will need to run Add Identifiers and Create Certificates again when updating to 0.5.x.

Summary of Tasks to Build Trio 0.5.x

These steps assume:

- You previously built Trio 0.2.x using GitHub Actions (Browser Build)
- You confirmed your Apple Developer license agreements are up to date

These are the new steps for you to follow:

1. Configure Browser Build Certificate Automation
2. Configure Fork with dev branch
3. Update Identifiers
4. Update Certificates
5. Build Trio 0.5.x

Browser Build Certificate Automation

Browser Build Certificate Automation was added to Trio 0.2.3 - if you have not added the ENABLE_NUKE_CERTS variable, you should add it now. Otherwise, skip ahead to Configure Fork with dev branch.

In order to utilize the new automatic certificate renewal feature, you'll need to add a new Variable. Variables are located in GitHub, in the same location as your Secrets. The exact location will depend upon whether you build using a GitHub organization or a personal account.

If you use a personal account, click on your Trio repository. If you have other repositories, just follow these same steps for each of them.

If you build using a GitHub organization and have already added this variable to your organization, there's nothing for you to do. All repositories in your organization are covered. Otherwise, click on your organization name.

The numbered steps correspond to numbers in the graphics below:


1. Choose settings
2. Scroll down to select Secret and variables
3. Choose Actions
4. Choose Variable

5. Tap on “Create new organization variable” or “Create new repository variable”

6. Add the variable `ENABLE_NUKE_CERTS`

7. Type in the value `true`

8. Tap on Add Variable.

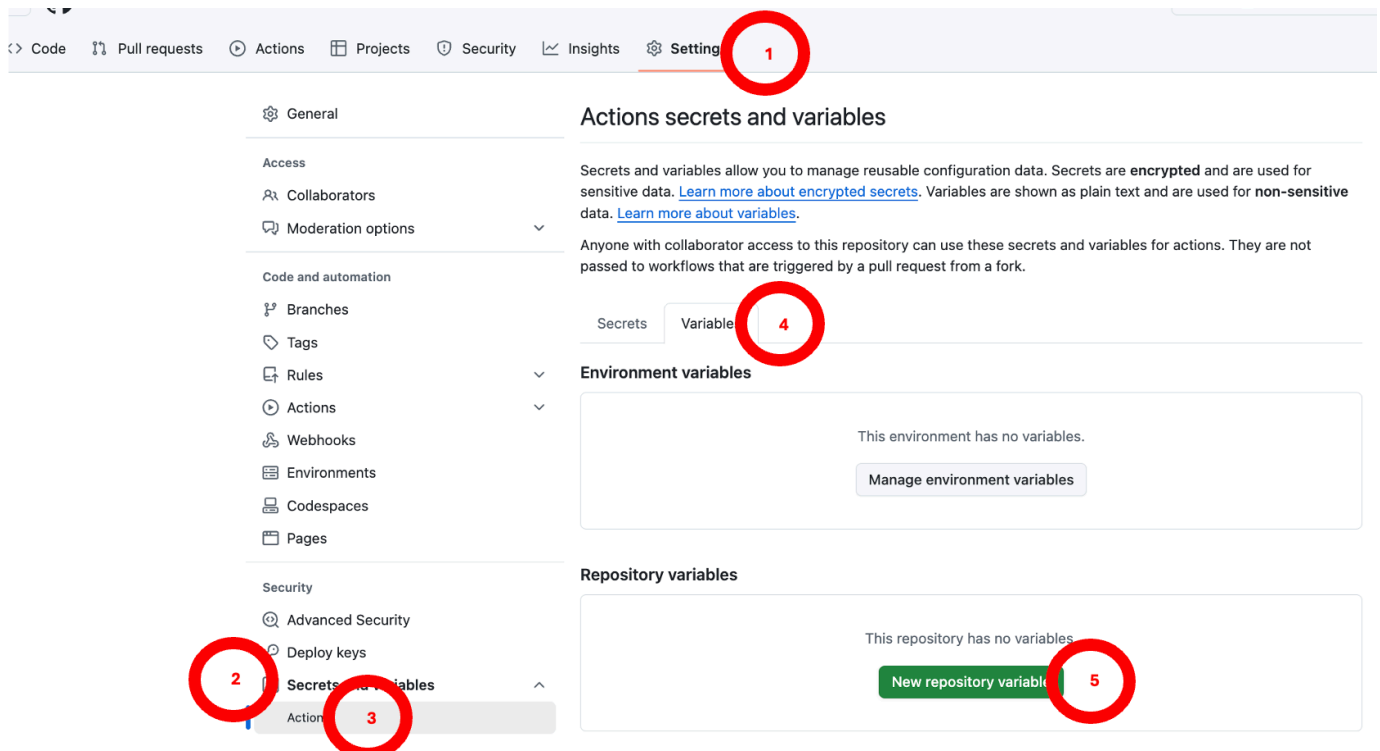

Actions variables / New variable

Note: Variable values are exposed as plain text. If you need to encrypt and mask sensitive information, [create a secret](#) instead.

Name *

`ENABLE_NUKE_CERTS`

- Variable names may only contain alphanumeric characters ([a-z], [A-Z], [0-9]) or underscores (_).
- Variable names cannot start with a number.
- Variable names cannot start with GITHUB_ prefix.

Value *

`true`

Add variable


Configure Fork with dev branch

The open beta testing for Trio uses the `dev` branch.


- If you do not have a `dev` branch you must first configure one following these directions.
- If you already have a `dev` branch, skip ahead to [Update Branch](#) and be sure to select the `dev` branch when you update your fork.

Add Branch

First, you need to get to your Trio repository and tap on the branch icon. (Do not worry about how many branches you have.):

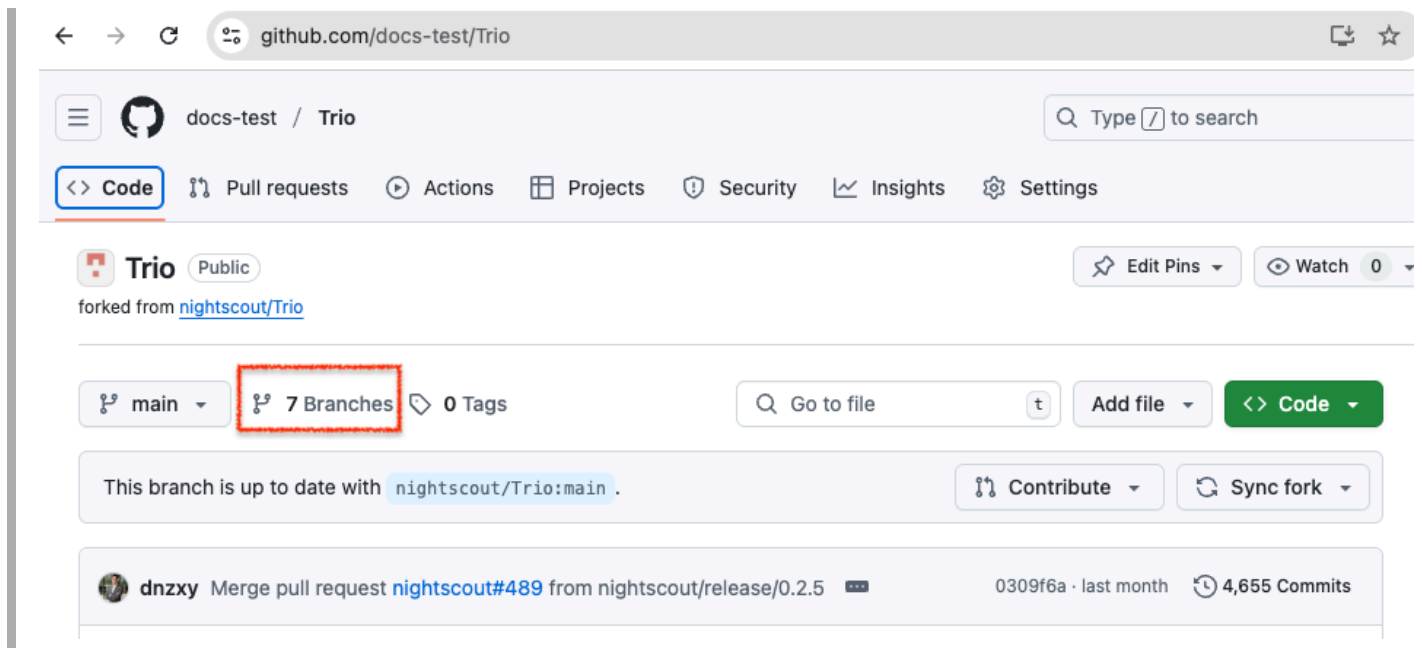

The next screen displays the branches you currently have. It shows a `New branch` icon in the upper right. If you don't see that, you are not logged in. Tap on the `New branch` button.

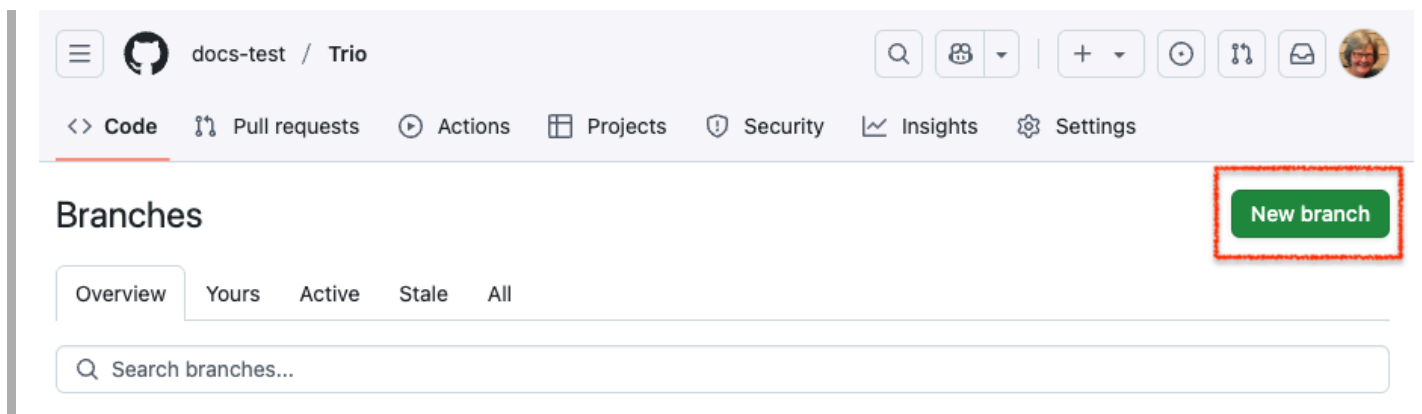


Each step in the list below matches with the number in the graphic. On the left half of the graphic is the default selections for your fork when `main` is the default branch. The right side shows the display after making the indicated selections.:

1. Click on the drop-down menu labeled 1 in the graphic and choose `nightscout/Trio`
2. Click on the drop-down menu labeled 2 in the graphic and choose `dev`
3. Click on the `Branch` name box labeled 3 in the graphic and type `dev`

- The branch name in your fork should always match the branch name you are adding; check that you type it correctly

4. Review the dialog items to make sure everything is correct, and then tap on Create branch

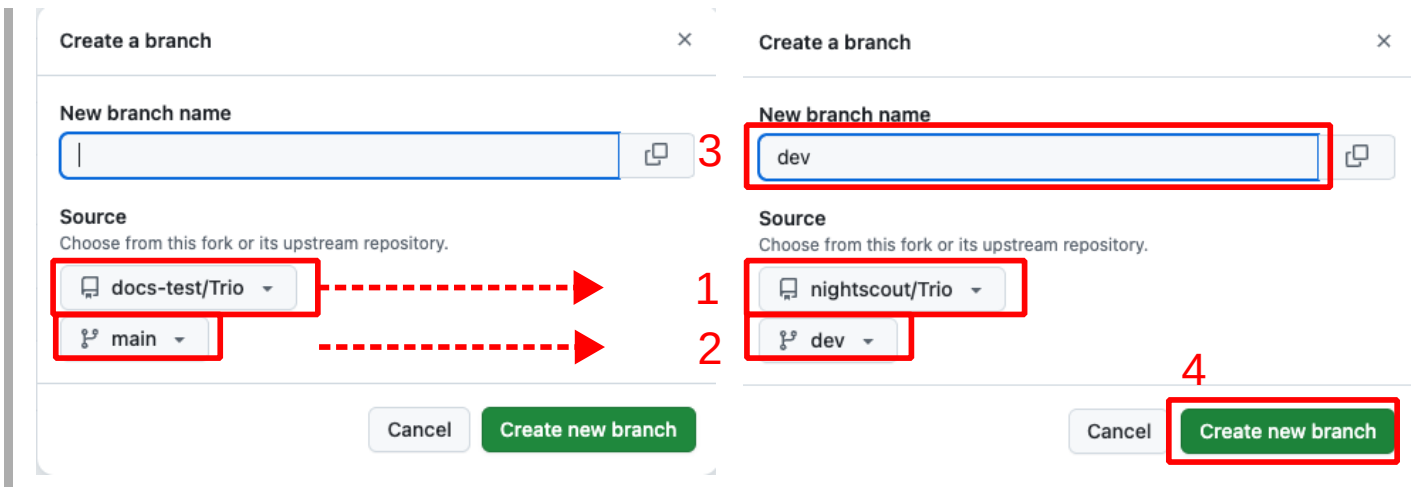


Update Branch

Be aware that building from the dev branch will upgrade to Trio 0.5.x. Once you have upgraded to 0.5.x, going back to 0.2.x is not supported.

Tap the Code button (upper left) and ensure this branch in your fork is up to date.

- Select the desired branch in the dropdown menu (this graphic shows main branch, to get 0.5.x, you must choose dev branch)
- If the message indicates this branch is "behind", tap on the sync fork button and then the Update branch button

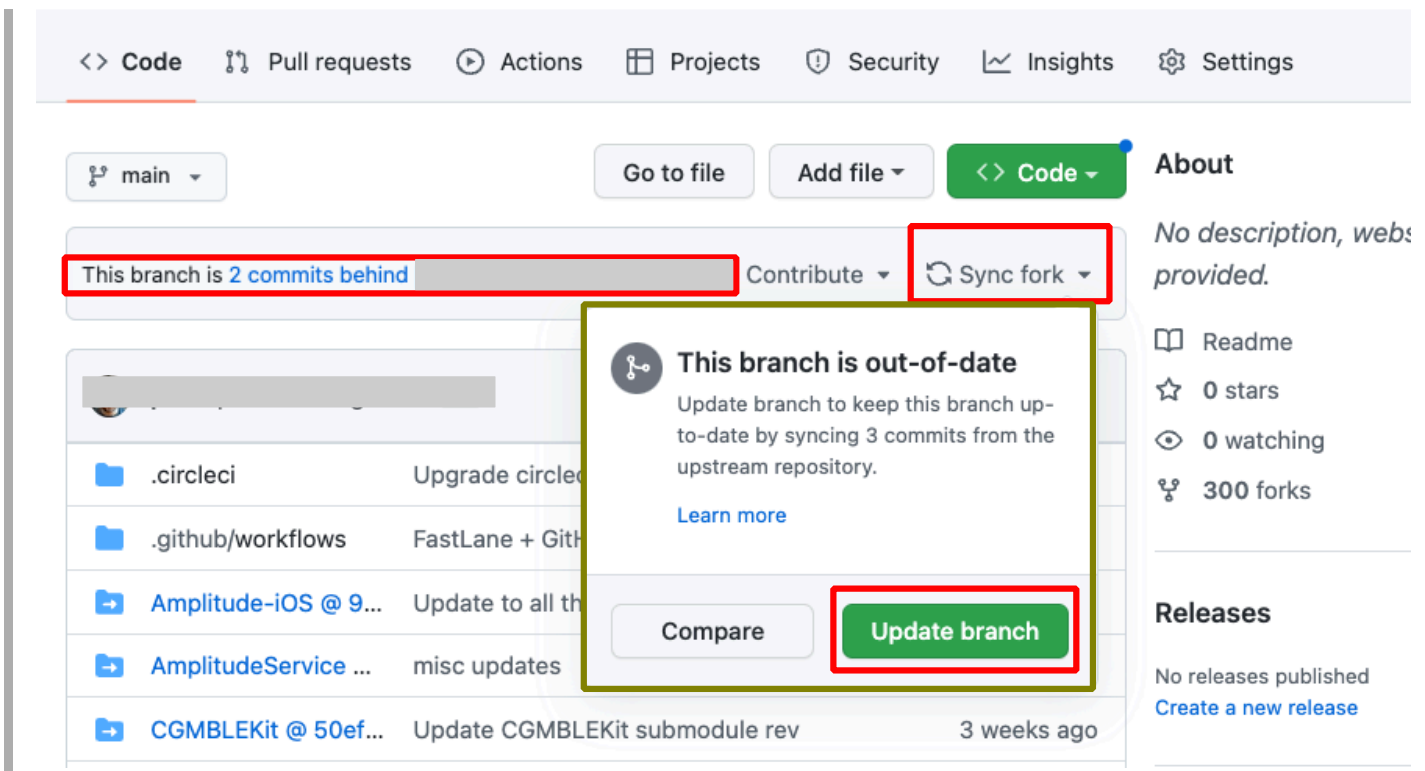


Update Identifiers

For Trio 0.5.x and newer, you must have *Apple* Push Notification enabled to build the app. This capability is added to the existing Trio Identifier by running the Action: `Add Identifiers` after you update your fork.

Refer to the graphic below for the numbered steps:

1. Click on the `Actions` tab of your `Trio` repository
2. On the left side, click on 2. `Add Identifiers`
3. On the right side, click `Run Workflow` to show a dropdown menu
 - You will see your default branch (typically this is `main`)
 - **To update the `Identifiers` so you can build the Trio 0.5.x, you must select `dev`**
4. Tap the green button that says `Run workflow`.

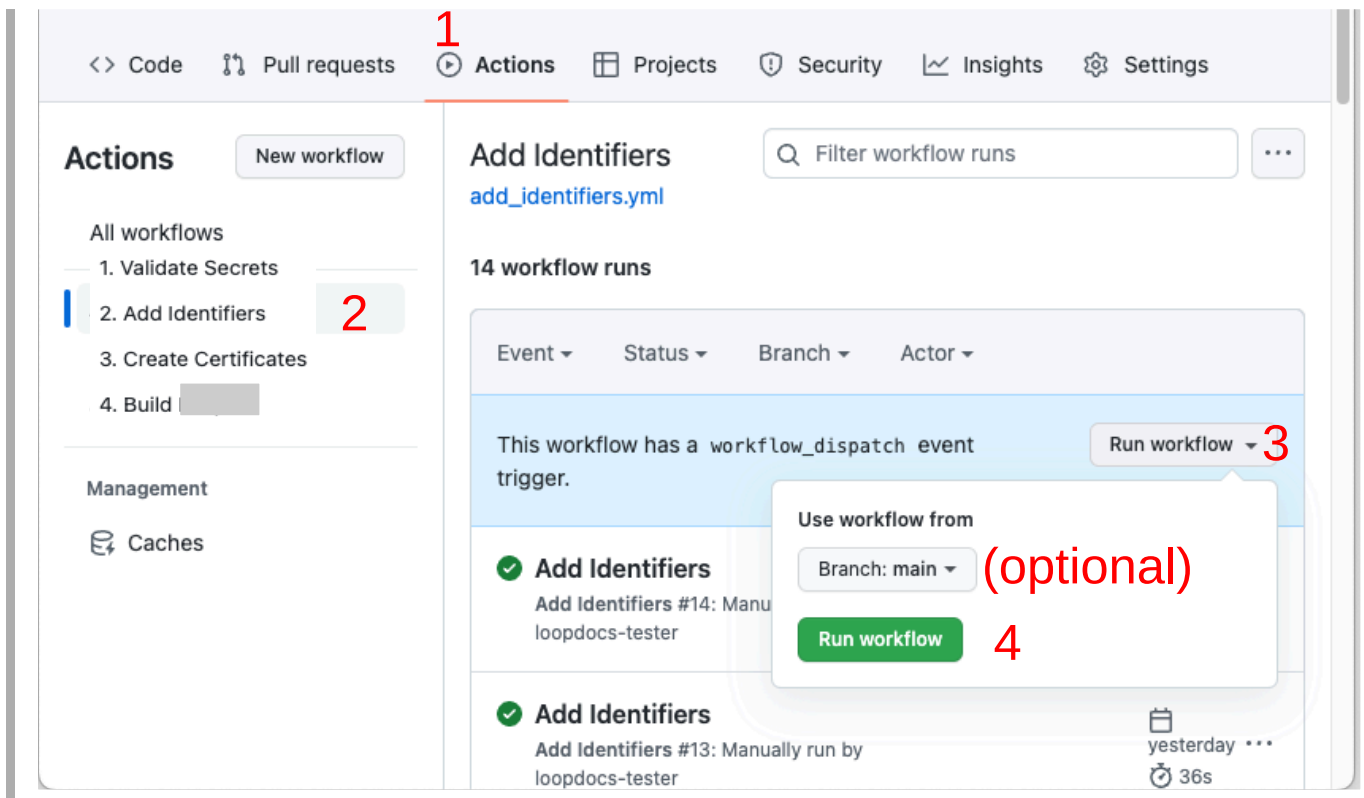

Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes


The `Add Identifiers` Action should succeed or fail in a few minutes. Do not continue to the next step until this one succeeds.

Update Certificates

For Trio 0.5.x and newer, you must have *Apple* Push Notification enabled to build the app. After you add this capability to the Trio Identifier, you must update the Certificate by running the Action: `Create Certificates`.

Refer to the graphic below for the numbered steps:

1. Click on the "Actions" tab of your Trio repository
2. On the left side, click on "Create Certificates"
3. On the right side, click "Run Workflow" to show a dropdown menu
 - You will see your default branch (typically main)
 - **Because you plan to build Trio 0.5.x, select dev**
4. Tap the green button that says "Run workflow".

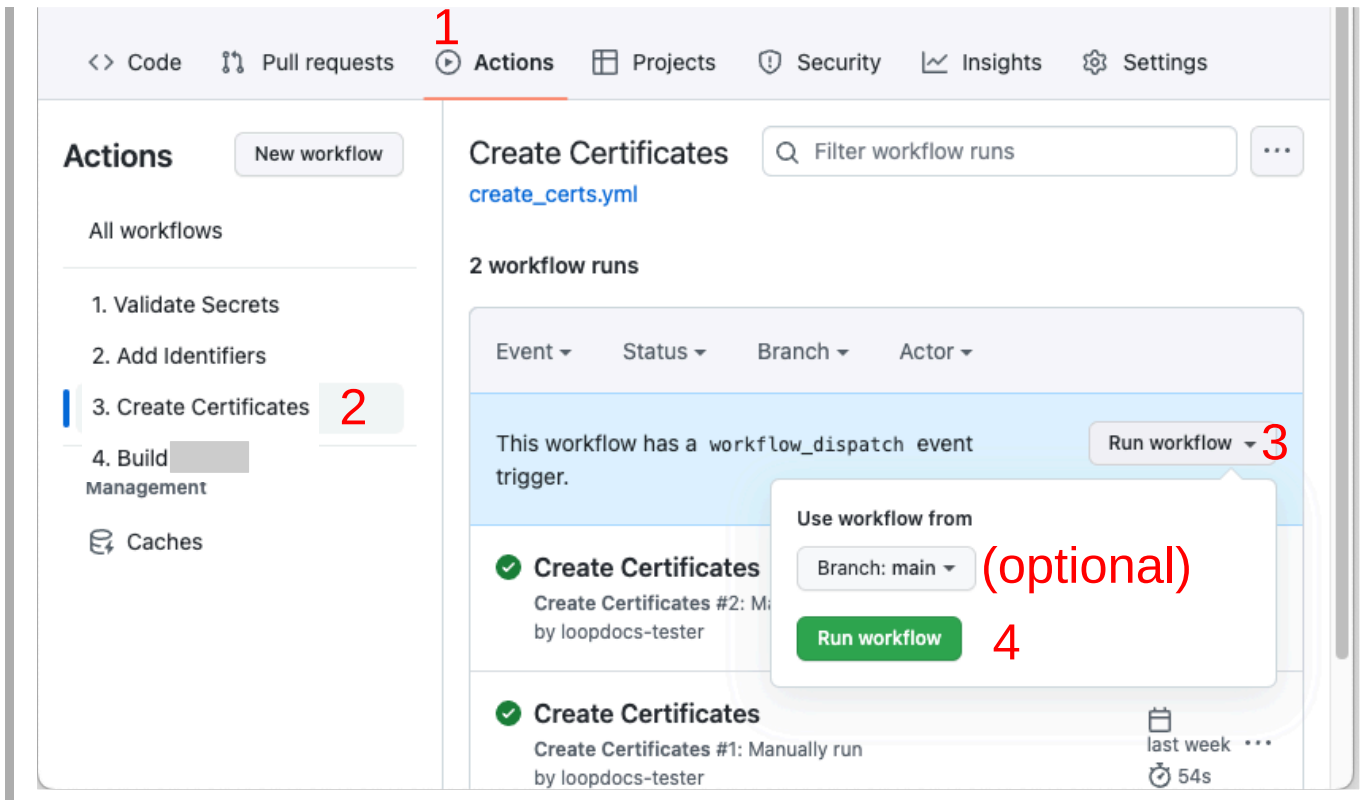

5. Wait a minute or two for the action to finish

Be Patient

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Once you see the green check mark by `Create Certificates`, the next step is to Build the app.

Build the App

If you completed all the steps on this page successfully (got a green checkmark ) , you are ready to run Action: Build Trio.

We recommend most users leave their default branch as `main` . That means you must remember to select the `dev` branch each time you build Trio 0.5.x

If you have one branch as default, for example `main` , and choose to build a different branch, there is an extra step when you `Build Trio` . In addition to the normal steps 1, 2, and 3 in the graphic below, you must also do the

(optional) step. Select the `dev` branch in the `branch` dropdown menu before continuing to step 4 and tapping on the green `Run workflow` button:

The screenshot shows the GitHub Actions interface. The top navigation bar has tabs for Code, Pull requests, Actions (highlighted with a red '1'), Projects, Security, Insights, and Settings. On the left, under 'Actions', there is a list of workflows: '1. Validate Secrets', '2. Add Identifiers', '3. Create Certificates', '4. Build Trio' (highlighted with a red '2'), '8. DONT RUN Add ... Disabled', and '8. DONT RUN close... Disabled'. Below this is a 'Management' section with 'Caches'. The main area shows the details for '4. Build Trio' (workflow file: `build_trio.yml`). It indicates '133 workflow runs'. A message states: 'This workflow has a `workflow_dispatch` event trigger.' There is a 'Run workflow' button (highlighted with a red '3'). A modal is open, titled 'Use workflow from', showing 'Branch: dev' (with a red '(optional)' next to it) and a 'Run workflow' button (highlighted with a red '4').

Refresh, Do Not Repeat

Hit Refresh


After you tap the green `Run workflow` button, *GitHub* can be slow to update.

- Refresh the browser if you are unsure if the action started
- Do not start a new action until the first one completes

Disable the DO NOT RUN Actions

You may notice some *GitHub* actions that have `DONT RUN` or `DO NOT RUN` in their names. Those are special actions that are used by the developers. They are skipped in your `Fork`, but you may see logs for them.

If they already have the notation `disabled` beside them, so won't see them run. Otherwise, click on each one, click on the three dots at the upper right, and then select `Disable workflow`. You will not be annoyed by seeing logs that say the action was automatically run and then skipped.

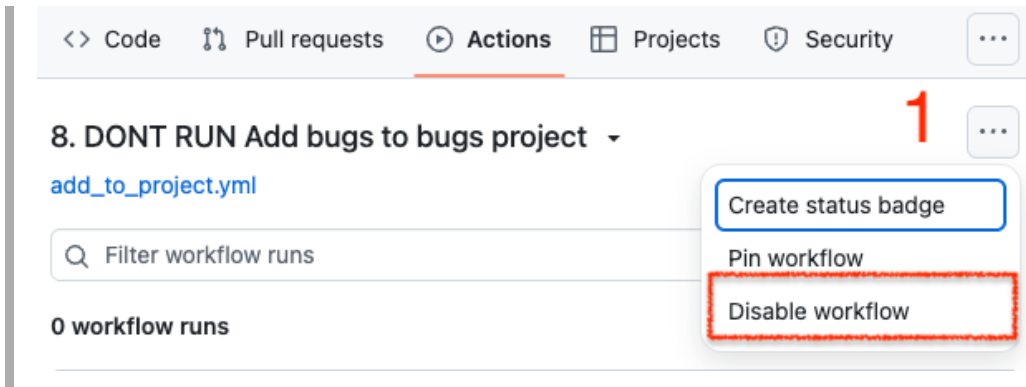

Update Build Errors

The most likely build error is that you did not [Update Identifiers and Certificates](#). In that case, you will see this **error** in the **build action annotations**.

Error



If you decide to look at the action **log**, instead of reading the annotations, you may see an **error** similar to the two graphics below.



```
from project 'Firebase'
518 [23:27:00]: ▶ ::error file=/Users/runner/work/Trio-dev/Trio-dev/Trio.xcodeproj::Provisioning profile "match AppStore org.nightscout.***.trio" doesn't
support the Push Notifications capability. (in target 'Trio' from project 'Trio')
519 [23:27:00]: ▶ /Users/runner/work/Trio-dev/Trio-dev/Trio.xcodeproj: error: Provisioning
profile "match AppStore org.nightscout.***.trio" doesn't include the aps-environment
entitlement. Profile qualification is using entitlement definitions that may be out of
date. Connect to network to update. (in target 'Trio' from project 'Trio')
520 [23:27:00]: ▶ note: Disabling previews because SWIFT_VERSION is set and
SWIFT_OPTIMIZATION_LEVEL=-O, expected -Onone (in target 'Trio' from project 'Trio')
521 [23:27:00]: ▶ note: Run script build phase 'Swiftformat' will be run during every build
because the option to run the script phase "Based on dependency analysis" is unchecked.
```

Solution

1. Update [Identifiers](#) - make sure you select `dev` branch
 2. Update [Certificates](#) - make sure you select `dev` branch
-

New Trio Builders

For experienced browser builders, this minimal set of instructions might be sufficient.

Warning - there is no hand holding in these directions.

Every app you build will use the same 6 Secrets.

- Fork from: <https://github.com/nightscout/Trio>
- You use the `Trio App Group` for Trio, see [Create the Trio App Group](#)
- [Identifiers](#) for Trio: see [Table of Identifiers](#)
- Add the `Trio App Group` to these [Identifiers](#):
 - `Trio`
 - `Trio Watch`
 - `Trio WatchKit Extension`
- In `App Store Connect`, the `Bundle ID` for Trio will be: `org.nightscout.TEAMID.trio`

Create the `Trio App Group`

If you already have a `Trio App Group`

- You can skip this step - your existing App Groups are found at this link: [App Group List](#)
- If your `Trio App Group` was created from a Mac with `Xcode`, you may choose to edit the Description to make the **NAME** match

If you do not have a `Trio App Group`:

- Go to [Register an App Group](#) on the Apple developer site and use the table below to help you create one.
- Replace `TEAMID` with your *Apple Developer* ID.

NAME	Xcode version (NAME)	IDENTIFIER
Trio App Group	group org nightscout TEAMID trio trio-app-group	group.org.nightscout.TEAMID.trio.trio-app-group


Table of Identifiers

These are the Identifiers created by running the GitHub action "Add Identifiers". Only 3 of them need to have the Trio App Group added. Be sure that you selected the dev branch when you start the "Add Identifiers" action. For more details, see Update Identifiers , but then return here when that action succeeds.

NAME	Xcode version (NAME)	IDENTIFIER
Trio	XC org nightscout TEAMID trio	org.nightscout.TEAMID.trio
Trio LiveActivity	-	org.nightscout.TEAMID.trio.LiveActivity
Trio Watch App	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp
Trio Watch Complication	XC IDENTIFIER	org.nightscout.TEAMID.trio.watchkitapp.TrioWatchComplication

Add Trio App Group to Identifiers

- Open the App IDs Identifier page for your Apple Developer Account.
- Click on the Trio Identifier and assign the Trio App Group to the Identifier - see graphic below.

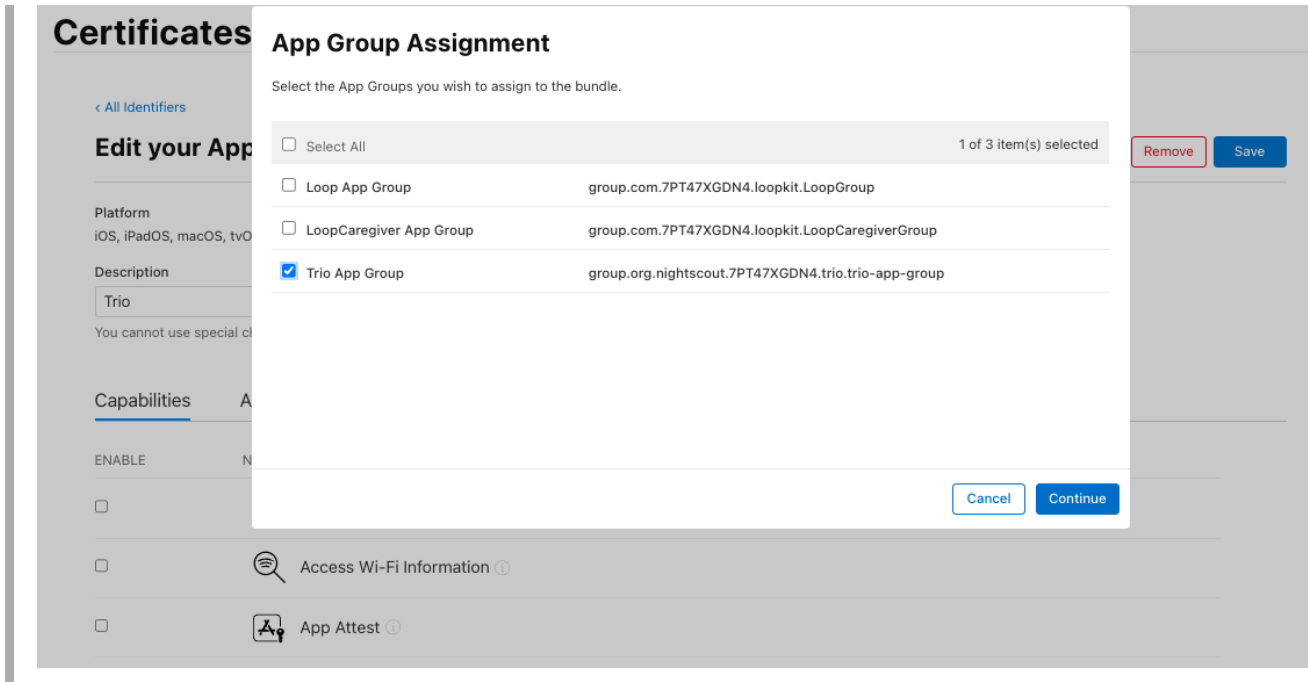

- Repeat this for the other 2 `Identifiers` that need to have an App Group assigned: `Trio Watch` and `Trio WatchKit Extension`

Configure Trio App

Follow the directions in [LoopDocs](#), but use the Trio `Bundle ID`

- In `App Store Connect`, the `Bundle ID` for Trio will be: `org.nightscout.TEAMID.trio`

Return to the main set of instructions on this page, [Update Certificates](#), and keep going until you have a successful build.